



HEAD OFFICE: E&M ACTIVE LIMITED, THE INNOVATION CENTRE, HERE EAST, QUEEN ELIZABETH OLYMPIC PARK, LONDON, E20 3BS,  
UNITED KINGDOM

[WWW.EMACTIVE.COM](http://WWW.EMACTIVE.COM)

MANUAL

**MAGICTOUCH**<sup>®</sup>  
O P T I C M A S S A G E R

## A BIG HELLO FROM E&M ACTIVE®

**THANK YOU FOR CHOOSING E&M ACTIVE®.** WE ARE DEVOTED TO ENHANCING WELL-BEING. OUR MISSION IS TO OFFER EVERYONE WORLDWIDE THE BEST POSSIBLE PRODUCTS THAT MEET THE HIGHEST STANDARD OF QUALITY, EFFICACY AND DESIGN.

**WE HOPE YOU ENJOY USING YOUR NEW E&M ACTIVE® PRODUCT.** PLEASE READ THE INSTRUCTION MANUAL THOROUGHLY BEFORE USE AND KEEP IT IN A SAFE PLACE FOR FUTURE REFERENCE.

**AT E&M ACTIVE® WE FRAME HEALTH&BEAUTY AND MATERIALISE LOVE.**

AS THIS IS OUR PHILOSOPHY,  
LIFE SHOULD BE LIVED WITH PLEASURE,  
LIKEWISE HEALTH&BEAUTY SHOULD BE INSPIRED BY LOVE.

**FIND OUT MORE ABOUT OUR BRAND.** VISIT US AT [WWW.EMACTIVE.COM](http://WWW.EMACTIVE.COM)

## LET'S GET STARTED

BEFORE YOU START .....	3
YOUR PRODUCT .....	4
SAFETY .....	13
TECHNICAL SPECIFICATIONS .....	14
CARE AND MAINTENANCE .....	15

## BEFORE YOU USE YOUR MAGICTOUCH® OPTIC MASSAGER, YOU NEED TO KNOW

- ✓ THE APPLIANCE IS FOR DOMESTIC USE ONLY.
- ✓ THIS IS NOT A TOY, IT SHOULD BE USED BY ADULTS ONLY.
- ✓ PLEASE KEEP OUT OF THE REACH OF CHILDREN AND PETS AT ALL TIMES, EVEN WHEN IN USE.

THE APPLIANCE COMES FULLY ASSEMBLED. BEFORE USING PLEASE ENSURE THAT ALL COMPONENTS AND FIXINGS ARE SECURE AND MAKE SURE THE APPLIANCE IS IN PERFECT CONDITION BEFORE OPERATION. PLEASE DO NOT USE ANY ATTACHMENTS THAT ARE NOT RECOMMENDED BY E&M ACTIVE®.



SOFT CLOTH INNER SURFACE

ADJUSTABLE FASTENING BELT



DISPLAY SCREEN

TWO ATTACHMENTS:




USB CHARGER



SOFT CLOTH BAG

### CHARGING THE BATTERY

- CONNECT THE USB CHARGER WITH A COMPUTER OR ADAPTER WITH THE SAME USB INPUT DOCK.
- THE LIGHT WILL TURN RED WHEN THE APPLIANCE BEGINS TO CHARGE. ONCE IT HAS FINISHED CHARGING, THE LIGHT WILL TURN GREEN.
- YOU MAY NOT USE THE APPLIANCE WHILE IT IS PLUGGED IN TO CHARGE. IF YOU PLUG IT IN DURING CHARGE, IT WILL AUTOMATICALLY TURN ITSELF OFF.
- PLEASE USE THE ADAPTER PROVIDED. DO NOT USE THE USB 2.0 PORT TO CHARGE THE BATTERY.
- THE LED LIGHT WILL BLINK TO ALERT YOU WHEN THE BATTERY NEEDS TO BE CHARGED.
- IT MAY TAKE 2 HOURS TO FULLY CHARGE THE BATTERY.

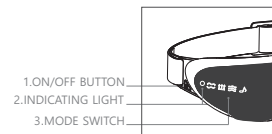

**MAGICTOUCH®**  
OPTIC MASSAGER

OPTION 1: PRESSURING MESSAGE, HEATING, VIBRATION AND BLUETOOTH.  
OPTION 2: PRESSURING MESSAGE, HEATING AND BLUETOOTH.  
OPTION 3: PRESSURING MESSAGE AND BLUETOOTH ONLY.



SMART AIR  
PRESSURE KNEADING



FAR INFRARED  
HEATING



VIBRATION



BLUETOOTH  
CONNECTION

## PREPARATION BEFORE USE

- ✓ OPEN THE BOX AND REMOVE THE APPLIANCE.
- ✓ REMOVE ALL OF THE PROTECTIVE PACKAGING.
- ✓ MAKE SURE THAT THE VOLTAGE COMING FROM THE OUTLET IS IDENTICAL TO THE ONE ON THE RATING LABEL.
- ✓ CHARGE THE APPLIANCE FULLY BEFORE USE.
- ✓ UNPLUG ONCE CHARGED AND PRESS THE POWER BUTTON TO SWITCH THE APPLIANCE ON.



ON/OFF SWITCH

USE THIS SWITCH TO TURN THE APPLIANCE ON. IT WILL THEN PRODUCE A SINGLE HIGH-PITCHED SOUND TO INDICATE THAT IT IS READY. YOU CAN TAP THIS BUTTON TO CHANGE THE SETTINGS OF THE APPLIANCE; THE DISPLAY LIGHTS WILL CHANGE ACCORDING TO WHICH LEVEL YOU ARE IN. HOLD THE BUTTON DOWN IN A LONG PRESS TO SWITCH THE APPLIANCE OFF.



MAGNETIC FIELD	TEMPERATURE	AIR PRESSURE	VIBRATION
1100*2 GS (INDUCTION STRENGTH)	39±3°C	MAX 45 KPa	13500±15% rpm

**TO CONNECT THE MAGICTOUCH® OPTIC MASSAGER BLUETOOTH SPEAKERS TO YOUR PHONE OR OTHER BLUETOOTH-ENABLED DEVICE:**

**BLUETOOTH IS ALREADY ENABLED ON THE OPTIC MASSAGER BY TURNING ON THE DEVICE.** PRESS THE ON/OFF BUTTON UNTIL YOU SEE THAT MUSIC ICON IS SELECTED.

**ACTIVATE BLUETOOTH ON YOUR DEVICE.** GO TO SETTINGS, FIND BLUETOOTH AND CLICK TO TURN IT ON.

**MAKE THE CONNECTION.** THE MAGICTOUCH NAME SHOULD NOW APPEAR ON YOUR DEVICE. (IF IT DOES NOT APPEAR, MAKE SURE THAT THE MAGICTOUCH® OPTIC MASSAGER ISN'T CONNECTED TO ANOTHER OF YOUR DEVICES).

**PAIR THE TWO DEVICES.** PRESS MAGICTOUCH NAME SHOWN ON YOUR DEVICE SCREEN TO PAIR. WAIT A FEW SECONDS... NOW YOU'LL HEAR 'PAIRING SUCCESSFULLY' FROM YOUR MAGICTOUCH. YOU CAN NOW RELAX AND ENJOY YOUR MAGICTOUCH® OPTIC MASSAGER!



**MAGICTOUCH®**  
O P T I C M A S S A G E R

### BENEFITS:

YOUR OPTIC MASSAGER COMBINES SMART PRESSURING MASSAGE, HEAT AND VIBRATION TO HELP RELAX THE EYE AREA.


YOUR OPTIC MASSAGER CAN HELP TO:

- ✓ RELIEVE TEMPLE PRESSURE AND STRAINS SURROUNDING THE EYE WHICH MAY CAUSE HEADACKES;
- ✓ IMPROVE CIRCULATION;
- ✓ REDUCE THE APPEARANCE OF DARK CIRCLES AND PUFFINESS;
- ✓ PROMOTE BETTER SLEEP;
- ✓ ELEVATE THE MOOD AND CALM THE MIND;
- ✓ AID DRY EYES;
- ✓ RELIEVE FATIGUE.

## SPINNING MAGNET MASSAGE

WE USE SPINNING MAGNET MASSAGE IN OUR OPTIC MASSAGER, WHICH COMBINES THE BENEFITS AND EFFICACY OF DEEP TISSUE MASSAGE WITH THE GENTLE FEEL OF A RELAXING MASSAGE.

WHAT DOES IT DO? WE ARE ALL MADE UP OF ELECTROMAGNETIC FREQUENCIES (OR VIBRATIONS) SO USING MAGNETIC MASSAGE AND THERAPY HAS QUICK AND POSITIVE EFFECTS ON THE BODY. HELPING CELLS TO REGAIN THEIR NATURAL ELECTROMAGNETIC FIELD CAN HELP TO ACTIVATE THE BODY'S NATURAL HEALING MECHANISMS.

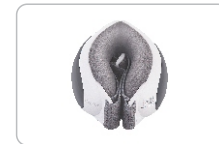


## SPINNING MAGNET MASSAGE

HOW DOES IT WORK? MAGNETIC FIELDS ARE MADE FROM THE NORTH AND SOUTH POLES OF THE MAGNET PUSHING AGAINST EACH OTHER AND EACH SIDE HAS A DIFFERENT USE.

NORTH POLE MAGNETS ARE USED TO DILATE BLOOD VESSELS, WHICH HELPS WITH ANY PROBLEMS CAUSED BY POOR CIRCULATION SUCH AS MIGRAINES AND STRESS HEADACHES.

SOUTH POLE MAGNETS CONSTRICT BLOOD VESSELS WHICH CAUSES THEM TO CREATE AN ANTI-INFLAMMATORY, REDUCING FLUID AND SWELLING.

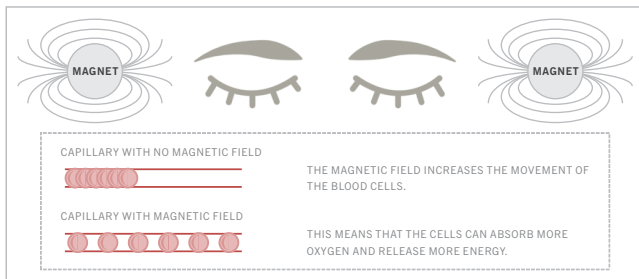

## SPINNING MAGNET MASSAGE



**MAGICTOUCH®**  
O P T I C M A S S A G E R

THE ELECTROMAGNETIC FIELDS IN OUR OPTIC MASSAGER BOOST THE PRODUCTION OF NEGATIVE IONS WHICH HELP YOU TO FEEL REFRESHED. NEGATIVE IONS ARE PRODUCED BY NATURE - THAT REFRESHING FEELING THAT YOU GET BY A WATERFALL, ON THE BEACH OR AFTER A STORM IS CAUSED BY YOUR BODY ABSORBING NEGATIVE IONS.

STEAM FROM HOT SHOWERS PRODUCES NEGATIVE IONS. THIS IS WHY A HOT SHOWER IS OFTEN THE ONLY EFFECTIVE WAY TO WAKE YOU UP IN THE MORNINGS.

NEGATIVE IONS HAVE SO MANY POSITIVE EFFECTS INCLUDING BUT CERTAINLY NOT LIMITED TO:

- ✓ INCREASING WELL-BEING AND MENTAL CLARITY
- ✓ IMPROVING ENERGY LEVELS AND FOCUS
- ✓ PROMOTING BETTER SLEEP
- ✓ REDUCING HEADACHES AND SICKNESS

**MAGICTOUCH® OPTIC MASSAGER  
CAN RELIEVE EYE STRAIN CAUSED  
BY EVERYDAY ACTIVITIES INCLUDING:**



USING COMPUTERS



WATCHING TV



USING SMARTPHONES



DRIVING



READING

## MAKING THE MOST OF YOUR PRODUCT:

THE APPLIANCE IS WIRELESS SO YOU CAN USE IT ANYWHERE AROUND THE HOUSE OR ON THE GO.

OUR OPTIC MASSAGER IS CUSTOMISEABLE – CHOOSE THE CONTROLS TO USE THE RIGHT SETTING FOR YOU.

THE SMART PRESSURISING MESSAGE ADJUSTS TO THE SHAPE OF YOUR FACE FOR THE ULTIMATE PERSONALISED COMFORT.

OUR OPTIC MASSAGER MASSES AROUND THE EYES SO THAT IT DOES NOT PUT ANY PRESSURE ON THE EYES THEMSELVES. THIS ALSO MEANS THAT YOU CAN SEE THROUGH THE GLASS WHILE WEARING THEM.

MOVING MAGNET MESSAGE PUSHES ENERGY MUCH DEEPER INTO AND THROUGH THE BODY HELPING TO CONTROL SOFT TISSUE PAIN AND GIVING AN ENERGISING EFFECT.

## BEFORE YOU USE THE PRODUCT:

PLEASE CAREFULLY READ THE FOLLOWING INSTRUCTIONS BEFORE YOU USE THE APPLIANCE.

THE PRODUCT MUST NOT BE USED WHILE DRIVING OR PERFORMING ANY OTHER ACTIVITY THAT REQUIRES FULL VISION.

PLEASE DO NOT USE THE APPLIANCE IN ANY ENVIRONMENT THAT IS EXTREMELY HOT OR HUMID.

DO NOT USE THE APPLIANCE IF YOU HAVE ANY EYE CONDITIONS OR HAVE JUST UNDERGONE EYE SURGERY. THIS INCLUDES GLAUCOMA, CATARACTS AND AMOTIO RENTINAE.

THE APPLIANCE IS NOT INTENDED FOR THE USE BY ANYONE WITH REDUCED PHYSICAL, SENSORY, MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE. THIS INCLUDES CHILDREN.

THE APPLIANCE SHOULD BE COMFORTABLE TO USE, DO NOT USE IF YOU FEEL ANY DISCOMFORT.

THE APPLIANCE IS NOT SUITABLE FOR USE WITH CONTACT LENSES; PLEASE REMOVE ANY CONTACT LENSES BEFORE USING.

THE APPLIANCE MUST BE HELD IN THE HAND DURING NORMAL USE.

TO CLEAN THE APPLIANCE, USE A SOFT, DRY CLOTH ON THE SURFACE, NEVER USE A WET OR DAMP CLOTH. WE ALSO RECOMMEND THAT YOU PERIODICALLY CLEAN THE APPLIANCE BY LIGHTLY BRUSHING OVER THE SOFT CLOTH INNER TO PREVENT ANY DUST FROM BUILDING UP.

BECAUSE THE APPLIANCE USES HEAT, ANYONE WHO IS SENSITIVE TO HEAT MUST BE CAREFUL WHEN USING THE APPLIANCE.

BATTERY SAFETY:

THE BATTERY MUST BE REMOVED FROM THE APPLIANCE IN ORDER TO DISPOSE OF IT SAFELY. PLEASE ENSURE THAT YOU DISCONNECT THE APPLIANCE FROM THE SUPPLY MAINS BEFORE REMOVING THE BATTERY.

THE BATTERY MUST BE DISPOSED OF SAFELY.

WHEN REPLACING THE BATTERY, ONLY REPLACE WITH THE SAME OR EQUIVALENT BATTERY TYPE. DANGER OF EXPLOSION IF THE INCORRECT BATTERY IS USED.

BATTERIES SHOULD NOT BE EXPOSED TO ANY EXCESSIVE HEAT, INCLUDING SUNSHINE, FIRE, OR OTHER HEAT SOURCES.

## TECHNICAL INFORMATION

BRAND	E&M ACTIVE/MAGICTOUCH®
PRODUCT	OPTIC MASSAGER
MODEL	EM01
INPUT	5V --- 0.8A
BUILT-IN BATTERY	3.7V 1150mAh 4.26Wh
RATED TIME	15 MINUTES
OPERATION FREQUENCY RANGE	2402-2480Mhz
MAX POWER	xxdBm
PRODUCT DIM	220*125*81(MM)
PACKAGING DIM	186*140*105(MM)
G.W	0. 7KG
N.W	0. 3KG

\*FOR MORE DETAILS, PLEASE REFER TO THE ON PRODUCT LABEL.  
DESIGNED & DEVELOPED BY E&M ACTIVE FROM LONDON, UK.



MADE IN CHINA.  
THIS APPLIANCE CONFORMS TO ALL RELEVANT EUROPEAN AND INTERNATIONAL DIRECTIVES.

FCC ID: 2AEHM-EM01 IC: 20096-EM01

THIS APPLIANCE CONFORMS WITH HEALTH, SAFETY AND ENVIRONMENTAL STANDARDS FOR THE EUROPEAN ECONOMIC AREA, INCLUDING WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT DIRECTIVE, AND HAS A CLASS III ELECTRIC SHOCK PREVENTION CERTIFICATION.

## CARE AND MAINTENANCE

1.ALWAYS DISCONNECT THE MAINS POWER SUPPLY BEFORE CLEANING THE APPLIANCE.

2.DO NOT POUR WATER OR ANY OTHER LIQUID ONTO OR INTO THE APPLIANCE.

3.DO NOT USE A WET CLOTH TO CLEAN THE APPLIANCE. USE A SOFT, DRY CLOTH TO CLEAN THE OUTER SURFACE. WE WOULD ALSO SUGGEST REGULARLY CLEANING THE APPLIANCE BY LIGHTLY BRUSHING OVER THE SOFT CLOTH TO PREVENT ANY DUST BUILD-UP.

4.WE SUGGEST THAT YOU STORE YOUR APPLIANCE IN THE ORIGINAL BOX WHILE NOT IN USE. IT IS IMPORTANT TO KEEP IT IN A SAFE AND DRY LOCATION, AWAY FROM ANY DUST.

## WARRANTY INFORMATION

THE PARTNER OF E&M ACTIVE® GRANTS STATUTORY WARRANTY PERIOD. THE APPLIANCE CARRIES A MINIMUM GUARANTEE PERIOD OF ONE(1) YEAR FROM THE DATE SHOWN ON THE ORIGINAL PURCHASE INVOICE OR SALES RECEIPT. IF YOU NEED HELP WITH USING YOUR APPLIANCE, SERVICE OR REPAIRS; PLEASE CONTACT THE ORIGINAL PLACE OF PURCHASE. ITEMS SHOULD USUALLY RETURNED TO THE ORIGINAL PLACE OF PURCHASE ALL YOUR COST.

THE WARRANTY GRANTED BY THE E&M ACTIVE® PARTNER IS LIMITED TO REPAIR OR REPLACEMENT OF THE PRODUCT. THE E&M ACTIVE® PARTNER UNDERTAKES WITHIN THE SPECIFIED PERIOD TO REPAIR OR REPLACE ANY PART OF THE APPLIANCE FREE OF CHARGE IF IT IS FOUND TO BE DEFECTIVE. THE WARRANTY DOES NOT COVER WEAR AND TEAR, USER ERROR, ACCIDENTAL DAMAGE, NATURAL DISASTERS, COMMERCIAL OR COMMUNAL USE, SECOND-HAND ACQUISITIONS, MISUSE OR REPAIR BY A PERSON OTHER THAN AN AUTHORISED SERVICE AGENT FOR THE DEALER. UNDER NO CIRCUMSTANCES SHALL THE APPLICATION OF THE WARRANTY GIVE RISE OF THE COMPLETE REPLACEMENT OF THE APPLIANCE OR ENTITLE THE CONSUMER TO DAMAGES.

## TROUBLE SHOOTING

TROUBLES	POSSIBLE REASON	SOLUTION
IT DOES NOT WORK.	THE SWITCH IS OFF.	TURN ON THE SWITCH.
	IT HAS RUN OUT OF POWER.	RECHARGE THE BATTERY.
	USB / MICRO USB IS NOT CONNECTED PROPERLY.	WHEN CHARGING, MAKE SURE YOU CHECK THAT THE USB/MICRO USB AND CHARGING COMPONENTS ARE PROPERLY PLUGGED IN AND CONNECTED.
	PRODUCT DAMAGE.	CONTACT PLACE OF PURCHASE.
THE PRESSURE IS SLOWING DOWN.	THE APPLIANCE IS RUNNING OUT OF POWER.	PLEASE CHARGE THE APPLIANCE AND USE AGAIN ONCE CHARGED.

# EM01 MAGICTOUCH OPTIC MASSAGER

FCC ID: 2AEHM-EM01 IC: 20096-EM01

E&M ACTIVE®

# EM01 MAGICTOUCH OPTIC MASSAGER

FCC ID: 2AEHM-EM01 IC: 20096-EM01

E&M ACTIVE®



## WARNING

### FCC

CAUTION: THE USER IS CAUTIONED THAT CHANGES OR MODIFICATIONS NOT EXPRESSLY APPROVED BY THE PARTY RESPONSIBLE FOR COMPLIANCE COULD VOID THE USER'S AUTHORITY TO OPERATE THE EQUIPMENT.

THIS DEVICE COMPLIES WITH PART 15 OF THE FCC RULES AND INDUSTRY CANADA LICENCE-EXEMPT RSS STANDARD(S). OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS:

(1) THIS DEVICE MAY NOT CAUSE HARMFUL INTERFERENCE, AND

(2) THIS DEVICE MUST ACCEPT ANY INTERFERENCE RECEIVED, INCLUDING INTERFERENCE THAT MAY CAUSE UNDESIRE OPERATION.

LE PRÉSENT APPAREIL EST CONFORME AUX CNR D'INDUSTRIE CANADA APPLICABLES AUX APPAREILS RADIO EXEMPTS DE LICENCE. L'EXPLOITATION EST AUTORISÉE AUX DEUX CONDITIONS SUIVANTES :

(1) L'APPAREIL NE DOIT PAS PRODUIRE DE BROUILLAGE, ET

(2) L'UTILISATEUR DE L'APPAREIL DOIT ACCEPTER TOUT BROUILLAGE RADIOÉLECTRIQUE SUBI, MÊME SI LE BROUILLAGE EST SUSCEPTIBLE D'EN COMPROMETTRE LE FONCTIONNEMENT.

NOTE: THIS EQUIPMENT HAS BEEN TESTED AND FOUND TO COMPLY WITH THE LIMITS FOR A CLASS B DIGITAL DEVICE, PURSUANT TO PART 15 OF THE FCC RULES. THESE LIMITS ARE DESIGNED TO PROVIDE REASONABLE PROTECTION AGAINST HARMFUL INTERFERENCE IN A RESIDENTIAL INSTALLATION.

THIS EQUIPMENT GENERATES, USES AND CAN RADIATE RADIO FREQUENCY ENERGY AND, IF NOT INSTALLED AND USED IN ACCORDANCE WITH THE INSTRUCTIONS, MAY CAUSE HARMFUL INTERFERENCE TO RADIO COMMUNICATIONS. HOWEVER, THERE IS NO GUARANTEE THAT INTERFERENCE WILL NOT OCCUR IN A PARTICULAR INSTALLATION. IF THIS EQUIPMENT DOES CAUSE HARMFUL INTERFERENCE TO RADIO OR TELEVISION RECEPTION, WHICH CAN BE DETERMINED BY TURNING THE EQUIPMENT OFF AND ON, THE USER IS ENCOURAGED TO TRY TO CORRECT THE INTERFERENCE BY ONE OR MORE OF THE FOLLOWING MEASURES:

-- REORIENT OR RELOCATE THE RECEIVING ANTENNA.

-- INCREASE THE SEPARATION BETWEEN THE EQUIPMENT AND RECEIVER.

-- CONNECT THE EQUIPMENT INTO AN OUTLET ON A CIRCUIT DIFFERENT FROM THAT TO WHICH THE RECEIVER IS CONNECTED.

– CONSULT THE DEALER OR AN EXPERIENCED RADIO/TV TECHNICIAN FOR HELP.

### FCC AND IC RADIATION EXPOSURE STATEMENT:

THIS EQUIPMENT COMPLIES WITH FCC AND CANADA RADIATION EXPOSURE LIMITS SET FORTH FOR AN UNCONTROLLED ENVIRONMENT.

THIS TRANSMITTER MUST NOT BE CO-LOCATED OR OPERATING IN CONJUNCTION WITH ANY OTHER ANTENNA OR TRANSMITTER.

### DÉCLARATION D'IC SUR L'EXPOSITION AUX RADIATIONS:

CET ÉQUIPEMENT EST CONFORME AUX LIMITES D'EXPOSITION AUX RADIATIONS DÉFINIES PAR LE CANADA POUR DES ENVIRONNEMENTS NON CONTRÔLÉS.

CET ÉMETTEUR NE DOIT PAS ÊTRE INSTALLÉ AU MÊME ENDROIT NI UTILISÉ AVEC UNE AUTRE ANTENNE OU UN AUTRE ÉMETTEUR.