

*make*ON

[www. makeonbeauty.co.kr](http://www.makeonbeauty.co.kr)

Product Name: Skin Light Therapy II Weight: 66g → 73g (based on design mock-up) Seller:

Amore Pacific Co., Ltd. Manufacturer: Amosense Co., Ltd.

100 Hangang-daero, Yongsan-gu, Seoul 90, 4sandan 5-gil, Jiksan-eup, Seobuk-gu, Cheonan-si,

Chungcheongnam-do

MADE IN KOREA 

Technology

Skin care device that cares from skin diagnosis to solution through light and microcurrent

LIGHT

Energy that light of a specific wavelength transmits into the skin



MICRO CURRENT

Activation of skin energy + promotion of product absorption



3min

Moisture supplementation care for dry skin

3min

Make-up care for well-applied makeup

3min

Elasticity up care that keeps the skin firm

3min

1min

1min

1min

Hypoallergenic care optimized for sensitive skin

Good to know before using!

1

Please recharge once after opening.

3 Page

2

Please select 'Recommended Mode vs Selected


4 page


Dry thoroughly with a towel and measure, do not apply skincare products for accurate results.


Featured Mode

- ① Turn on the power.
- ② Ready when the color of the button blinks
- ③ A beep sounds at the center of the ball. Measurement is complete after a long time

Measurement results

lack of moisture →  (Blink)

In case of lack of oil →  (Blink)

In case of lack of elasticity →  Blink

Use Featured Mode (3 minutes)

the sensor is placed on the face Automatically run the mode when touched

selection mode

- ▷ Use after selecting a mode by pressing the button



Sweating or prolonged exposure to the sun may not produce the right results.

*Only the beep is repeated if the measurement is not successful. Wait until the sound stops without removing it from the measurement site.

Product Introduction(

Makeon Skin Light Therapy III


A skin care device that provides a solution for skin change that is renewed by finding out the lack of skin

■ Touch sensor that recognizes skin

- Works only when in contact with non-dry skin → Use after applying a sufficient amount of cream formulation.

■ One button for easy use

- Press the button once to proceed in the following order: mode is changed.

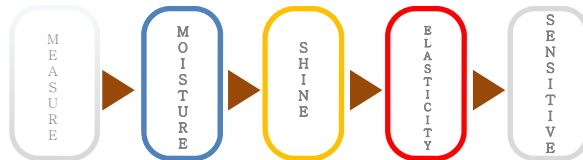

■ A sensor that measures the condition of the skin on your face

- If you measure the center of the ball for 3 seconds in the measurement mode, The most necessary mode is automatically If you want to use a mode other than the recommended mode Press the button to change to the desired mode and use it.


■ Long lasting use on a single charge

- 60 minutes of use with a full charge within 3 hours (can be used 20 times)
- After opening, the power does not turn on because the battery is low. Please use it after charging.

How to use


Select

1 Start

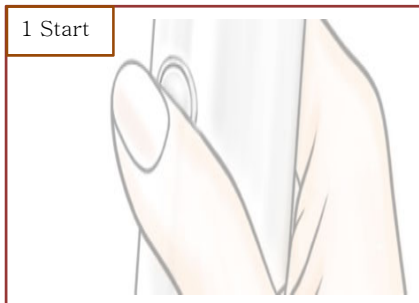


until you hearPress the button to turn on the power..

2 Measure

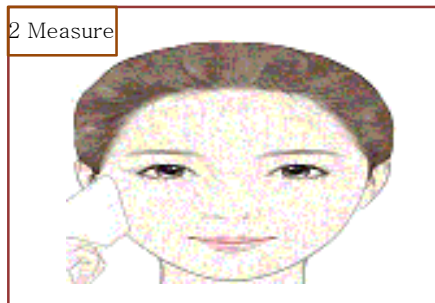

Point the device at the center of the cheek and hear a beep.wait until the end

3 Check

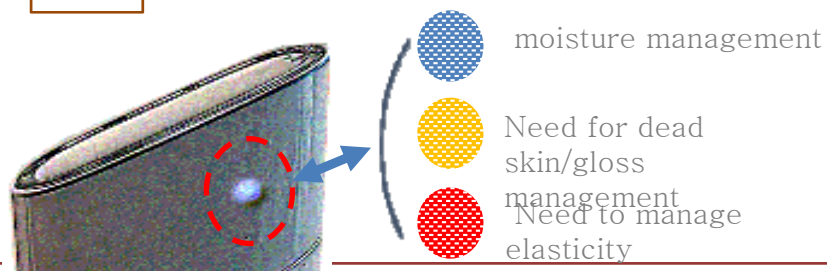


Check the required mode by the color of the side or button.

4 Care



After sufficiently applying the product to enhance the effect, both cheeks (1 minute),Take care of the rest of the area, such as around the eyes, forehead (1 minute),

5 Result

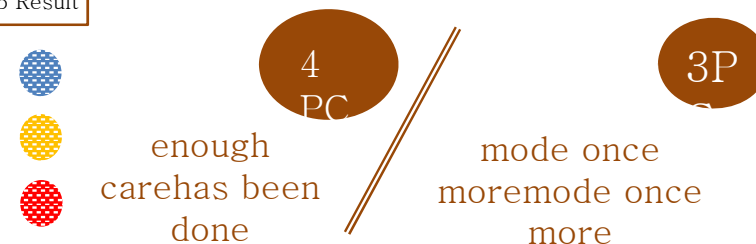

** Measurement results are visible only when used for 3 minutes in measurement mode*

After 3 minutes of use, check the effect with the side LED.

If the measurement is not performed properly, the alarm will continue to sound, so please wait until the measurement is finished.


How to use

select

After the skin care step of toner and lotion, use after selecting the desired mode while applying enough cream to the face.



2. Intensive care for the area to enhance the effect, and finish with essence for moisturizing/soothing after use.

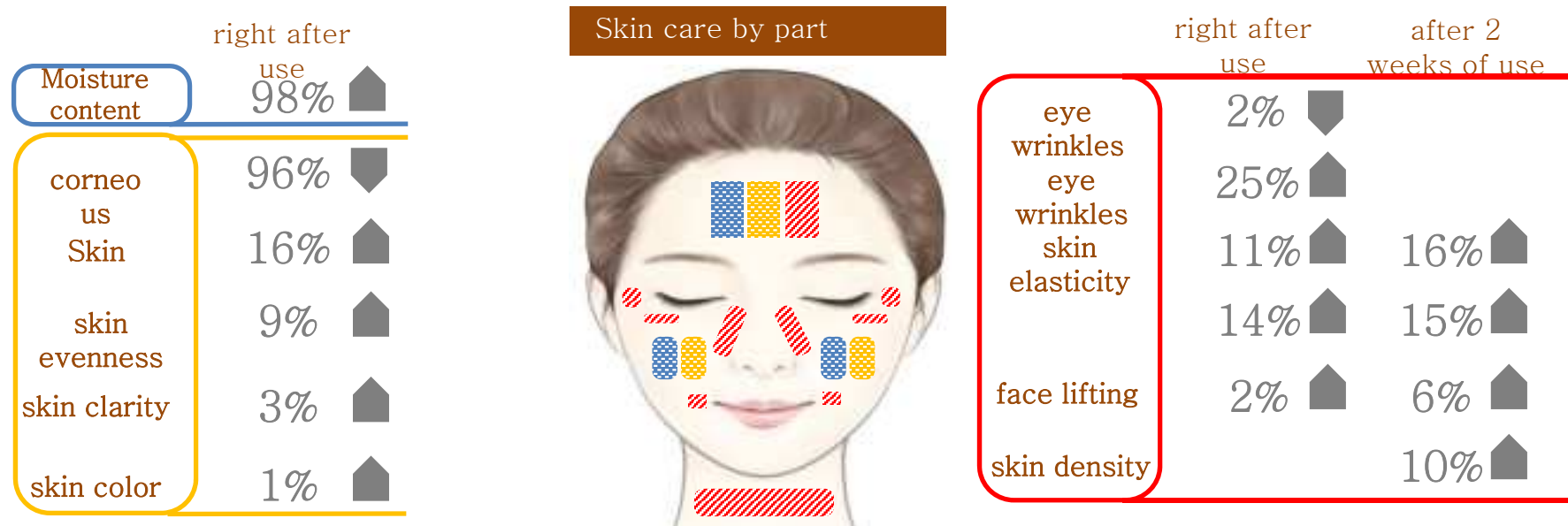

Questions about use


Can the measurement mode be affected by the season?

- Spring/Autumn: Due to dry weather, the frequency of moisture and oil care is highly recommended.
- Summer: Due to the intense heat, you sweat a lot, so firming care may be intensively recommended..
- Winter: Due to the cold and bitter wind, moisture and oil care are highly recommended.

Can't you check the measurement result with the LED on the side after using the mode selection?

- After the normal measurement in the measurement mode and the use of the automatically set mode is completed, the side LED You can see how much better it is. (shows up to 4)

Having a hard time checking if the use mode is working?

- You can check which mode is active by checking the color of the button..
((Example: LED color above the button  → BLUE mode)
- Use it so that it comes in contact with your skin, and if you feel the light and vibration, you are using it correctly.

Questions about use

Do you feel a stinging sensation on your skin?

- Safety was verified by conducting a skin safety evaluation test for skin irritation¹
- If you feel uncomfortable with the stinging sensation, you can get a good effect when using it with a mask sheet attached.

Do you feel the light is too bright?

- Completed certification of 'photobiological safety test' for eyes
- However, be careful not to look directly at the light for a long time.
- When using around the eyes, close your eyes to minimize glare. If you feel any discomfort around your eyes, such as eye fatigue, pain, foreign body sensation, etc., please stop using it temporarily.

Precautions when using

when using

If you press it hard on your face, it may irritate your skin. Please rub lightly.
Do not use on areas with severe redness or swelling due to skin irritation during use. Please discontinue use until it subsides. • Do not use on infected, burned or cut skin as skin damage may occur. • If you have a wound on your skin, use it after treatment..

when storing

• Do not use if pregnant or nursing mothers.
Do not drop the product or subject it to external impact. • Do not allow children or pets to put the product in their mouth or bite it. • Store at room temperature away from direct sunlight, and do not store at too high or low temperatures. • When maintaining the device, wipe lightly with a soft cloth once a week. • Do not leave in direct sunlight for a long time. It may cause discoloration.

when charging

• Using a charger that is not suitable for the product may shorten the life of the product or cause a malfunction. • Do not excessively bend the charger cord or use the product while charging..

• The charging adapter is not designed to be waterproof, so do not immerse it in water or charge it in a wet place.

do not separate

Free repair is not available if the user arbitrarily modifies or separates the product.

problem solving

I wonder	I'll give you an answer
Is it okay to use every day?	<ul style="list-style-type: none"> It is designed to be gentle on the skin even if used daily, but it may cause severe stinging depending on the skin type. In the case of sensitive skin, use only 1 set (3 minutes) after applying enough skin care products. If skin trouble occurs after use, stop using it for a while before using it.
A tingling sensation during use Are you feeling uncomfortable?	<ul style="list-style-type: none"> Depending on your skin type or condition, there may be slight irritation. When using for the first time, it is recommended to use MODE 1 (BLUE), which has the weakest microcurrent, for 3 minutes. Gradually increase the usage time step by step, and apply enough skin care products such as cream or massage gel before care.
LED wavelength during use Can't you see?	<ul style="list-style-type: none"> Negative ions and vibration can be checked only when the proximity sensor touches, and the proximity sensor can be checked with a gray dot on the top. The sensor does not work on dry skin, so please use it after applying a skin care product..
Even in combination with dermatological treatment it's okay	<ul style="list-style-type: none"> Since it is a beauty device, not a medical device, it is recommended to use it according to the doctor's opinion related to skin disease. Use it after determining whether to use it or not.
wearing jewelry Is it available??	<ul style="list-style-type: none"> It is possible to use it because there is no direct effect on surfaces that come into contact with the skin, such as earrings and necklaces. Skin care products used after cleansing are recommended to be used without precious metals.
Do you hear a 'beep beep' sound while using it?	<ul style="list-style-type: none"> It is a notification sound that guides you to use it by moving it according to the time easily for each part. A beep sounds when one minute elapses after the power is turned on. Both cheeks (1 minute), eyes & T zone (1 minute), It is convenient to use in the order of forehead, chin & philtrum (1 minute).

problem solving

I wonder	I'll give you an answer	
A “beep beep beep beep” sound from the main unitDo you see a red notification light blinking?	This feature notifies you when charging is required. After connecting the main body and the charging cable After charging is complete When the white full charge notification light on the main unit flashes, charging is complete. The time required to fully charge may vary depending on the charging environment.	
Even if you press the power button while chargingNot working?	For safety in use, the main body is designed not to operate while charging is in progress. After charging is complete, separate the charging cable from the main body before using it.	
Even if you connect the charging adapterNot charging??	Please check if the product is properly connected	Make sure the product is well connected with the cable. Please check that the outlet and adapter are disconnected. Please check once more that it is not
	Check that the outlets in your home are working properly.	lease connect to another outlet
Does the product not turn off?	It could be that you let go of the power button too quickly.	Press and hold the power button for 1-2 seconds then release
Product won't turn on?	Please use it after charging whether the product is charged	Please use it again after charging for more than 30 minutes.

quality assurance

product warranty

product name	Skin Light Therapy III		
warranty period	1 year from date of purchase		
	Y	M	D
Purchase date	※ Be sure to indicate the date of purchase If not listed, it may be invalid		
where to buy			

In the condition of normal use according to the precautions such as the user manual within the warranty period In case of a breakdown, we provide free A/S according to the warranty regulation, and the paid service of the warranty regulation If applicable, A/S will be provided for a fee. If you want A/S, the warranty certificate and warranty regulations After checking the contents, contact the AMOREPACIFIC Customer Service Center
Customer service center (receiver pays): 080-023-5454 If there is something wrong with this product, we will compensate you in accordance with the Consumer Dispute Resolution Standards announced by the Fair Trade Commission.

Warranty provisions

Free Service

The warranty is for the main body of the device.

Free service is available for malfunctions that occur during normal use within one year (warranty period) after purchase.

Even during the warranty period, the following cases will be charged for service.

1. In the event of a failure due to careless handling by the consumer (dropping, submersion, breakage, or damage after purchase)
2. Changes in appearance (scratches, deformation, discoloration, etc.) that occur during normal use
3. 'In case of failure due to not observing the 'precautions'
4. In case the consumer arbitrarily disassembles or modifies
5. If this warranty is not presented or the date of purchase is not indicated in this warranty

FCC Compliance Statement

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to correct the interference by one of the following measures

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter

FCC Radiation Exposure Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 5 mm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

FCC Compliance Statement

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation

Product Specification

Frequency Range	2400 MHz to 2483.5 MHz
Charger	USB-C
Operating Voltage	DC 3.7 V (Battery)
Temperature Range	-20 °C to 50 °C
Wireless	BLE(1 Mbps) / 40 ch
H/W Version	1.0
S/W Version	1.0