


Step 6: Review User Guide

The User Guide provides key information about Libris, including Important Safety Information. Please review it carefully and keep it nearby for reference.

If you have questions, call your response team.

By activating your Numera Libris you agree to the Terms of Use and understand the Important Safety Information in the Numera Libris User Guide.

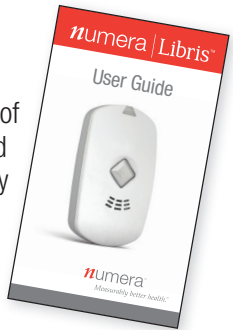

Quick Reference

Turn Libris ON: Press the Call Button and hold until the Battery Indicator flashes.

Turn Libris OFF: Press the Call Button twice quickly; hold the second press for 7 seconds or until Libris plays a message confirming your request to turn off. Press the Call Button again to confirm.

Place an emergency call:

1. Press and hold the Call Button.
2. Listen for the alert tone, followed by the audio message: “Calling your response team now. Please hold the device and position the microphone near your mouth to complete your call”.
3. Your response team will answer and notify emergency services for you.

IMPORTANT:
Pendant Calling Devices will only operate within a limited distance from the Charger. Maximum recommended distance is 150 feet from the Charger.

Battery & Cellular Indicators

Battery Indicator around Call Button shows the level of battery charge.

Battery Indicator	Status*	Action Required
Green	Good	Ready for use
Amber	Moderate	No action needed
Red	Low	Charge soon

*Battery is designed to hold a charge for 36 hours but may be reduced as a result of user activity, cellular coverage, talk time, device-specific settings, and device life.

Cellular Indicator on top of Libris shows cellular signal strength.

Cellular Indicator	Signal Strength*
Green	Strong
Amber	Moderate
Red	No Signal

*A cellular signal may not be available in all areas.

About Libris



1. **Cellular Indicator** shows cellular signal strength.
2. **Battery Indicator** shows battery status.
3. **Call Button** initiates a two-way call to response team.
4. **Program Indicator** shows when optional Pendant devices are added.
5. **Power Connector** for power cord.
6. **Program Button** used to add optional Pendant devices.

numera | Libris™

Quick Start Guide



numera™
Measurably better health.™

REGULATORY INFORMATION

We, Nortek Security & Control, LLC of 1950 Camino Vida Roble STE 150, Carlsbad, CA 92008, declare under our sole responsibility that this device complies with Part 15 of the FCC rules.

FCC & IC Notice

This device complies with Part 15 of the FCC Rules and Industry Canada license exempt standard(s). Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference received that may cause undesired operation. Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

This Class B digital apparatus complies with Canadian ICES-003. Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician to help.

WARNING: Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate the equipment.

To satisfy FCC/IC RF exposure safety requirements, a separation distance of 20 cm or more should be maintained between this device and person's body (excluding extremities: hands, wrists, feet and ankles).

S30-xxxxx



For more information about Libris visit

www.numera.com

Numera, Numera Libris, "Measurably Better Health" are trademarks of Nortek Security & Control LLC. All other company and product names used herein may be trademarks, or registered trademarks of the respective companies.
Copyright © 2016 Nortek Security & Control LLC 10009954 X2

Enjoy new freedom

This Quick Start Guide will help you set up and activate Libris. **Please read this entire guide before using Libris to ensure it is set up correctly and ready for use.**

What's in the Box



Step 1: Activate Libris

1. Insert the Power Cord into the Charger. Plug the large end of the Power Cord into a standard wall outlet.
2. Place Libris in the Charger.
3. Libris will automatically turn on and begin activation.
4. Follow the voice prompts until directed to connect to your response team.



Step 2: Call to Set Up

1. Press and hold the Call Button.
2. Your response team will answer shortly to customize Libris and complete set up.
3. After the call, charge Libris in the Charger for 4 hours. Then it is ready for use.

IMPORTANT:

- Libris is **not ready for use** until you have:
- Charged Libris in the Charger for 4 hours
 - Completed a call with the response team




Step 3: Adding Pendants

The Libris Charger can be programmed to recognize optional Pendant Calling Devices.

To add a Pendant Calling Device:

1. Be sure the Libris is in the Charger and the Charger is powered on.
2. Press and release the Program button on the rear of the Charger. The Program indicator will light amber.
3. Within 20 seconds, activate a Pendant Calling Device.
4. The Program indicator should light green when the Pendant is recognized.
5. Repeat steps for additional Pendants.

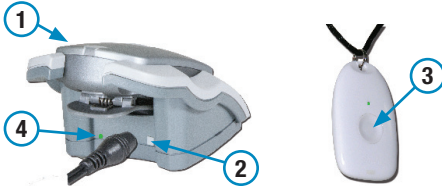

Step 4: Testing Pendants

If optional Pendant Calling Devices are added, they must be tested.

To test Pendant Calling Devices:

1. Be sure the Libris is in the Charger and the Charger is powered on.
2. Activate a Pendant Calling Device by pressing the button on the Pendant.
3. Your response team will answer shortly to verify your testing.

IMPORTANT:

Pendant Calling Devices can only place an emergency call if the Libris is resting in the Charger and the Charger Power Cord is plugged into AC power.

Pendant Calling Devices do not have a speaker or microphone, talk to the response team through the Libris unit.

Step 5: Wearing Libris

Libris can be worn on a belt using the Belt Clip, or as a pendant around your neck using the Lanyard/Clip. If worn with the Lanyard, Libris should be worn under clothing to prevent excess movement of the device.

To attach the Belt Clip or Lanyard/Clip:

1. Align the notched side of Clip with the notched area on the back.
2. Press inward until the Clip notches are nested within the Libris notches.
3. Firmly push the Clip toward the top of Libris. It is secure when you hear a click or you can no longer see the notches on the Clip.

