



Occean Smart Sonic  
Electric Toothbrush

User Manual

Thank you for using Occean SE  
Smart Sonic Electric Toothbrush.

### Important Safety Tips

- Due to the high-speed vibrations, it is normal to experience a tingling sensation while using the toothbrush for the first time. This feeling will gradually go away by regularly using the toothbrush.
- If the toothbrush is not used for a long period of time, the battery will run out due to the power consumption from standby mode. This can be solved by charging, the charging indicator on the brush will not light up instantly. Please wait for 6 - 10 minutes for the toothbrush to go back to its normal state.
- It is normal for the toothbrush to heat up when it is charging. If the toothbrush becomes extremely hot, immediately stop charging and/or turn it off and take it to Occean's authorized customer service for checking.
- The toothbrush can only be used with the Occean special head.
- The toothbrush is for oral use only. Do not use it for other purposes.
- The toothbrush is a personal item. The product is not intended for use in a dental clinic and should not be used for multiple patients.
- The toothbrush should not be placed under direct sunlight for a long period of time. It should be stored in a cool and dry environment.

### People with Sensitive Oral Cavity

The following people should use this product with great caution.


People who have undergone oral peridental surgery or are being treated.

People who have mouth diseases, loose false teeth or gums with haemorrhage.

People with sensitive and susceptible oral cavity, they should consult a professional dentist before using this product.

The toothbrush is not suitable for people with reduced mental and physical capabilities. Not suitable for infant and children under the age of 12. Children who use this product should be strictly supervised by an adult.

Please regularly clean the power plug to prevent the accumulation of dust and dirt which can reduce the charging effect.

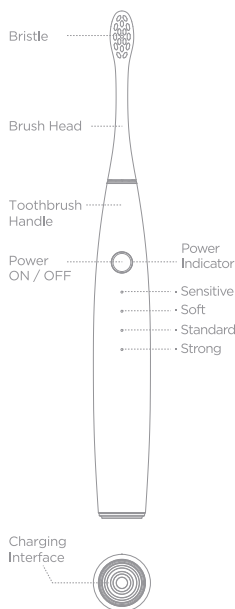

### Instructions

#### • Direct Brushing Mode(Default)

Can be enabled through the app. Short press on the Power button will skip the Awake mode and it will directly start the cleaning process. Press the Power button for 2 seconds to go to the next level of cleaning while brushing.

#### • Anti-splash (App Enabled)

Turn on the anti-splash function through the Occean App, which effectively prevents liquid from splashing when the toothbrush is not operating in direct contact with the tooth.

#### • Awake Mode (App Enabled)

Shortly press the Power button to turn on the toothbrush and it will play the welcome sound. It will automatically enter into dormant state if not operated for more than 15 seconds.

Shortly pressing the toothbrush after turning on the power will make the toothbrush carry out a specified cleaning plan.

Following the cleaning plan, the toothbrush will automatically stop working after a specified time. If you want to stop brushing before reaching the recommended time, short press the Power button to stop brushing.

Press the Power button for 2 seconds to change the cleaning intensity while brushing.

### 4-S Level Cleaning

- L1 ( · ) ... 1 Light On ... Sensitive
- L2 ( · · ) ... 2 Lights On ... Strong
- L3 ( · · · ) ... 3 Lights On ... Standard
- L4 ( · · · · ) ... 4 Lights On ... Soft

#### **Sensitive**

Suitable for first time electric toothbrush user with sensitive teeth.

#### **Soft**

Suitable for first time electric toothbrush user.

### Standard

Suitable for experienced electric toothbrush user.

#### **Strong**

Suitable for user with high requirement cleaning process.



### Smart Tooth Brushing Solution

By learning from many experiments and thanks to suggestions from dentists, the Occean SE is an innovative product in the industry. Our new brand of smart and sonic toothbrushes meet the need of every different individual. Download and log into our App for more information.

### Recommended Tooth Brushing Method

#### • Tooth Brushing Time

Oral experts (Dentists) recommended that each brushing should be 2 to 3 minutes long.

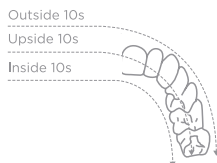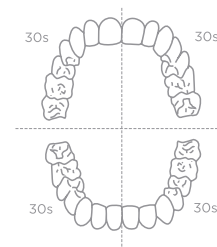

### Steps



Insert the brush head carefully into the toothbrush. A small gap should be left between the head and the main body for work space so the brush head can work normally.



To prevent toothpaste from splashing, smear a small amount of toothpaste on the head before using it. (A sonic electric toothbrush saves more toothpaste than a manual toothbrush.)



Soak the toothbrush head and put it into your mouth, then turn on the toothbrush to start brushing.



Slowly move the bristles at 45° to the gum line to brush teeth properly.



Hold the toothbrush handle and slowly move it in circular motion. Do not press too hard.



Please rinse with clean water after brushing.

### Specifications

Name	Occean Smart Sonic Electric Toothbrush
Type	Occean SE
Magnetic Power	270 gf.cm
Power of Subject	1.8 W
Charger	Occean USB Cradle
Input	5V / 1A
Service Time	60 Days ( * )
Battery Capacity	2600 mAh ( ** ) Li-ion
Charging Time	-6 Hours ( *** )
Bluetooth	Bluetooth 4.0 above


\* Validated by 2 minutes brushing period, twice a day.

\*\* Theoretically, battery has capacity range from 2,250mAh to 2,950mAh.

\*\*\* It is suggested to use a USB charger of 5V/1A or 5V/2A.

### Magnetic Charging

When the OcClean SE touches the special charging cradle, the bottom of the toothbrush will light up, which shows the toothbrush is being charged. In case of incorrect charging (such as short circuit and overcurrent), the warning LED will light up.

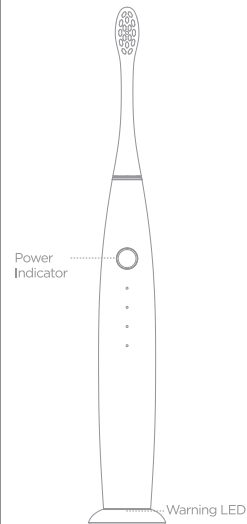

### Charging Process

0% - 95%, the Power Indicator light flashes.

95% - 100%, the Power Indicator light is constantly on.

### Use Process

If the toothbrush battery is at less than 5%, the Power Indicator will flicker in red as a reminder to charge the toothbrush.

### Notes

The battery can't be replaced by users themselves and can only be replaced by the staff in OcClean Authorized Client Center. Unauthorized bodies are forbidden to disassemble this product.

The USB cable can't be replaced by unauthorized body. In case of damage, please contact the OcClean After-sale Center for maintenance.

### Cleaning and Maintenance



#### Cleaning Brush Head

Clean the brush head after use and air dry it. Do not pull the brush head for this act as it may result in deformation of the toothbrush.



#### Cleaning Body

Use clean water to wash the surface of handle and scrub it with soft cloth and air dry it. Do not press the rubber sealing ring with sharp item or metal object, otherwise it will be damaged.



#### Cleaning Charger

Scrub the dirty parts of the charger with a soft cloth. Pull out the battery plug from the socket prior to cleaning the charger. Scrub the surface of charger with wet cloth.



#### Change the Head

It is suggested by dentists that the head should be changed every 2-3 months to maintain optimal brushing.



### Remark

X The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;

O The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;

Toxic and harmful substance or element of parts	Environment Protection				
	Pb	Hg	Cd	Cr6+	PBB
Vibrating Electroacoustic device	X	O	O	O	O
Accessories	X	O	X	O	O
Battery	X	O	O	O	O

Name of parts	Environment Protection				
	Pb	Hg	Cd	Cr6+	PBB
PCBA	X	O	O	O	O
Plastic parts	X	O	O	O	O
Metal Hardware	X	O	O	O	O
FPC	X	O	O	O	O

### Guarantee Certificate of OcClean Smart Sonic Electric Toothbrush

#### User Information

Name: \_\_\_\_\_

Tel: \_\_\_\_\_

E-Mail: \_\_\_\_\_

Address: \_\_\_\_\_

Sales information

Retailer: \_\_\_\_\_

Invoice No: \_\_\_\_\_

Tel: \_\_\_\_\_

Sales Date: \_\_\_\_\_

Production Place: \_\_\_\_\_

SN

Type: \_\_\_\_\_

Product Information

#### Warranty Card

Delivery Date

Fault description/Reason/Disposal

Employed equipment

Acceptance date

Serviceman

Please ask for a maintenance report from the maintenance center and properly keep it with the warranty certificate. Shenzhen Yunding Information Technology Co., Ltd. in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers guarantee for customers who purchase OcClean Smart Sonic Electric Toothbrush. In case of inconsistency between the service articles and the provisions of guarantee, the provisions of guarantee and the contents in the guarantee certificate shall govern.

Shenzhen Yunding Information Technology Co., Ltd.

Floor 28, Tower C, Dachong Commercial Center, No. 9676 Road, Shenzhen, China  
Zip Code: 518057

If there is fault occurring in the product you purchase, please contact email: service@occlean.com

### OcClean Facebook Page



[www.facebook.com/occleanofficial](http://www.facebook.com/occleanofficial)

Follow the  
OcClean Official Facebook Page

### Download App



[www.occlean.com/download](http://www.occlean.com/download)



[www.occlean.com](http://www.occlean.com)

### FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: -- Reorient or relocate the receiving antenna.

-- Increase the separation between the equipment and receiver. -- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. -- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID: 2AN5D-OCLEANSE