




TUBELET
Water Flosser

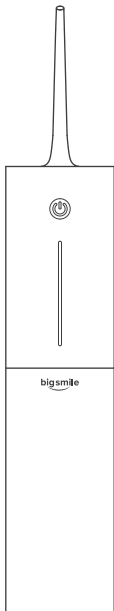

Instruction Manual

For your safety, please read
the instruction book carefully before using.

FCC STATEMENT

1. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference.

(2) This device must accept any interference received, including interference that may cause undesired operation.

2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- * Reorient or relocate the receiving antenna.
- * Increase the separation between the equipment and receiver.
- * Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

* Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. This equipment should be installed and operated with a minimum distance of 20 cm between the radiator and your body.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be taken, including the following.

CAUTION: READ ALL INSTRUCTIONS BEFORE USE




To reduce the risk of flares, electrocution, fire, or injury:

1. The cable provided fits into a 5V DC USB Type A charger.
2. Never operate the appliance if it has a damaged power cord, if it is not working properly, if it has been dropped or damaged or has been submerged in liquid for an extended time.
3. Keep the power cord away from hot surfaces.
4. Close supervision is necessary when this product is used by, or near children or individuals with disability.
5. Use this product only for its intended use as described in this manual or recommended by your dental professional.
6. Do not directly point the water jet under the tongue, into the ear, nose or other delicate areas. This product operates with high pressure that can damage sensitive areas. See For Your Safety for proper usage.
7. Only use attachments recommended by the manufacturer.
8. Never operate where oxygen or aerosol sprays are being used.
9. The appliance can not be operated while charging. Wait until it is fully charged and separated from the charger before use.

WARNING - FOR YOUR SAFETY

1. If you had oral or periodontal surgery within 2 months, contact your dentist before use.
2. When starting any new oral hygiene program, your gums may have some bleeding. Contact your dentist if excessive bleeding occurs or if the bleeding continues after two weeks of use.
3. Do not use mouthwash or solutions that contain iodine, bleach or tea tree oil.
4. Do not use while wearing dental jewelry.
5. Do not use if you have any open wound or sores on your tongue or in your mouth.
6. If you are taking any antibiotics or other medications, consult your dentist before using this oral hygiene aid.
7. Closely supervise children and individuals with special needs to ensure proper use, and prevent them from any injury from this appliance.
8. Be sure that only the water tip touches your teeth, gums or top of the tongue to avoid any injury or damage to those areas.
9. Use only for oral hygiene cleaning. Do not use for any other purpose.
10. As with any new oral hygiene program, individuals with heart problem should check with their physician prior to using the water flosser.

PRODUCT ITEMS AND INTRODUCTION

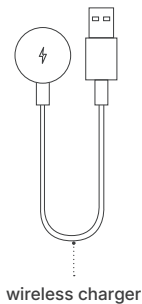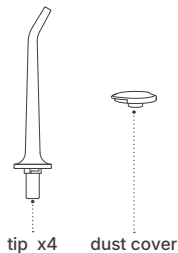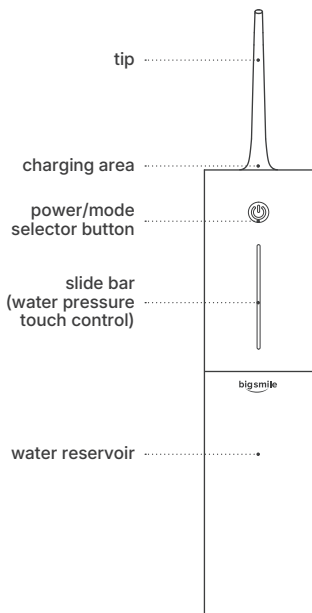


CHARGING

Note: The recommended safety charging voltage is 5V DC. Before charging, make sure the output voltage of the power adapter matches the recommended power rating. A higher voltage may cause damage to the water flosser.

1. Fully charge the water flosser for 4-5 hours using the included wireless charger. **The lights will flash rapidly to indicate low battery and as a reminder to recharge it.
2. Use the provided wireless charger with a USB Type A power adapter to charge the water flosser (the charging voltage is 5V DC).
3. Place the charger gently on top of the appliance until the magnets are aligned and affixed for charging.
4. The light on the SlideBar will stay fully lit when the appliance is at full charge.




READ ALL INSTRUCTIONS BEFORE USE

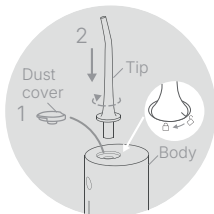

Attaching the Tip

1. Remove the dust cover.
2. Insert the tip into the water outlet hole, align the target line on the "skirt" of the tip to the unlock icon □ on the top of the body, then rotate clockwise to secure it aligned with the lock icon 🔒.

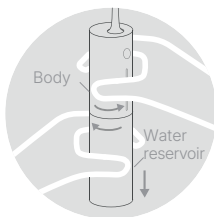


Remove the Water Reservoir

Counter twist the water reservoir and the body then pull apart to open.

-NOTE: Align the arrow with the printed icons from lock 🔒 to unlock □.


* When new, some material scent may be present. It will fade after several uses.

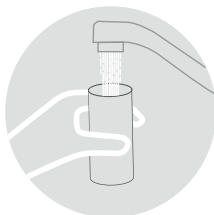

Fill with Water

Fill the reservoir with clean water.


* For longer lifespan of the device, it is best to avoid using mineral rich water that can build up hard deposit in the pump.

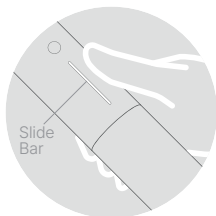

Start the Device

Quickly push the ON/Off button or hold the device to activate it.

Touch and slide down the Slide Bar to adjust to the lowest pressure setting.

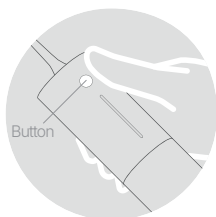


Mode Selection

(While on Standby)

Double Click to cycle Strong/Gentle Press and Hold to cycle Strong/Pulse [Strong/-Pulse can switch while in use] then Quick Click: ON/OFF

* To prevent accidental intensity changes, the Slide Bar will automatically lock, if no adjustment within 3 seconds of power-up. If you need to adjust the intensity, please press the button to unlock the slide-bar at any time.

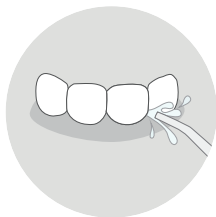


Start Cleaning

Place the tip inside the mouth pointed towards a tooth or a gap.


Keep mouth slightly closed then press the ON/Off button to start.

Slowly move the tip around and clean between each teeth.

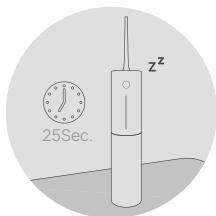

After Cleaning

Press the ON/Off button to stop the pump. Set the water flosser on a stable flat surface. It will go into Sleep Mode after 25 seconds.

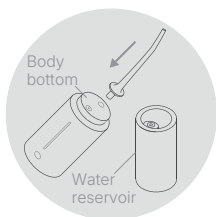

Storing the Tip

Open the water reservoir. Insert the tip into the hole at the bottom of the device body. Close the water reservoir.



Storing the Dust cover

1. Open the water reservoir.
2. Insert the dust cover into the hole at the bottom of the device body.
3. Close the water reservoir.

IMPORTANT NOTES

- * You may experience slight discomfort or bleeding when you use the product for the first time, which will likely disappear in 1 to 2 weeks. If you experience excessive bleeding or bleeding continues to occur after 2 weeks, please contact your dentist.
- * When the power is low, the status light will blink twice in rapid succession, and will repeat every 5 seconds. Please follow the charging instruction to safely recharge the water flosser.

MAINTENANCE

- * Use only low mineral (distilled is recommended) water and mild detergent to wash.
- * Do not wash in a dishwasher. or with water over 40 ° C strong chemicals such as bleach or hydrogen peroxide to clean the product.
- * Remove the tip and clean it with running water and wipe with damp cloth after each use. Do not bend, stretch or twist the tip.
- * Detach the water reservoir and clean it with running water and gently brush to get rid of any residue.
- * If the product is not being used for an extended period, please store the product in a cool and dry place away from direct sunlight.
- * For best results, it is recommended to change your tip every 4 to 6 months.

PRODUCT SPECIFICATION

Model No. : PB-01

Net Weight : 265g


Size : 219x53x53mm

Adapter : DC 5V $\overline{=}$ 1A

IP Level : IPX7

VISIT OUR WEBSITE

www.bigsmilecare.com

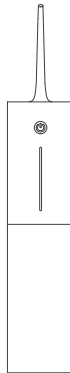

bigsmile

WATER FLOSSER