


**FCC ID: 2BOVE-HK-712**

1. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.
2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.
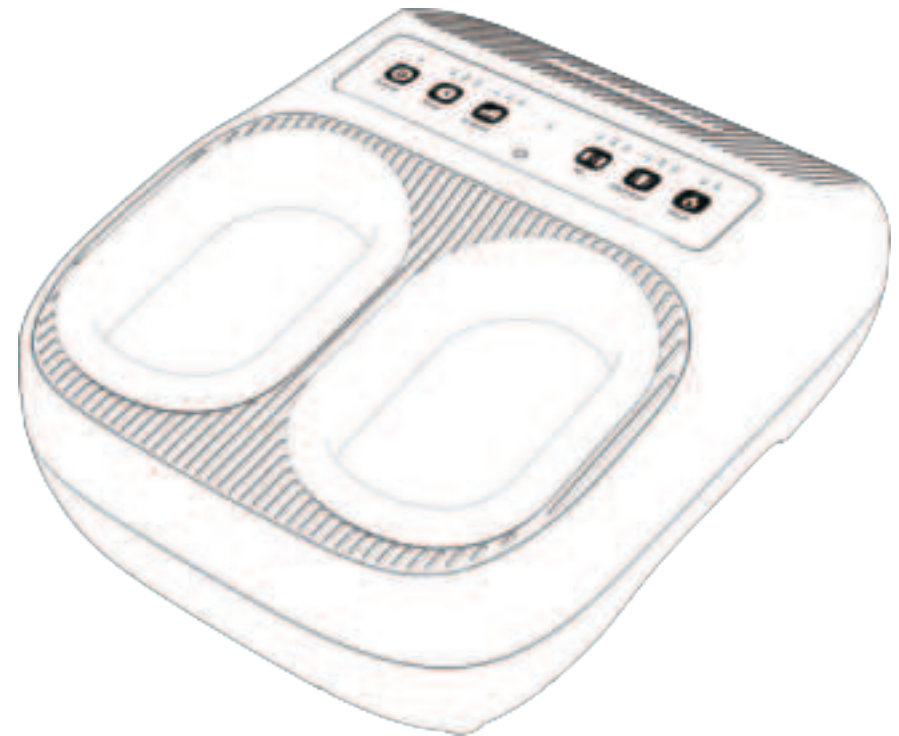
## SHIATSU FOOT MASSAGER USER MANUAL



○ BLACK

○ WHITE

○ PINK





Please read this manual before use, if have any issues, contact us by Amazon E-mail or **veicomtech@hotmail.com**

## USE CONTRAINDICATIONS:

1. Check that the voltage matches the rated specifications on the massager. Don't unplug the power supply with wet hands and not leave the massager unattended after it has been plugged into the power supply.
2. People with cardiac pacemaker, intolerant to electric impulse, with venous thrombosis, tetanus, local metallic foreign body or bleeding tendency must consult their health care practitioner before using this appliance.
3. Pregnant women should consult their health care practitioner before using this appliance.
4. If you suffer from any problems with your legs or feet such as varicose veins, please consult your health care practitioner before using this appliance.
5. Stop using this appliance immediately if you feel discomfort or any pain and irritations occurs.
6. Avoid using the massager when your feet are wet
7. Do not sit or step hard on the massager.
8. Do not use this appliance consecutively for over 30 minutes.
9. Do not use the appliance while sleeping or if you are feeling drowsy in case of over use for safety concerns.
10. Always turn the power off before you insert or remove a plug. Please remove by grasping the plug and do not pull on the cord.
11. Make sure that the power applied on this appliance matches the rated power.
12. Close supervision is necessary when this appliance is used by children. Children under 8 years old are prohibited to use this appliance.
13. This appliance should be used for its intended purpose only. Users should be responsible for consequences caused by improper operation on this appliance.

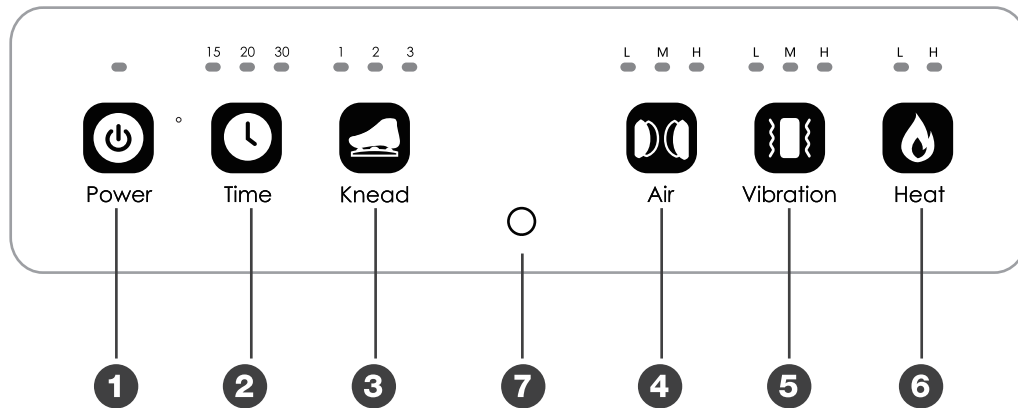
## 4. FUNCTION SELECTION

Customize your massage experience using these functions:

Function	Operation	Note
Knead		Activates combined kneading + air compression
Air		Air compression only When knead function is on, air can't be turned off alone
Vibration		Leg vibration only
Heat		Heating function

**Note: All functions can operate simultaneously. Press corresponding buttons repeatedly to adjust intensity levels.**


## INSTRUCTION:



- 1** Power: Start/ Stop the massager
- 2** Time: Set operating time
- 3** Knead : Activate/ Deactivate the kneading and air function
- 4** Air : Activate/ Deactivate the air compression function
- 5** Vibration: Activate/ Deactivate the vibration function
- 6** Heat: Activate/ Deactivate heat function
- 7** Infrared receiver

- 8** One-touch start four functions simultaneously
- 9** Lock enabled functions and the function keys can't be adjusted

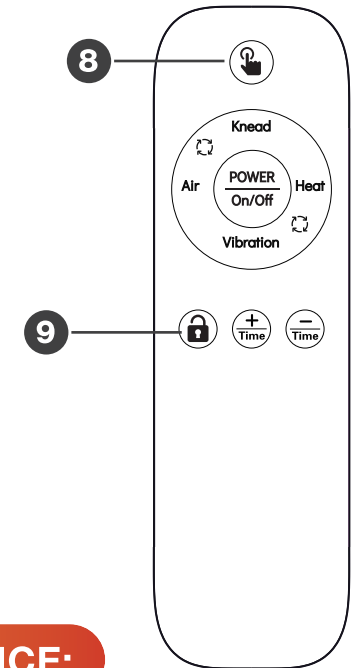Note: The receiving range of the remote control is limited to 1 Meter. Please use it within this range.




## CLEANING AND MAINTENANCE:

1. Before cleaning the massager, should be unplugged from the power cord
2. Do not use benzene, thinner or other corrosive solvents to clean your massager
3. Do not use a wet cloth to clean the massager, the power switch or alcohol. You must use a dry cotton cloth
4. If you do not use it for a long time, it is recommended that you roll up the power cord of the massager, and the massager is stored in a dry and clean place
5. Do not store the instrument in high temperature, open fire or direct sunlight
6. The use of the instrument should be gently held and placed, do not fall, pressure, impact or strong impact!

## FOOT MASSAGER OPERATION GUIDE


### 1.Initial Setup

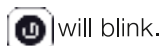1.Connect the foot massager to the power supply.




2.Press  to turn on the massager

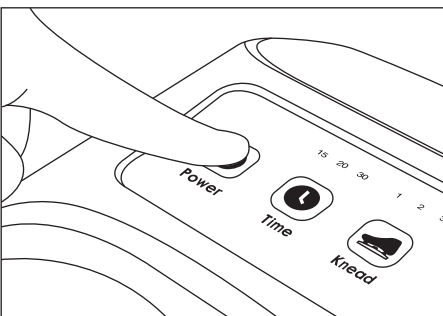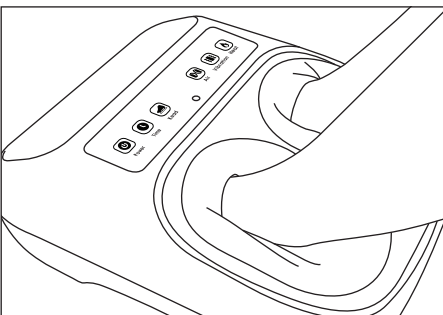Default setting: Time indicator lights up with 10-minute duration




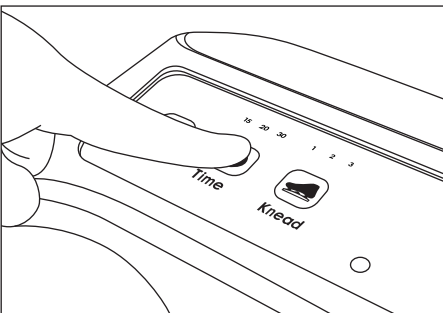
### 2.Preparation

Place your feet in the massage area



### 3.Time setting

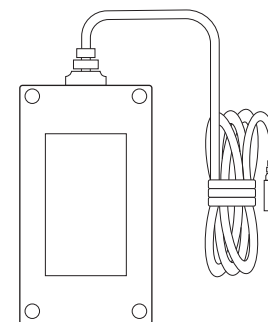
Press  to select your preferred session duration (15/20/30 minutes) .Default time is 15 mins




## WHAT'S IN BOX



Foot massager



Power supply

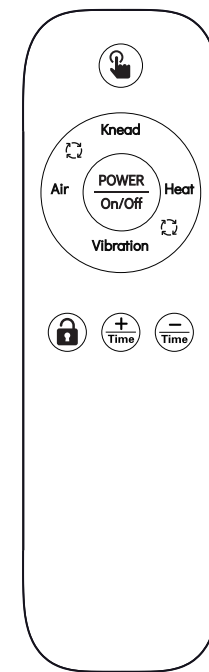

Remote control

## SPECIFICATIONS

Product name	Foot massager
Power	48 W
Frequency	50/60HZ
Input Voltage	AC 100-240 V
Output	12 V; 4A
Timer	15-20-30 min
Temperature	95-107.6 °F