

# Qure



[www.qureskincare.com](http://www.qureskincare.com)  
[support@qureskincare.com](mailto:support@qureskincare.com)

**P19-0023/MJ-03**

**Q-REJUVALLIGHT PRO FACEWEAR  
USER GUIDE AND MANUAL**

---

Includes important safety information.  
Read all instructions prior to use.

## **Q-Rejuvalight Pro Facewear Gives You Clinical-Grade Facial Treatments. For a fraction of the cost. Without You Leaving Your Couch.**

Bringing the latest science-based proven light therapy technology to your home in a revolutionary, hands-free portable, wearable design.

Q-Rejuvalight Pro Facewear is your clinical-grade skin treatment that you can do easily, all by yourself, in the comfort of your own home. Like going to a beauty clinic 5 times a week—without leaving your couch.

Proven skin treatments, tailored to your skin. Qure's customizable anti-ageing and anti-breakout technology allows you to treat your skin, everyday, and see immediate, longer-lasting results.

Qure's anti-ageing and anti-breakout technology, scientifically proven, repairs fine lines, boosts collagen production whilst simultaneously destroying and preventing acne. Completely improving the texture and quality of skin, leaving you with youthful, smoother and brighter complexion.

Now, you can get clinical results in your own home, on your own time.

Comes with a compatible iOS/Android app that allows you to fully customize, track and plan your treatments. You have full control of your skincare treatments, with a personalized program that's unique to your skin profile. Download the Qure Skincare app and step into the world of personalized skincare.

## GET STARTED


Follow these steps to make sure the strap is well connected to your Q-Rejuvalight Pro Facewear mask.

Clip 1 Clip 2




### STEP 1

Position the mask so the LED lights face towards you.

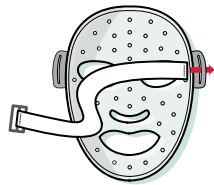


### STEP 2

Slide the strap through clip 2. Make sure that the Qure logos on the strap are facing down.

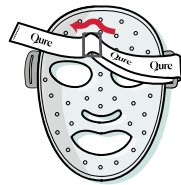

### STEP 3

Insert the strap through the slide adjuster. You should now see the Qure logos on the strap.

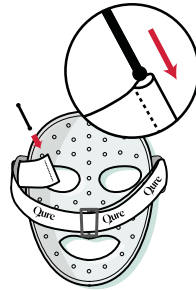

### STEP 4

Slide the strap through clip 1, push the pin through the end of the strap, and secure the pin by pulling the strap.

## WHAT'S INCLUDED



Q-Rejuvalight Pro Facewear



Detachable, adjustable head strap for universal fit



Universal USB charging cord (5V)



Qure Storage bag

### INDICATOR LIGHT

Lights up when the mask is turned on

### CHARGING PORT

Charge once, complete up to 10 treatments

### UNIVERSAL POWER BUTTON

Turns device on/off, puts device in Bluetooth pairing mode & selects treatment modes.

### SIDE ADJUSTER

Effortlessly adjust the size of the strap for a comfortable fit.

### STRAP


Adjustable and removable strap that holds mask in place.

### CLIP

Easily remove and re-attach the strap

### MASK

Hands-free portable, wearable design. Universally fits all facial shapes and features, both men and women

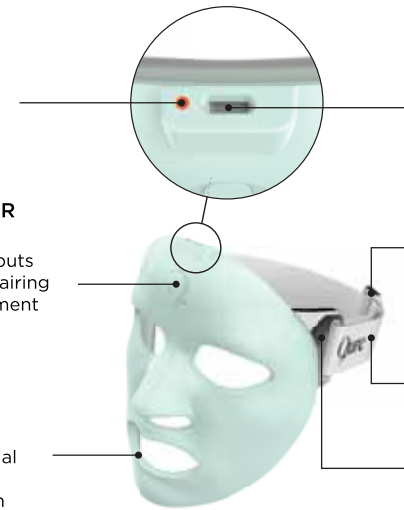

## SPECIFICATIONS MATERIALS

Color: Light Blue  
Size: 18.2 cm (h) x 14.8 cm (w) x 9.6 cm (d)  
Weight: 200g  
Interface: One button operation

## IMPORTANT SAFETY INFORMATION

- Do not operate the Q-Rejuvalight Pro Facewear where it could come into contact with water.
- Do not use or charge the Q-Rejuvalight Pro Facewear with power adapters that are damaged and/or faulty.
- Do not harshly bend or abuse the wire of the USB charging cord.
- Do not attempt to take the Q-Rejuvalight Pro Facewear apart. It is not serviceable by the user. This will void the warranty.

### Disclaimer

This product is not intended to diagnose, cure or prevent any disease. Qure Skincare makes no claims, representations, or warranties regarding the ability of this product to cure any physical, skin, or mental condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

### Product Representation

Qure is constantly improving products with the latest technology. As a result, the actual product may differ slightly from the product shown in images on the packaging and in this user guide.

## WARNINGS

- Do not harshly bend the device.
- Do not sit or stand on the device.
- Do not aim or look directly into the lights. Doing so may cause temporary vision discomfort. Placement of the device should ensure that no part of the device is covering the eyes.
- This device is designed for use on adults over 18 years of age and not recommended for use on infants or children.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant.
- Please consult your doctor prior to use if you have, or suspect you have cancer.
- Do not use the device if it is wet or if it could fall into water while it's plugged in. If the device has fallen into water, do not touch it until you first unplug it.
- Do not continue treatment if your skin is red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than 2 hours you may have light-sensitive skin, stop using the device and contact a healthcare professional.
- Do not use this device to treat skin that is painful to touch, raised (higher than surrounding skin), crusting or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- This device contains a Lithium Ion Battery. Before disposal of this device, please review your local laws and requirements surrounding Lithium Ion Battery disposal. The preferred method of disposal is through recycling.

## CAUTIONS AND CONTRAINDICATIONS

### Cautions

- Discontinue use of the device if it overheats or becomes uncomfortably hot.
- Always unplug device from outlet after charging is complete.
- Do not store devices in direct sunlight or location with temperatures over 100° r-Fahrenheit (38°C) as this may cause damage to the Q-Rejuvalight Pro Facewear.
- Use only the supplied Universal USB Power Cord (5V 50/60hz). Use of incorrect power cord could damage the device.
- Do not attempt to open or repair the device or power supply.

### Contraindications

*Ask your doctor prior to using if you:*

- Are taking prescription or over-the-counter drugs known to cause light sensitivity.
- Have any suspicious or cancerous lesions.

### Test yourself for light sensitivity prior to using on your face:

1. Power on your device.
2. Place the device on your arm and conduct a 3 minute treatment.
3. After treatment is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Discontinue use and contact Qure support.


## Q-REJUVALIGHT PRO FACEWEAR OVERVIEW

Q-Rejuvalight Pro Facewear is your customized, science-backed light therapy treatment. You can now get clinical-grade results in the comfort of your own home and see immediate, longer-lasting results.

Q-Rejuvalight Pro Facewear is a safe, natural, and effective treatment for rejuvenating your skin, reducing wrinkles & fine lines, lessening the appearance of scars & blemishes, treating acne, and boosting collagen. Completely improving the texture and quality of skin, leaving you with youthful, smoother and brighter complexion.

Backed by a huge base of clinical research, it has been studied and tested across hundreds of peer-reviewed clinical trials with overwhelmingly positive results for skin health, collagen production, improved appearance, reduction of inflammation, and helped people look younger.

### HOW QURE'S LIGHT THERAPY TECHNOLOGY WORKS

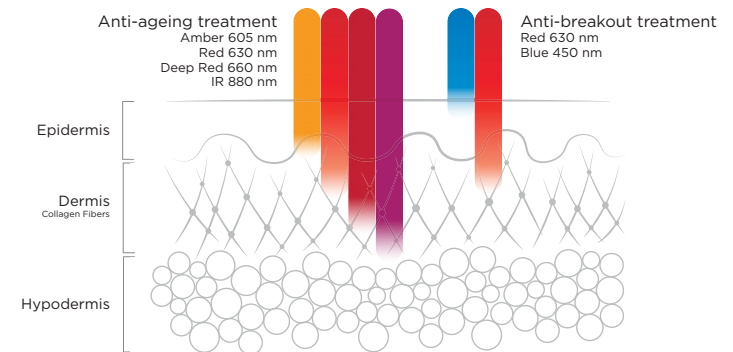

Q-Rejuvalight Pro Facewear delivers safe, concentrated wavelengths of natural light directly to the underlying layers of your skin, with no chemicals, UV rays, or excess heat, with clinical-grade LEDs to power up your cells.

These wavelengths of natural light stimulate the mitochondria in our cells to produce energy more efficiently. That increased energy allows your fibroblast cells to create more natural collagen.

Collagen protein makes up 30% of the total protein in the body, and 70% of the protein in the skin! As a person ages, the epidermis (outer layer of skin) becomes thinner and loses elasticity through a process called elastosis. As this happens, a person tends to show more signs of aging and acquires more wrinkles.

Q-Rejuvalight Pro Facewear light therapy restores more healthy cellular function, stimulating the production of collagen—which is why so many people have experienced the skin & anti-ageing benefits of a clearer and smoother complexion, with less fine lines and wrinkles.

Qure's Anti-Ageing Red Light Therapy stimulates your skin cells that other non-red wavelengths cannot reach and thus cannot affect the underlying tissue in the way that Qure's technology does. The benefits are numerous;



Brightens  
Complexion



Reduce Wrinkles  
and Fine Lines



Increase Collagen  
Production



Increase Healing  
and Recovery



Decrease  
Inflammation

Qure's Anti-Breakout light therapy technology treats mild to moderate acne, killing the acne-causing bacteria. It eliminates the bacteria found in the oil glands of the skin and causes the oil glands to shrink in size. The result is less oil in the skin and less bacteria, leading to less acne in the treated areas.

Q-Rejuvalight Pro Facewear Anti-Breakout Light Therapy not only destroys acne bacteria, but reduces redness, inflammation, and prevents future breakouts.

Q-Rejuvalight Pro Facewear is programmed with a 3 minute automatic shutoff function for professional anti-ageing and acne-fighting results with short, simple, day by day treatments in the comfort of your own home producing in-clinic results.

Q-Rejuvalight Pro Facewear contains 2 different treatment modes:

- Anti-ageing
- Anti-breakout

#### What To Expect:

Users typically see initial results within 2 weeks, with the best results occurring after a 8-week course of treatment

## Q-REJUVALIGHT PRO FACEWEAR TREATMENTS

Customize your treatment via the Qure Skincare app. Choose between two treatment modes (anti-ageing and anti-breakout) and tailor your treatments on 5 facial areas to suit your skin's needs.

### **Reduce fine lines and wrinkles:**

Increases production of collagen and elastin.

*Select Mode 1: Anti-Ageing*

### **Acne breakouts:**

Kills topical bacteria causing breakouts and speeds up recovery time, giving skin a fighting solution to acne.

*Select Mode 2: Anti-Breakout*

### **Even skin tones:**

Helps even skin complexion for naturally beautiful skin.

*Select Mode 1: Anti-Ageing*

### **Reduce redness:**

Alleviates Inflammation and rosacea.

*Select Mode 1: Anti-Ageing*

### **Revitalization:**

Promotes a brighter glow for more radiant skin.

*Select Mode 1: Anti-Ageing*

### **Soothing:**

Supports in reducing the size of swollen capillaries and alleviates inflammation and rosacea.

*Select Mode 1: Anti-Ageing*

### **Skin brightening:**


Helps break down deposits of melanin, hyperpigmentation brown spots, within the skin for even and more youthful-looking skin.

*Select Mode 1: Anti-Ageing*

### **Increases hydration:**

Increased blood flow brings more blood and nutrients to the area, for more hydrated and glowing looking skin.

*Select Mode 1: Anti-Ageing*

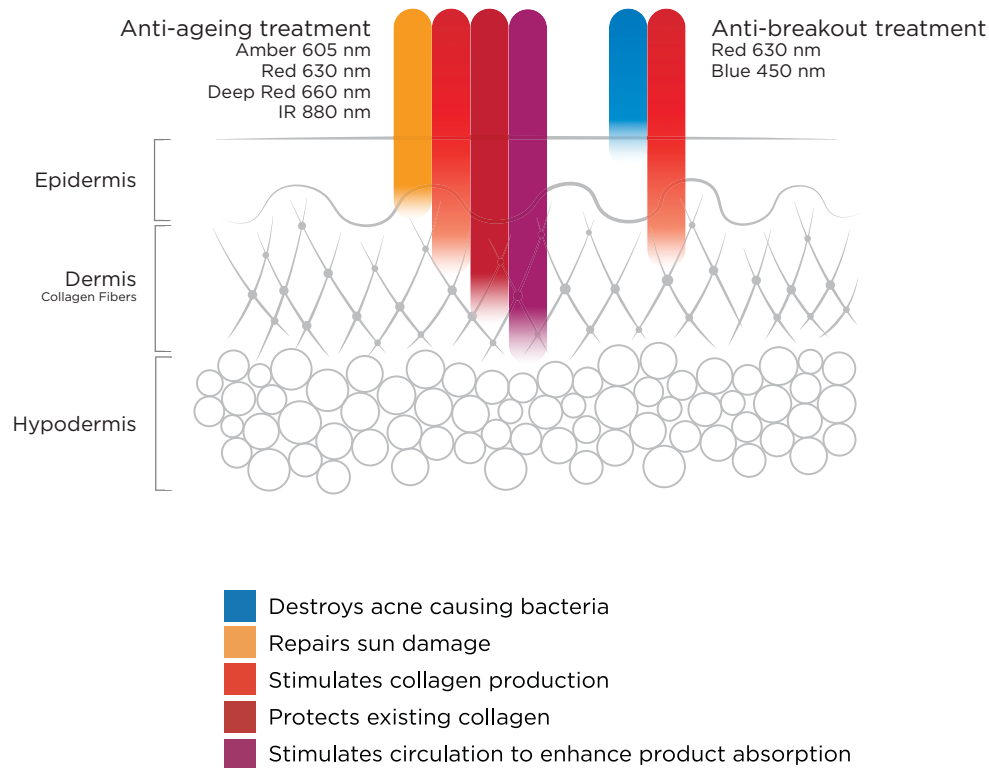

## HOW TO OPERATE YOUR Q-REJUVALIGHT PRO FACEWEAR

### Set up:

1. Unpack your Q-Rejuvalight Pro Facewear and follow these steps below.
2. Unpack the USB-C cable and plug it into a universal USB power adapter (not included in the box) and into outlet for charging (use only UL or other locally certified and compliant USB power adapter).
3. Plug USB-C power/charging cord into the port connector into the Q-Rejuvalight Pro Facewear.
4. Charge your Q-Rejuvalight Pro Facewear for 2 and a half hours. Device is fully charged when the Power RED indicator light next to the charging port is turned OFF.
5. Gently unplug the USB-C power/charging cord from the Q-Rejuvalight Pro Facewear after fully charged.
6. Place the Q-Rejuvalight Pro Facewear in the protective storage bag when not in use. Do not store device in direct sunlight or locations with high temperatures.

## INSTRUCTIONS FOR USE:

### Phase One:

1. Download the Qure Skincare app.
2. Open the app on your mobile device and follow the instructions to set up your profile and pair your Q-Rejuvalight Pro Facewear with the app.
3. Turn on your Bluetooth on your phone.
4. Turn on Q-Rejuvalight Pro Facewear by holding the center button on the forehead until the light is on.
5. Once your Q-Rejuvalight Pro Facewear appears on the Bluetooth available devices list, click on it to connect.
6. White indicator light appears when Q-Rejuvalight Pro Facewear is paired.



### Phase Two:

1. Wash and fully dry your face prior to use, removing all make-up .
2. Assemble the strap to the Q-Rejuvalight Pro Facewear by sliding the ends of the strap behind the hooks of the clips.
3. Place and position the mask onto your face and adjust the strap (use the slide adjuster) to fit comfortably around the head and make sure the mask stays in place.
4. Follow the instructions in the app to start and customise your treatment. You can choose from two treatment options: Anti-Ageing and Anti-Breakout.
5. To choose/change the treatment without using the app, follow the instructions below:
  - Power on the device by holding the Q button on forehead until the red indicator light is on
  - Press the button once for Anti-Ageing Treatment
  - Press it again for Anti-Breakout Treatment
  - Switch off the device by holding down the Q button until the indicator light turns off
6. Keep Q-Rejuvalight Pro Facewear on your face. Unless treatment time is set differently on the app, the device will automatically turn OFF after completion of 3-minute treatment.
7. Remove Q-Rejuvalight Pro Facewear from face and wipe device clean (see cleaning instructions).
8. Place Q-Rejuvalight Pro Facewear into the provided protective storage bag.
9. After treatment, apply creams or topical treatments as desired to skin.

### Helpful Tips for the Best Results:

- Use Qure Skincare App to customize your treatments, set reminders and to get informed about Qure's light therapy technology.

- Use your Q-Rejuvalight Pro Facewear daily.
- Clean your skin before use for an optimal result.
- Do not use the device after the application of fatty creams or products that contain SPF as this could block light penetration.
- Clean Q-Rejuvalight Pro Facewear after each use.
- Don't look into the light directly, this might cause an inconvenient temporary blindness.

### Maintenance:


1. It is important to keep the Q-Rejuvalight Pro Facewear clean and to clean after each use (see cleaning instructions).
2. This device has no serviceable parts. Battery is not accessible and not serviceable. Users should not attempt to access or service the battery as this will damage the device and void the warranty.
3. Q-Rejuvalight Pro Facewear includes only an ON/OFF button. All settings are set to the factory setting.
4. Charge Q-Rejuvalight Pro Facewear only with the supplied USB-C power/charging cord and a universal USB-C power adapter (not included).
5. Power indicator light will flash red when battery levels are low. This indicates it is appropriate to begin recharging the Q-Rejuvalight Pro Facewear.
6. Compatible Power Supply: 5V 50/60hz)



"white indicator" light  
mask turned on



"red indicator" light  
battery indicator

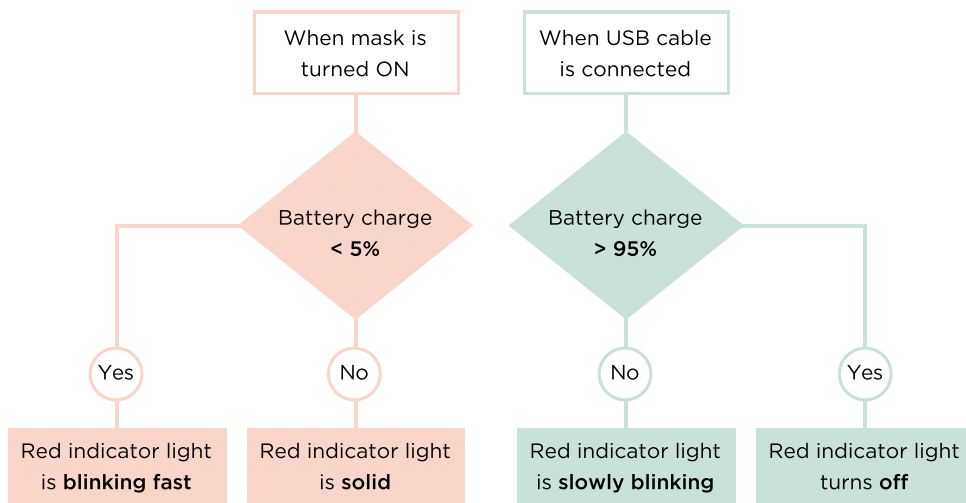

### Cleaning:

1. Turn Q-Rejuvalight Pro Facewear OFF and unplug from any power source.
2. Wipe the inside of the mask with a damp cloth (make sure the cloth is not too wet as it might damage the electrical components) to remove any residue from skin oils after each use. This is especially important if the device is being used by more than one individual.
3. Clean Q-Rejuvalight Pro Facewear body by using a damp or microfiber cloth. Do not use chemical or alcohol-based cleansers on the exterior.
4. Allow Q-Rejuvalight Pro Facewear to dry completely prior to use.

### Troubleshooting:

1. Q-Rejuvalight Pro Facewear is not turning on:
  - a. Is Q-Rejuvalight Pro Facewear fully charged? If the white indicator light is blinking then please charge the Q-Rejuvalight Pro Facewear.
  - b. Is the power cable connected to the Q-Rejuvalight Pro Facewear and the RED indicator light still turned off? - Please contact Qure customer support.
  - c. If the Q-Rejuvalight Pro Facewear still does not turn on, please contact Qure customer support.
2. Q-Rejuvalight Pro Facewear is not charging:
  - a. Is the USB-C power cord firmly inserted into the Q-Rejuvalight Pro Facewear? (Does the power indicator light turn on?)
  - b. Try a different USB-C Power Adapter to ensure it is working properly.
  - c. Try a different wall outlet to determine if there is power to the outlet.
  - d. If the Q-Rejuvalight Pro Facewear still does not charge then please contact Qure customer support.

3. Q-Rejuvalight Pro Facewear shuts off abruptly, prior to completion of 3 minute treatment:
- Q-Rejuvalight Pro Facewear is designed to shut off automatically after 3 minutes.
  - Q-Rejuvalight Pro Facewear may shut off after multiple repeated uses if it exceeds desired temperature. Allow the Q-Rejuvalight Pro Facewear to cool for several minutes and try again.
  - If the Q-Rejuvalight Pro Facewear continues to shut off prior to 3 minutes, please contact Qure support.

## LIMITED WARRANTY INFORMATION

We value you as a client and strive for the highest quality products to earn your satisfaction. Please use this as your reference for any customer service or warranty needs.

Warranty is limited to the terms set out below:

### ***What is covered by the limited warranty?***

This 1-year limited warranty covers defects in materials and workmanship in your Q-Rejuvalight Pro Facewear when purchased through an authorized partner or directly from Qure Skincare. Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred.

If the product was a gift then please obtain a copy of the original receipt or have the original purchaser process the warranty claim on your behalf.

### ***What is not covered by this limited warranty?***

This limited warranty does not cover: damage, defects, failure caused by abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, or loss of device.

### ***How long does this limited warranty last?***

The limited warranty begins from the date of purchase and expires after one year. The warranty period is not extended if your device is replaced by Qure during the warranty period.

### ***What will Qure do if your device is defective?***

If you discover a defect in your Q-Rejuvalight Pro Facewear covered by the warranty, before the 1-year warranty period expires, contact Qure customer support for assistance. If we confirm a defect covered by this limited warranty, Qure will decide whether to replace your device with either a new or refurbished device. If we find that your device is not defective, we will notify you and return the device to you. All Limited Warranty services must be authorized and performed by Qure or a service provider authorized by Qure.

## WARRANTY LIMITATIONS & DISCLAIMERS

This warranty gives you specific legal rights, and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, Qure disclaims all implied warranties. If applicable law does not permit the disclaimer of implied warranties,

then the duration of implied warranties are hereby expressly limited to the duration of the express limited warranty set forth above.

In no event shall Qure be liable for any special, incidental, consequential, punitive or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability) or otherwise, even if Qure is advised or aware that such damages may arise. Such damages include, but are not limited to: loss of Qure, cost of any substitute equipment, the claims of third parties or damage to property.

If the applicable law does not allow limits on warranties, or on remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between Qure and you, the purchaser. This allocation is recognized by both parties and is reflected in the price of Qure. By using Q-Rejuvalight Pro Facewear you acknowledge that you have read these terms, conditions, and limitations, understand them and are bound by them.

### **FCC Statement**

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.