


# iMagic™

MODEL# 17BT0302

## DIGITAL SMART SCALE USER MANUAL

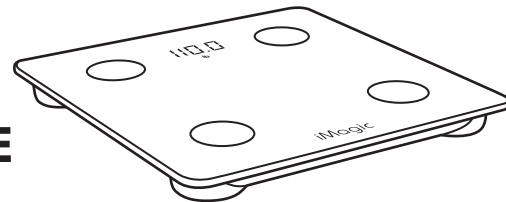Thank you for choosing the iMagic™ Digital Smart Scale and we hope that you enjoy using it for your health tracking and monitoring needs. Using Bio-electrical Impedance Analysis (BIA) technology to offer a full range of data in personal Weight, BMI, Body Fat, Body Water, Basal Metabolic Rate(BMR) and Metabolic Age will provide assistance in achieving goals and a healthier lifestyle.
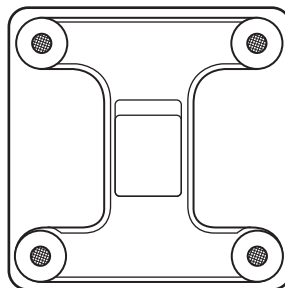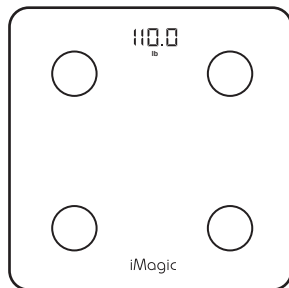
### PACKAGE CONTENTS:

Digital Smart Scale  
Batteries: 1.5V AAA x 3  
User Manual

### SPECIFICATIONS:

Size: 11" x 11" x 1"  
Units: lb / kg  
Max. Capacity: 396 lbs /180 kg  
Increments: 0.2 lbs/ 0.05 kg

### IMPORTANT SAFETY INSTRUCTIONS:

The iMagic™ Digital Bluetooth Scale is a highly accurate precision measuring device and certain basic safety precautions should be followed when using:

- Never jump or step aggressively onto scale. Never place heavy objects on scale
- To avoid damaging tempered glass always handle and move scale carefully
- Avoid stepping directly on the edge of scale
- Feet and scale surface should always be dry before using
- Do not attempt to disassemble
- Do not use device or electrode if it has been damaged. The continuous use of a damaged unit might result in injury or serious harm
- Keep scale in dry conditions. Avoid exposure to moisture, heat, dust, lint or direct sunlight
- Do not use scale to diagnose or treat any medical condition. You should consult your physician(s) when undertaking any diet or exercise program
- Anyone that is pregnant or has a pacemaker should consult their physician(s) before using
- Do not mix different types of batteries such as alkaline, carbon-zinc, or rechargeable batteries
- Do not mix old and new batteries. Always replace the whole set of batteries at one time
- When not in use for long durations, please remove batteries

### OPERATING INSTRUCTIONS:

#### 1) BEFORE OPERATING


- To ensure accuracy always place scale on even and hard surface when in use. If using for the first time or if scale is ever moved to a new location, wait 10 seconds before stepping on scale. We recommend that measurements be taken at the same time of day

- When stepping on scale the iMagic™ logo should be closest to you. This will allow the LED display to be read properly
- Confirm phone is connected to network
- The iMagic™ Digital Smart Scale can connect to your mobile phone via "iMagic Home" app and bluetooth. Download "iMagic Home" app from App store for Iphone or Google Play for Android phone. Make sure your phone is Apple iOS 7.0 or above, or Android 4.3 or above, and with bluetooth 4.0 or above

#### 2) OPERATING

- Download "iMagic Home" app from App store for Iphone or Google Play for Android phone
- Register your own account
- Add your scale device "17BT0302" on the app per steps
- Input your personal data to start your health tracking experience
- For Weight Measurement use (Note: bare feet are not required), stand on scale evenly and wait for results to be displayed on scale or sent to phone via "iMagic Home" app
- For Body Fat and other analysis (Note: dry bare feet are required), stand on scale allowing feet to make contact with (4 silver) electrodes
- Listen for notification on phone to view Health Report
- For multiple users repeat above steps

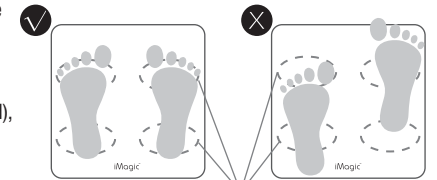Electrodes

### TROUBLE SHOOTING:

| NO DISPLAY                              | NO BLUETOOTH CONNECTION                        | LOW BATTERIES POWER | OVERLOAD                                              | CALIBRATING                                    |
|-----------------------------------------|------------------------------------------------|---------------------|-------------------------------------------------------|------------------------------------------------|
| No display                              | 6800<br>kg                                     | Lo                  | EE-                                                   | CAL                                            |
| Check batteries for proper installation | Make sure Bluetooth is turned on in your phone | Change batteries    | Exceeding maximum weight capacity of 396 lbs / 180 kg | Wait 10 seconds to allow calibration to finish |

### WARRANTY:

- The scale has been checked for accuracy prior to leaving the factory and has a guarantee against any defect in workmanship and parts for a period of one (1) year.
- Within one (1) year if you wish to make a claim, please call US & CANADA Toll Free #1-833-4-iMagic (8:30 am - 5:00 pm E.S.T.) Monday - Friday for assistance.

#### Caution:

This device complies with Part 15 of the FCC. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

\* Apple and Iphone are registered trademarks of Apple Inc.  
\* Google and Android are registered trademarks of Google Inc.  
\* Bluetooth is a registered trademark of Bluetooth SIG Inc.

**SAVE THESE INSTRUCTIONS FOR FUTURE USE!**

Distributed by:  
Great Star Industrial USA, LLC  
Huntersville, NC 28078 USA

For assistance, call:  
1-833-4-iMagic  
www.imagicdevice.com  
MADE IN CHINA

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This device complies with Industry Canada's licence-exempt RSSs. Operation is subject to the following two conditions:

(1) This device may not cause interference; and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication.

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MPE Reminding

To satisfy FCC / IC RF exposure requirements, a separation distance of 20 cm or more should be maintained between the antenna of this device and persons during device operation.

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\* Apple et iPhone sont des marques de commerce déposées d'Apple Inc.
\* Google et Android sont des marques de commerce déposées de Google Inc.
\* Bluetooth est une marque de commerce déposée de Bluetooth SIG Inc.