

User Manual



FOOT MASSAGER

360°ALL DIRECTIONAL FOOT MASSAGER

Table of contents

Safety Instructions

Warranty

Product Components

Features and Use Diagram

Setting

Specifications

Troubleshooting

Maintenance

Safety Instructions

➤ Important Safety Instructions

- 01) Be sure to read the user manual in a detailed way prior to use to ensure safety and proper use of the product
- 02) Make sure to use the product by following the instructions in this manual
- 03) After using or before cleaning the product. please disconnect to the power
- 04) Always cautiously use the heating function to avoid serious burn
- 05) Do not apply the product to too insensitive and over sensitive skin or to body parts of poor blood circulation
- 06) Never allow children or those of physical and mental disorders to use the heating function or improper operation that may be dangerous
- 07) To prevent fire,electric shock or injuries caused by overheating, never use it under a pillow or a blanket
- 08) Do not cover it during use so as to avoid malfunction or damage resulting from motor overheating
- 09) Children, those of physical and mental disorders and the disabled should be under guide when using the product
- 10) Please stop using it immediately if you feel uncomfortable
- 11) When the power cords or plugs are damaged,inappropriate operation is performed and the product has been dropped or damaged or dropped into water.do not use this product please contact the service center tor check and repair
- 12) Do not use spare parts that are not recommended by the manufacture
- 13) Do not insert or drop objects into any opening of the product
- 14) Please close all operation control before disconnecting the product to the power if you want to cut off the power
- 15) This product cannot be used for self-diagnosis or treatment and please visit professional medical care if in need of diagnosis or treatment

➤ Use Environment

- 01) Do not use the product in a wet or dusty environment in order to avoid product failure
- 02) Do not use this product at room temperature over 40 °C
- 03) Do not use the product outdoor
- 04) Do not use the product next to an electric heater or scores or in direct sunlight

Safety Instructions

➤ People Not recommended Using the Product

- 01) Please consult your doctor before using if you have the following situations
- 02) Receiving medical treatment or experiencing physical discomfort
- 03) Suffering from hypertension.heart disease. cerebral arteriosclerosis.and cerebral thrombosis
- 04) Suffering from cancer
- 05) Pregnant or in a period
- 06) Suffering from osteoporosis
- 07) Using a pacemaker or implantable medical devices
- 08) Suffering from acute diseases
- 09) Minors


➤ Buttons

- 01) Do not operate the button with wet hands
 - 02) Keep the buttons from water or other liquid
 - 03) Do not fall asleep when using the product
 - 04) Do not use this product when you are drunk or feeling uncomfortable
 - 05) Do not directly contact the massager pad with your hands, feet joints, knees chest, and abdomen
 - 06) Do not use this product over thirty minutes in a row
 - 07) Do not use if the product's cloth cover has been damaged
 - 08) Do not use the product one hour before and after eating so as to avoid adverse effects
 - 09) Do not enable the heating function when the product is being used in a low room temperature It is recommended that you should gradually adjust the temperature to a normal range before enabling heating
 - 10) If the product is used immediately after being transferred From a cold environment to a warm environment, the accumulation of water vapor inside the product may affect the normal operation of the product or even cause mechanical failure
- You are suggested to put the product in normal room temperature for 1 hour before starting to use it after the transfer

Warranty

- 01) Warranty is not transferable
- 02) Users can enjoy free warranty from the date of purchase during the warranty period(The surface is not covered under warranty)
- 03) The product warranty is valid for 12 months from the invoice date
- 04) If we write that supplier shall provide a warranty card
- 05) The purchasing date shall be based on the effective receipt or, invoice issued by the retailers
- 06) This warranty shall not apply to leasing, trade, ecommerce, institutions and other non-buyer users
- 07) The following situations are not covered the warranty
 - a) Destruction caused by human factors, including the non-normal using environment usage and storage that are not compliant to the instructions, and damage caused by improper use, etc
 - b) Disassemble, repair and modification of the product without our authorization
 - c) Damage caused by users poor transportation after Purchase
 - d) Damage resulting from any other force majeure (such as flood, fire, lightning and earthquakes)

Product Components

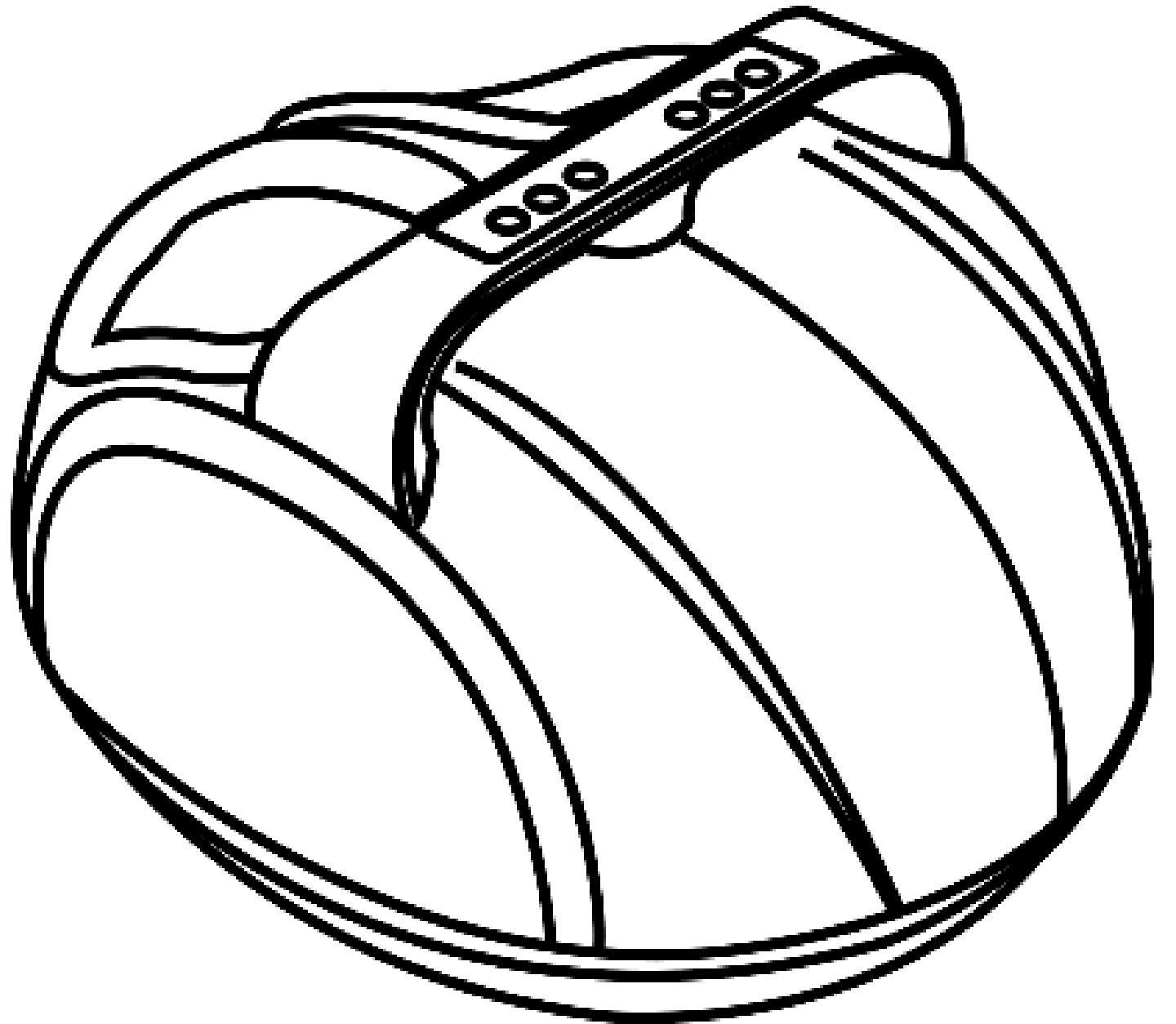

- **Central Surface:** Six independent buttons – controlling the product's settings and functions
- **Central Massager Pad:**Scratching-Acupressure on the arches with heating
- **Front Massager Pad:**Roller Kneading-Acupressure on the front soles
- **Rear Massager Pad:**Push Kneading-Acupressure on the heels
- **Interior airbag:**Air pressure on the insteps
- **Removable and washable fabric cover:**cleaning
- **Bottom:**Cord Wrap-keeps things tidy

Features and Use Diagram


➤ Product Features

- 01) Egg-shaped design.fashionable and beautiful.Lightweight and portable
- 02) A combination of a variety of function including scratching.roller kneading and push kneading acupressure, air pressure, and heating
- 03) Massager Pad that provides all-around massages and in -deep acupressure to clear the body's meridians and foot reaction zones
- 04) Removable and washable cover, allowing you to keep clean at any moment
- 05) Cord wrap. bringing convenience to life

➤ Use Diagram



Settings

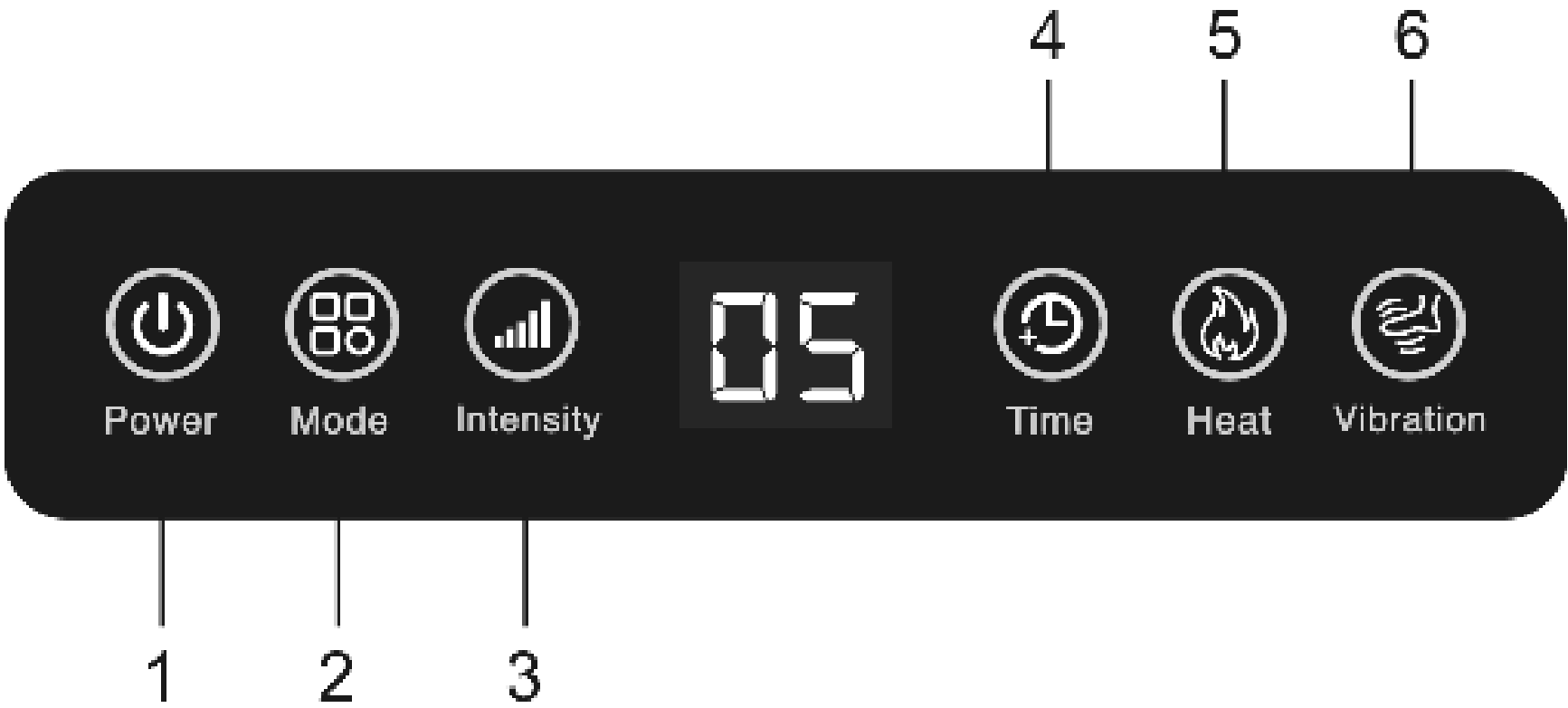

1.Power: Turns on/off the machine

2.Mode: Switch between three massage modes (default mode 1)



3.Intensity: Controls the massage intensity,to your preference (9 levels,default intensity 5)

4.Time: Set massage time (05-30mins,default 15mins)

5.Heat: Press the button to control the heating independently (default on)

6.Vibration: U1 U2 two levels (default on)

Specifications

Model	TMF-01SB
-------	----------

Voltage	110V 60Hz
---------	-----------

Power	50W
-------	-----

Troubleshooting

Please refer to the following several tips and possible solutions if problems occur during usage of this product

The product may be returned to the service center or undergone inspection and repair by authorized technicians

Problems	Possible Causes	Remedies
Not working	1.Unplugged 2.Power button not pressed	1.Connect the product to power or check whether or not the product is well-connected 2.Press the power button
Suddenly stops running	1.Power disconnection 2.Operation duration over 15 mins 3.Product overrun and overheat 4.Too much operating pressure that the massage balls stop rotation	1 .Check whether or not the product is plugged 2.Restart the product 3.Restart the product after cooling 4.Restart the product

Maintenance

➤ Storage

- 01) The product should be put into the color box and stored under safe.dry well-ventilatcd and cool condition when it is not in use
- 02) Do not twist switching the power button and the cord
- 03) Do not hang the product with the cord
- 04) Avoid scratching by sharp objects on the surface of the product
- 05) Do not put the product indirect sunlight or under condition of high temperature

➤ Cleaning

- 01) Make sure the power button is off and unplug the power cord before cleaning
- 02) Use a dry cloth to clean the surface of the product
- 03) Do not hang the product with the cord

➤ Note

- 01) Do not use corrosive detergents.such as benzene or thinner.to clean the product

➤ Maintenance

- 01) Keep the product clean and away from high humidity conditions
- 02) Please use a dust cover(not include in the box)to cover the product to prevent dust from coming into it if it is not used for a long period

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.