

# BLUETOOTH EYE MASSAGER

Operating instruction



Please read the instruction carefully before use

## PLEASE READ THIS INSTRUCTION BEFORE USE

Thank you for buying this eye massager. Please read the operation manual carefully before use, and pay special attention to safety matters, so as to operate and use this product correctly

This manual provides users with a safe and effective way of operation. Please keep it properly after reading it for future reference

## ORDERS TO RECORD

1. Safety precautions	-----	2
2. Product profile	-----	2-3
3. Please confirm the commodity content	-----	3
4. Name and function of each part	-----	4
5. product usage	-----	5
6. Cleaning, maintenance and storage	-----	6
7. Faults and solutions	-----	7
8. Product specifications	-----	7

## SAFETY PRECAUTIONS

If you wear contact lenses, please remove contact lenses before using this product.

If there is any dirt or foreign matter in the mask, please wipe it with a wet cloth before using it to avoid accidental eye.

During the charging process, the product is not supported to start up (charging is conducted under the working state, the product starts the self-protection program and stops working automatically).

Do not use the eye massager during outdoor walking or driving.

Do not use eye massager in hot and humid places.

## PRODUCT INTRODUCTION

The eye massager uses the microcomputer chip control technology for the entire human nature simulation control, the application of air pressure that a whole new concept to the eye care, under the control of air pressure in the computer to a full range of muscles around the eye pressure, vibration and coupled with the internal function and the function of the far infrared heating, the eye acid bilges, the eye pain, relieve eye fatigue.

### 1. FEATURES:

The eye massage device can be connected to mobile phone to play music via Bluetooth;

Portable wireless usage, USB charging mode;

Remote infrared heating microcomputer chip control technology;

5 massage techniques;

15 minutes automatic timing;

Air bag massage (massage the muscles around the eye area by proper matching between air filling and air releasing);

5V voltage output is more secure, one key operation is simple and convenient;

New vibration mode;

## 2. SCOPE OF APPLICATION:

For the eyes

## 3. CONTRAINDICATIONS:

Disables eye injuries or inflammation.

If there are symptoms of eye diseases can not be used, general implementation of eye surgery, glaucoma, cataract, retinal detachment and other major eye diseases patients should not be used.

For patients with deep stroke, hypertension and hypotension, please consult a professional physician first. If you feel any discomfort during use, please stop using immediately.


Please avoid using the eye massager during drinking or physical discomfort

PLEASE CONFIRM THE COMMODITY CONTENT





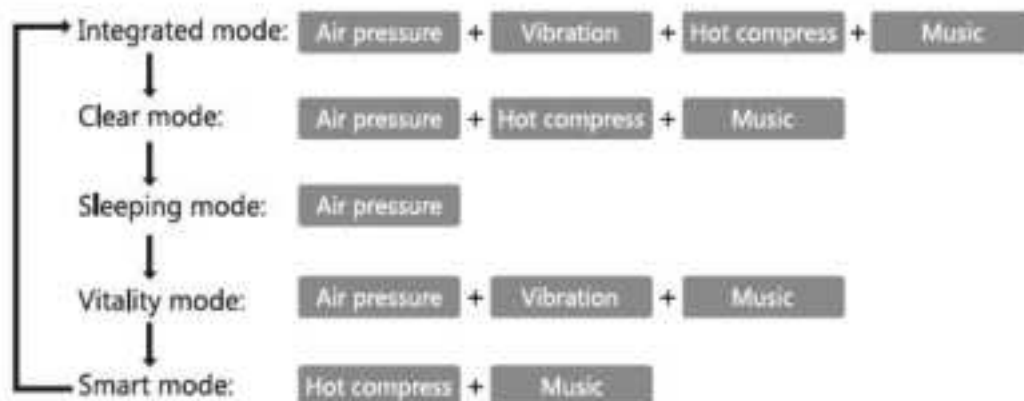

Eye massage device, USB charging cable, Instruction manual



## ONE KEY CONTROL LOOP



## NAME AND FUNCTION OF EACH PART

- (1) Long press the on/off key for 3 seconds to start the integrated mode ( air pressure + vibration + hot compress + music)
- (2) Press the  key to start clear mode (air pressure + hot compress + music)
- (3) Continue to press the  key to start sleeping mode (air pressure)
- (4) Continue to press the  key to start vitality mode (air pressure + vibration + music)
- (5) Continue to press the  key to start smart mode (hot compress + music)
- (6) Press the on/off key twice to pause/play the music
- (7) Press the  key to adjust the temperature (low temperature, medium temperature and high temperature)
- (8) Press the  button to adjust the volume.
- (9) Press the on/off key for 3 seconds to turn off the machine function
- (10) Open the Bluetooth of your mobile phone and search the Bluetooth name: YBY.  
The devices will automatically connect.




## PRODUCT USAGE

### Method of use

#### 1. Preparation before use

- (1) take out the eye massager gently
- (2) put the eye massage device on the eye position

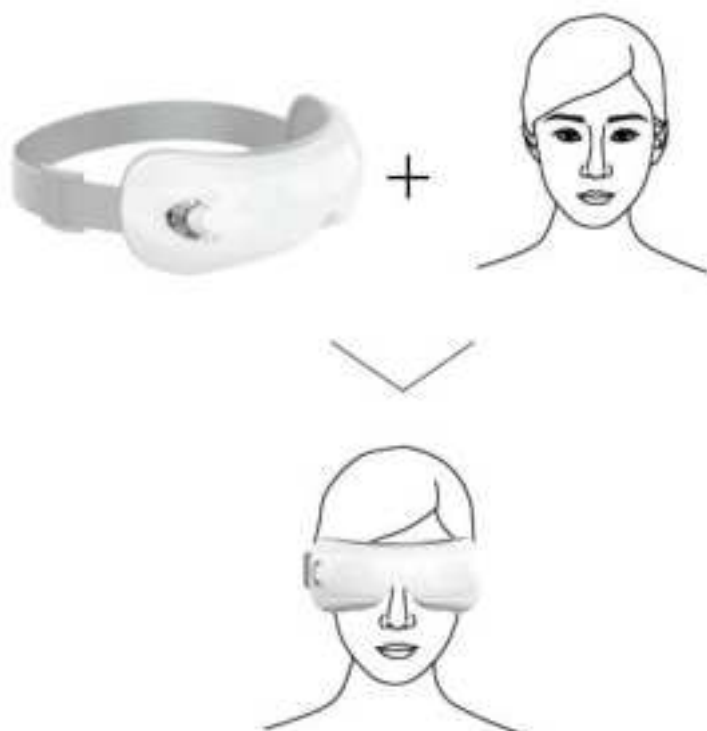


## PRODUCT CLEANING, MAINTENANCE AND STORAGE

Product cleaning, storage and maintenance

### 1. Product cleaning and maintenance

If the product is dirty, please wipe it with a dry soft cloth.

If it is particularly dirty, you can wipe it with medical alcohol.

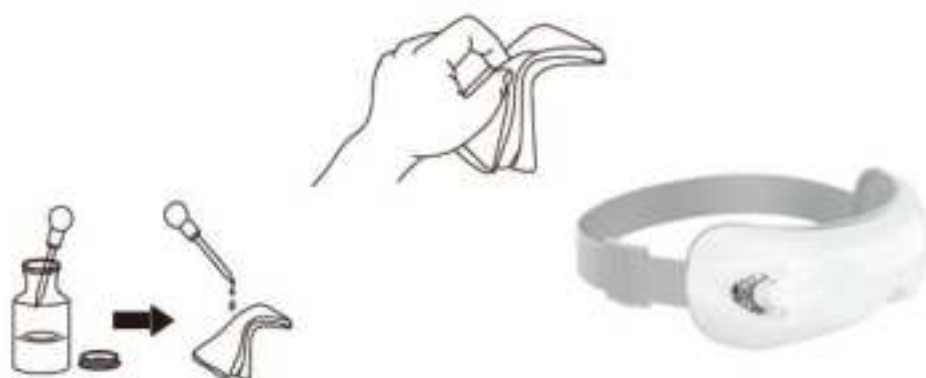

### 2. Keep

Keep it dry, ventilated and clean

Please keep out of reach of children

Do not disassemble the main engine without permission

## COMMON FAILURES AND SOLUTIONS

### Unable to boot:

1. If the battery power is too low, please charge immediately
2. Host protection (with electricity) connect the charger, charge, and allow the host to automatically remove the protection (charge for 5 seconds)
3. The battery is not charged and the charger is damaged
4. The pressure intensity is too weak

**Note:** all maintenance work must be done by our authorized maintenance personnel

## PRODUCT SPECIFICATION

The product name	Eye massager
Automatic timing	15 minutes
The rated voltage	DC 5V
Rated power	3.4W
Power adapter	AC 100-240V DC5V



## COMMON FAILURES AND SOLUTIONS

The fault	Possible cause of failure	The solution
Unable to boot	Low battery	Please charge immediately
	Motherboard protection (with electricity)	Connect the charger to charge, and allow the host to automatically remove protection (5 seconds after charging)
	The battery is not charged	Replace charger
	In charge	Unplug charger
The pressure is too weak	Low power	charging
Eye patch not hot	Low power	charging
	Model is wrong	Tune to the corresponding pattern

## FCC Statement

This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: 1) this device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

## FCC RF Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.