


EQL
essence

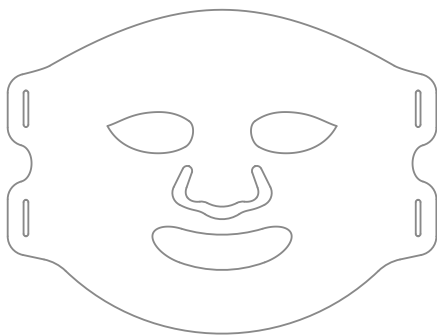


AURO Mask

User Manual

- 1 Product & Accessory
- 2 AURO Mask - Classic
- 3 Mask Appearance
- 4 How To Use Controller
- 5 How To Start
- 6 APP Download & Connection
- 7 APP Function
- 8 AURO Mask - Premium
- 9 Mask Appearance
- 10 How To Use Controller
- 11 How To Start
- 12 APP Download & Connection
- 13 APP Function
- 14 Notice
- 15 Maintenance, LED Indicator
- 16 Product Specification
- 17 Federal Communication Commission
Interference Statement
- 18 Warranty

1 Product & Accessory

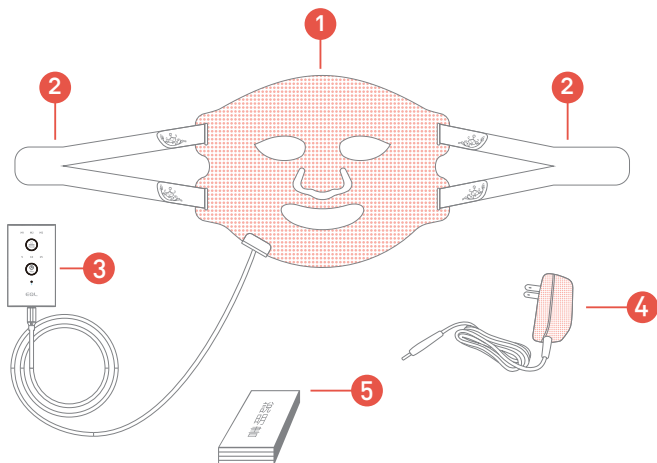

① AURO Mask x1

② Strap x2

③ Controller x1

④ Charger x1

⑤ User Manual x1



AURO Mask - Classic

RUBY | FOR WRINKLES

GOLD | FOR BETTER SKIN TONE

3 Mask Appearance



- 1 Venetian style mask
- 2 Strap with unique pattern
- 3 Velcro

How To Use Controller

4



1 Power switch button

2 Time indicators

3 Timer switch button

4 Bluetooth indicator


5 Power adapter port

6 Reset hole

7 Mask-connecting port

Power On	Long press power button for 3 seconds, Bluetooth indicator will keep blinking
Select Time	Short press timer switch button
Shut down (Manual)	Long press power button for 3 seconds
Shut down (Auto)	The device will shut down automatically after idling for 5 minutes
Factory Reset	Use a small point object or needle to poke in the reset hole, the device will be back to factory settings

5 How To Start

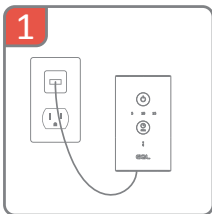


1. Charge the device
(please charge for 4 hours
before 1st time use)




2. Clean your face



3. Apply daily skincare products

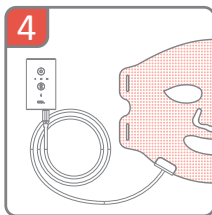

4. Plug mask into controller

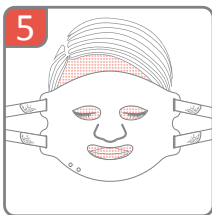

5. Put on the mask,
and adjust the
strap to comfort



6. Set up timer

★ if you want to use with AURO APP, please
check the guide for set-up on next page



APP Download & Connection

6

1 APP download


Search "AURO Mask" from Google Play or App Store.

Or scan QR code:

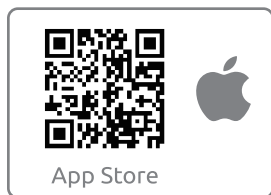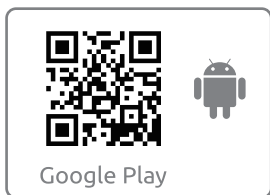

2 Connect to the AURO Mask


Turn on your smart phone Bluetooth and you will find the AURO Mask, press connect

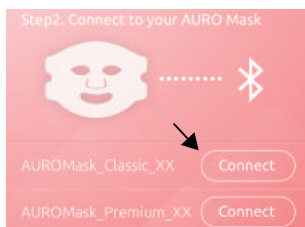

3 Name the AURO Mask




4 Start to Use

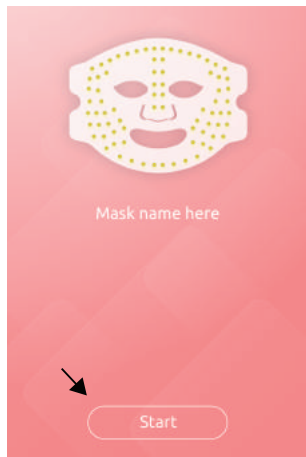


7 APP Function

(Please go to eqltech.com.tw to download complete user guide for AURO)

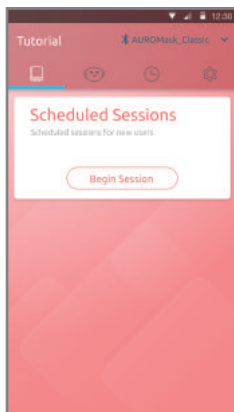


Tutorial program

A 28-day tutorial program designed for new users

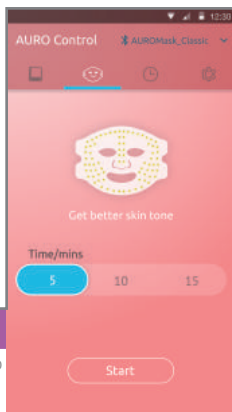


AURO Control

Replace your controller to operate the AURO Mask

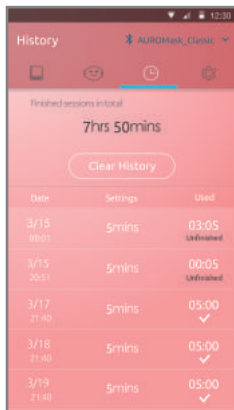

History

Track your skincare treatments

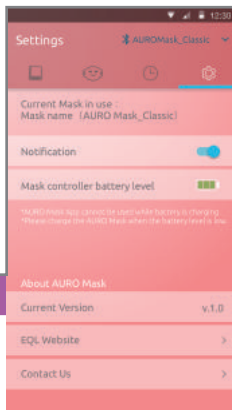

Setting

Turn on/off notification, check battery status, other supports from EQL

AURO Mask - Premium

M1 | 紅黃光

立體緊緻

水嫩勻亮

M2 | 紅黃藍光

無暇美肌


光澤淨膚

M3 | T字藍光加強

長效控油

平衡收斂

9 Mask Appearance

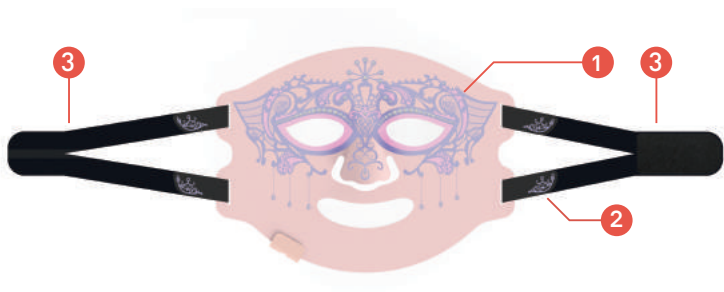

- 1 Venetian style mask with diamond
- 2 Strap with unique pattern
- 3 Velcro

How To Use Controller


10



- | | |
|----------------------------|------------------------|
| 1 Mode indicators | 5 Bluetooth indicator |
| 2 Power/Mode switch button | 6 Power adapter port |
| 3 Time indicators | 7 Reset hole |
| 4 Timer switch button | 8 Mask-connecting port |

Power On	Long press power button for 3 seconds, Bluetooth indicator will keep blinking
Select Time	Short press timer switch button
Select Mode	Short press Mode switch button
Shut down (Manual)	Long press power button for 3 seconds
Shut down (Auto)	The device will shut down automatically after idling for 5 minutes
Factory Reset	Use a small point object or needle to poke in the reset hole, the device will be back to factory settings

11 How To Start

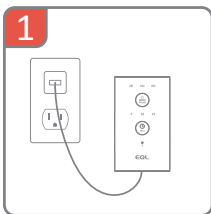


1. Charge the device
(please charge for 4 hours
before 1st time use)




2. Clean your face



3. Apply daily skincare products

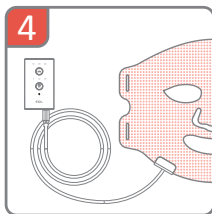

4. Plug mask into controller

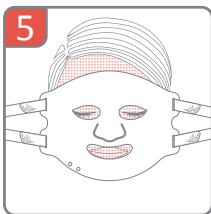

5. Put on the mask,
and adjust the
strap to comfort





6. Set up timer and
mode

★ if you want to use with AURO APP, please
check the guide for set-up on next page

1 APP download


Search "AURO Mask" from Google Play or App Store.

Or scan QR code:

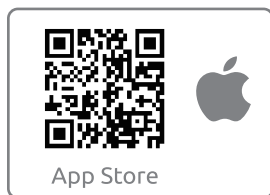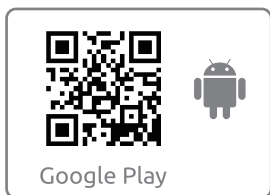

2 Connect to the AURO Mask


Turn on your smart phone Bluetooth and you will find the AURO Mask, press connect

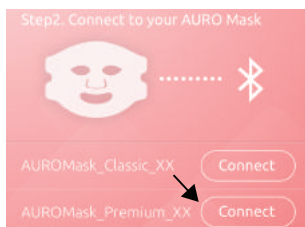

3 Name the AURO Mask




4 Start to Use

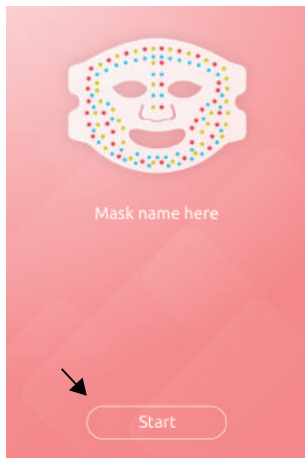

13 APP Function

(Please go to eqltech.com.tw to download complete user guide for AURO)


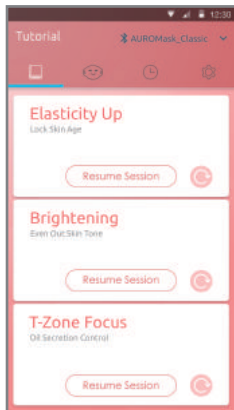


Tutorial program

A 28-day tutorial program designed for new users

AURO Control

Replace your controller to operate the AURO Mask

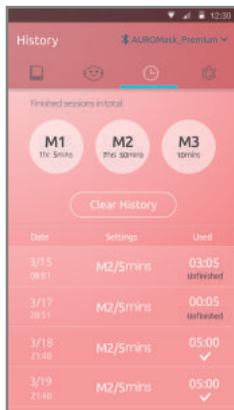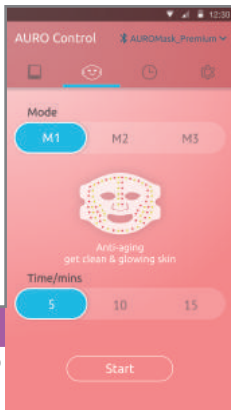


History

Track your skincare treatments

Setting

Turn on/off notification, check battery status, other supports from EQL

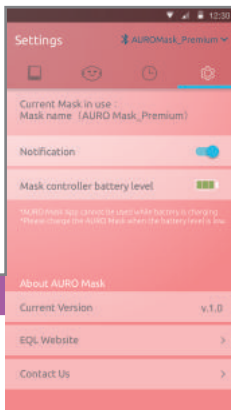

1. Please charge 4 hours before 1st time use.
2. DO NOT USE the device while charging, APP function will be suspended while the device is charging.
3. DO NOT USE if the AURO mask has fallen in water or is damaged in anyway. Water can damage the parts inside.
4. DO NOT USE tissue mask when you are using the AURO mask, it may reduce light therapy effects.
5. DO NOT USE more than 30 minutes per day.
6. USE ONLY the power supply and controller that came with your AURO mask.
7. DO NOT use the AURO if you think there maybe something wrong or defective.
8. DO NOT look directly at the light. Please turn on the power after wear the mask.
9. DO NOT USE the AURO mask if your skin is sensitive to light. You may know if you are sensitive to light as some medicines you take, or products you use, can cause light sensitive skin.
10. DO NOT USE the AURO mask if you are pregnant or suspect that you may be pregnant without consulting a doctor first.
11. CALL a doctor right away if you think that you may have had a negative reaction from using the device.
12. Controller's function will be suspended when doing tutorial skincare program, please use AURO APP to operate the device.

15 Maintenance

1. Use a soft, damp but not wet cloth to clean the silicone surfaces. You may also use alcohol on a soft cloth to clean the silicone surface.
2. Store the AURO Mask in a safe, dry and cool place.
3. Keep the AURO Mask away from fire.
4. DO NOT bend the AURO Mask when not in use.
5. DO NOT wash the AURO Mask.

LED Indicator

Bluetooth Indicator	Blink	AURO Mask is NOT connected to smart phone
	Stay on	AURO Mask is connected to smart phone
Time Indicator	Blink (Single LED)	Low battery
	Circulate	The device is charging
	Stay on	Full battery

Mask Size	270mm x 204mm x 15.6 mm
Mask Size	50mmx 86mm x29.8mm
Power Supply	Input : AC 100 ~ 240V Output : 50~60Hz Daya : DC 5V/ 2A
Connectivity	Low Power Bluetooth
Operation Temperature	37 °C +/- 2 °C

Federal Communication Commission

Interference Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

FCC Caution:

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Warranty

EQL warrants that the product will be free from defects in materials and workmanship under normal use for a period of 24 months after its original purchase date. If the product proves defective during this warranty period, EQL, at its option, will either repair the defective product or will provide a replacement in exchange for the defective product.

A service fee may be charged in the following situations:

- Improper handling of the product without referencing the product manual.
- Damages due to human error, humidity, or self-disassembly/assembly.
- Damages due to unforeseen disasters and circumstances, such as fires, floods, typhoons, and earthquakes, etc.
- Damage to software, consumables, and natural wear-and-tear.
- If the product is not purchased from EQL authorized channels

EQL



Facebook



LinkedIn



HDH5A001010

EQL
essence