UREVO SMART SCALE

User Manual

IMPORTANT: Read all instructions carefully before using this product.
Retain this owner's manual for future reference. The specifications of this product may vary from this photo, subject to change without notice.

**Manufacturer: Shenzhen Yile
Dynamic Technology Co., LTD.**

Address: 1205-1206, Tianliao Building,
Tangchang Road, Tanglang Community,
Taoyuan Street, Nanshan District,
Shenzhen, China

Version: 1.0

MADE IN P.R.C

Activate UREVO App Now
For Longer Warranty &
Smart Features Access

Mode.URWS001

 **UREVO**

Table of Content

02	EN	DISPOSAL
03		FCC STATEMENT
04		PRODUCT SPECIFICATION
04		DISPLAY INSTRUCTIONS
04		SAFETY PRECAUTIONS
05		DIRECTIONS FOR USE
06		APP DOWNLOAD AND MEASURE
08		SMART SCALE WARRANTY POLICY

Disposal

The packaging material is recyclable. Please dispose of it in an environmentally friendly manner and take it to a recycling center.

Dispose of the device in an environmentally friendly manner.

Dispose of it at a recycling center for used electrical and electronic devices.

You can obtain more information from your local authorities.

The batteries should be removed prior to disposal of the device and disposed of separately from the device.

To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but must be taken to suitable collection points.

Also note the relevant statutory regulations governing the disposal of batteries.

FCC Statement

CAUTION: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

CAUTION: This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Statement:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Product specification

Size	278*278*24mm
Weight	1.18kg
Battery	3*AAA
Division value	50g
Measurement Range	0.3-180kg
Display	LED
Material	ABS, Safety tempered glass
PS:Specific information for each model, According to the outer packing	

Display instructions

1. Lo(low battery)
2. C (test error)
3. OL(Overweight prompt)
4. Unable to boot or display (replace the battery)

Note: After the product appears "OL", please take out the battery and reinstall it, which can be measured within the weighing rang

Safety precautions

Be sure to read safety precautions before use.

- This product should not be used by those equipped with pacemaker,artificial lung or wearing electrocardiograph.
- Do not use this product on slippery surfaces

Directions for use

The correct way to measure

- Please put the product on the hard flat floor, do not put on the carpet or others soft ground, or it will cause inaccurate measurement data.
- Clear before each measurement.
Reset method:After lightly stepping on the scale, the screen number will automatically return to 0.0lb or 0.00lb


place the scale on level ground

Do not place on carpets or soft floors

The correct way to stand

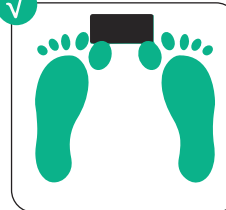APP download and measure

Step 1: APP Download

- Scan QR code to download UREVO SmartCoach App.
- iOS can be downloaded through the App Store and Android can be downloaded through Google Play


Step 2: How to connect the scale for the first time

- Click "Add Device", enter the search device page, click "weight scale" to search;
- Activate the scale, search for the scale device and connect it;
- You can start taking weight data;


Step 3: Weigh & View data

Go to the Scale home page to view measurements and more.

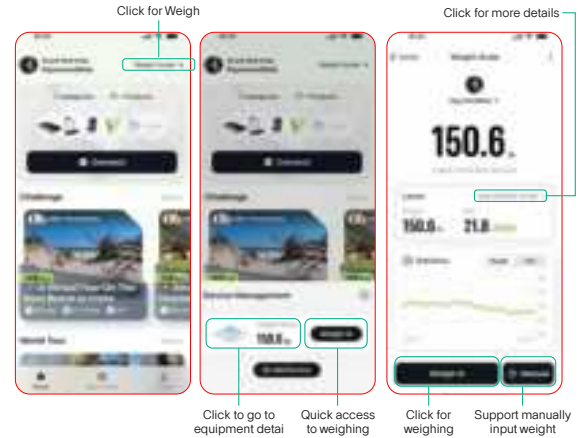Units Switch:

The scale is designed with two kinds of weighing units to choose, kg/lb. Press the unit key at the back of scale to choose the weighing unit you'd like to use.

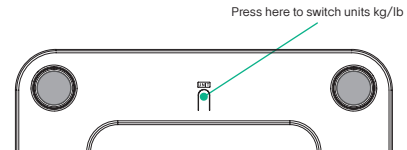SMART SCALE WARRANTY POLICY

1. From the date of purchase, you can enjoy a year of replacement service.
2. Within 15 days of purchase, the product quality problems (non-man-made damage), the merchant bears the back and forth freight, product quality for a year, within a year to provide new.

UREVO, using its sole discretion, will either repair or replace free of charge any part(s) proven to be defective under normal home use. Any repair or replacement shall provide no new warranty coverage, but shall retain only the remaining portion of the original product's warranty. This warranty is offered only to the original purchaser and is not transferable. Proof of original purchase is required.

Name	
Cellphone number	
Purchasing date	
Business name	
Product model	

Maintenance record	The fault phenomenon	
	Repair content	
	The date of repair	
	Maintenance personnel	

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