


CARETAKERSENTRY with two-way voice pendant

Quick Install Guide

Instructions for Installing Your New Personal Emergency Response System

Base Bottom View

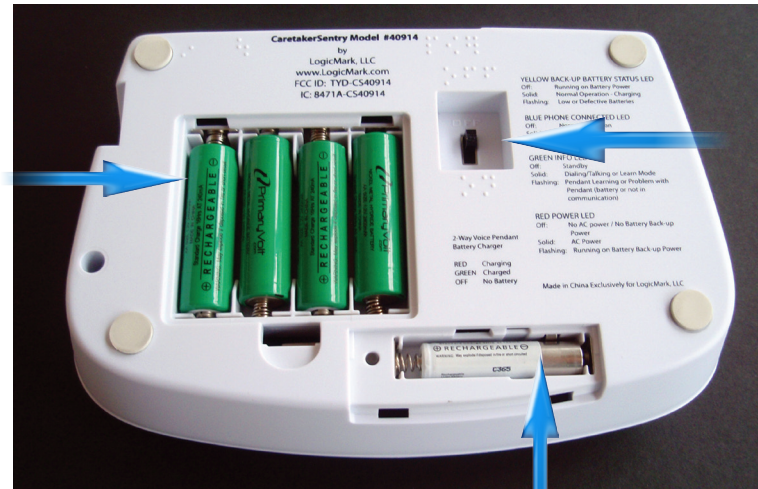


Backup Batteries
* **Note Orientation of batteries**

Master
ON/OFF
Switch



Two-way Voice Pendant Battery
Charging Port (for optional Two-way
voice pendant accessory)
* **Note Orientation of battery**

Base Back View

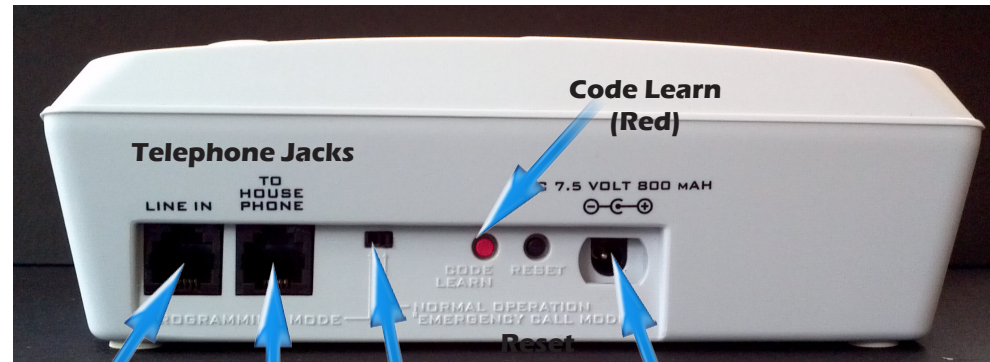

Line IN

To House
Phone

Programming Switch

Programming Mode "Emergency Call Mode"
Normal Operation


Power DC 7.5
800 mAh

Base Side
View



Cellular Accessory Jack

Pendant
Part #41915

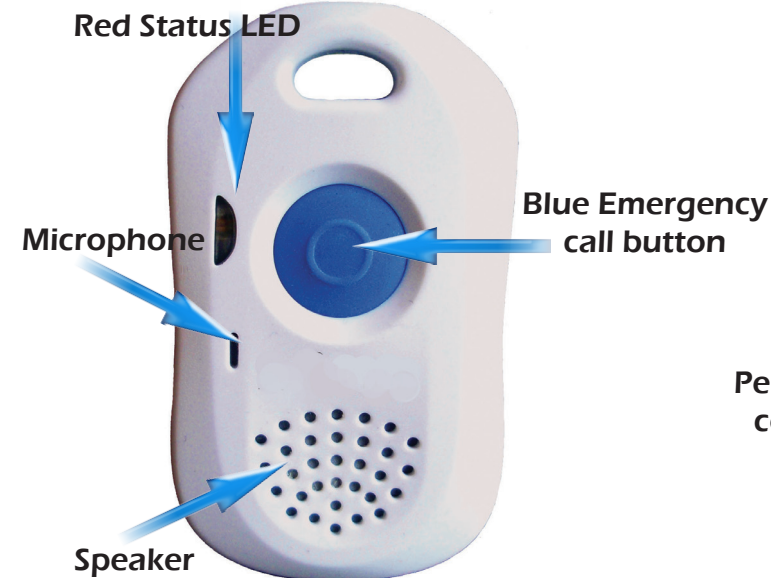


Red Status LED

Microphone

Blue Emergency
call button

Speaker

Pendant battery
compartment

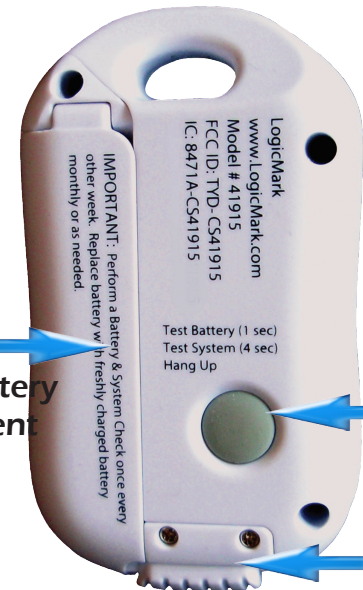

Gray System Check/
Battery test button

Pendant battery
cover latch

Base Top View
Part #40914



Speaker

Cancel (1 second)

Blue CareTaker
button

Red Emergency
button

Microphone

LogicMark Disclaimers and Limited Warranty

COMMUNICATION AND RESPONSE LIMITATIONS: Purchaser acknowledges that signals which are transmitted over telephone lines, or other modes of communication pass through communication networks wholly beyond the control of LogicMark LLC and are not maintained by LogicMark LLC, and, therefore, LogicMark LLC shall not be responsible for any equipment or communication failure which prevents transmission signals from reaching your contact list including emergency 911 operators or damages arising therefrom. **Purchaser acknowledges that LogicMark LLC provides no response to the System's equipment.** The equipment is designed to communicate with the a central station of your choice and LogicMark LLC is not and shall not be responsible for ambulance, police or other emergency response time or that any response will be provided by the central monitoring station.

TESTING AND SERVICE OF THIS EQUIPMENT: The equipment, once installed, are in the exclusive possession and control of the Purchaser, and it is Purchaser's sole responsibility to test the operation of equipment and request warranty service if the equipment is under warranty.

PURCHASER'S EXCLUSIVE REMEDY: Purchaser's exclusive remedy for LogicMark LLC's default hereunder is to require LogicMark LLC to repair or replace, at LogicMark LLC's option, any equipment or part of the personal emergency alert system which is non-operational during LogicMark LLC's warranty period.

LIMITATION OF LIABILITY: This equipment is not designed or guaranteed to prevent any loss or injury. This Limited Warranty and Disclaimer of Liability constitutes the terms of sale and use of the equipment, and if there should arise any liability on the part of LogicMark LLC as a result of any cause whatsoever, regardless of whether or not such loss, damage, or personal injury was caused by or contributed to by LogicMark LLC's negligence to any degree or failure to perform any obligation or strict products liability, such liability will be limited to an amount paid by the Purchaser to LogicMark LLC for the product, or to the sum of \$350.00, whichever is greater.

For warranty information, see your dealer for details.

© LogicMark, LLC v4.1

CaretakerSentry Base Unit - Model 40914
FCC ID: TYD-CS40914
IC: 8471A-CS40914
CaretakerSentry Two-Way Pendant - Model 41915
FCC ID: TYD-CS41915
IC: 8471A-CS41915
REN 3.7

CARETAKERSENTRY Step-by-Step Installation

1

Confirm Accessories in Kit

2 Lithium Ion
Pendant batteries



4 Rechargeable AA
Batteries



Telephone Cord



AC Adapter



Lanyard



*Wriststrap



*Beltclip



*The holster is the same for both Wriststrap and Beltclip. Simply unscrew the strap and replace with the clip or vice-versa

Model #41911

3

Plug Telephone Line and Power Cord into Wall Jacks



Plug the AC Adapter power cord into the power jack and twist 1/4 turn clockwise to lock.




(OPTIONAL) Plug the telephone into the "To House Phone" jack.

2

Place Batteries in Unit


Place the Green AA Rechargeable Batteries in the compartment on the bottom of the base unit. **Make sure to turn ON/OFF switch to ON. You will hear the announcement "Running on Battery Power." Unit will do so until the AC Adapter power supply is plugged in.**

Install 1st pendant battery in pendant.
(When battery is properly installed, pendant will automatically perform a battery check - listen for these voice prompts everytime a battery is installed.)

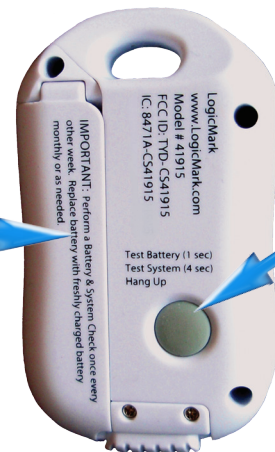

**Battery should standby in the pendant for 4-6 months (depending on use). We recommend to exchange/charge the batteries monthly.

Install 2nd pendant battery into Two-way Voice Pendant Battery Charging port
RED LED = Charging
GREEN LED = Battery Charged
No LED = defective battery or no battery installed

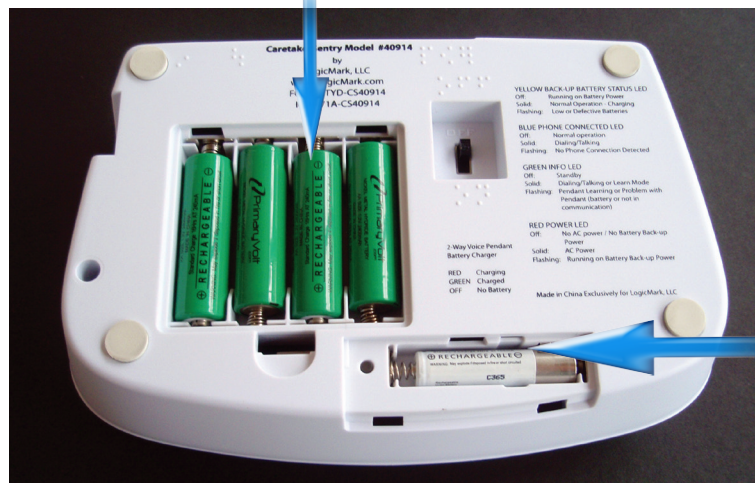

4

Test the Unit

Performing a System Test

- Push GRAY EMERGENCY button on back of pendant for **4 seconds**
- Listen for voice prompt from base unit:
 - "All Systems are OK – Battery is OK" is heard if the system is correctly installed
 - If any other voice prompt is heard, the unit will describe the problem
- We recommend that you perform this System Test weekly**

Place call to Emergency Operator:

- Push BLUE EMERGENCY button for 2 seconds on Pendant (RED LED on the pendant should turn ON)
- You should hear a dial tone, dialing and then communication tones between the CaretakerSentry and central station. Several seconds later the operator should answer the call
- Talk normally into the pendant– just like speaker phone
- At the end of the conversation, the CaretakerAlert should automatically hang-up – if not, push the CANCEL button for 1 second (either on back of pendant or on base unit.)

