FLOSSOM 花至

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FLOSSOM Wearable Intelligent Flexible Beauty Apparatus

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Thank you for choosing our FLOSSOM beauty apparatus.

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For correct and safe use of the product, please read this manual carefully before use. Keep it properly for future reference.

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Contents

Product specification	01
Interface and buttons	01
Who and where the product is not intended for	02
Functions	03
Directions for users	04
Usage	05
Cautions and safety warnings for use	06
Troubleshooting	09
Q & A	10
Disposal of built-in rechargeable lithium battery	11
Hazardous substances or elements	12
After-sales service information	13
Warranty card	14

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Product specification

Product name: FLOSSOM wearable intelligent flexible beauty apparatus
Model: WE200-M
Dimensions: Approx. W 165mm x D 155mm x H 48mm (skinpilot)
Battery capacity: 7.4V 250Ah
Rated input: 5V/2A (skinpilot)
Charging duration: 60min. (depending on the charging conditions)
Battery life: 60min. (depending on the operating conditions and the level selected)
Standard: GB4706.1-2005/GB4706.15-2008
Place of origin: China

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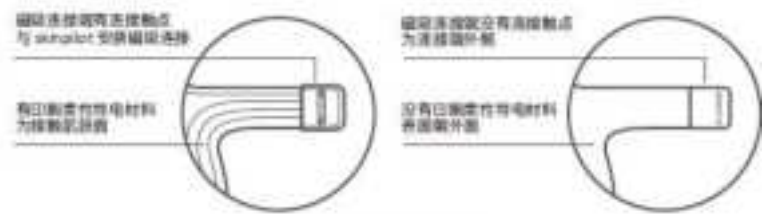
Who and where the product is not intended for

The product is a small household appliance instead of a medical device.
Please read the following content about the forbidden behaviors and the behaviors needing attention carefully before use and use the product as instructed herein. The special circumstances below are not exhaustive. If it is difficult to judge whether the product is suitable for you, please consult a doctor or the customer service of our flagship store or the authorized store where the product was purchased or call our customer service hotline (400-848-0056) directly.

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For the following people, circumstances and body parts, the product is forbidden to use.

- Any person using electronic medical implants such as a heart pacemaker, any electronic medical device for life support such as an artificial heart-lung machine, or any wearable electronic medical device such as an electrocardiographic recorder;
- Any person using contact lens. (He/she can use the product after removing the contact lens.)
- Any person with a heart disease, cardiac dysfunction or tending to have any problem of this sort;
- Any person who cannot feel stimulation;
- Any person under 18 years old (including any infant or child). (The product and its accessories should be put in a place unreachable to children.)
- Any body part not designated for use in this manual (including the following ones):
 - Any area around the eye, such as the eyeball or upper eyelid;
 - Throat or prominentia laryngea;
 - Oral cavity or mucosa;
 - Private part;
 - Any skinpilot with damaged skin, rash, skin inflammation or cancer;
 - Raised nevus or viral warts;
 - Body part with a metallic implant;
 - Fracture site

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09

For the following people or circumstances, please use the product with caution or after consulting a doctor.

- Any person with any acute disease, tuberculosis, fever, hypertension, infectious disease, malignant tumor, blood disease, diabetes or other visceral disease, or currently receiving a treatment. (Consult a doctor before using the product.)
- Any person with sensitive or fragile skin, any current or previous skin disease, atopic dermatitis, allergic dermatitis, contact dermatitis, or a worry about the adverse effect on skin. (Consult a doctor before using the product.)
- Any person with angiotelectasis due to hepatosis or long-time use of any steroid hormone medicine. (Consult a doctor before using the product.)
- Any person with high myopia (above -6D). (Consult a doctor before using the product.)
- Any person in or after an eye or facial surgery. (Consult a doctor before using the product.)
- Any person with problems such as hemorrhage, swelling or scabbing after a medical cosmetic surgery. (Consult a doctor before using the product.)
- Any person with a neurological disorder (Consult a doctor before using the product.)
- Any person hospitalized for treatment and taking any medicine. (Consult a doctor before using the product.)
- Any person currently receiving dental treatment or with dental filling). (Consult a doctor before using the product.)
- Any person after excessive exposure to sunlight or with sensitive or fragile skin. (Use the product with caution.)
- Any person drunk or shorter after drunkenness. (Use the product with caution.)
- Any female in her menstrual period. (Use the product with caution.)
- Any female in her pregnancy or lactation period (Consult a doctor before using the product and use it with caution.)

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Safety function

Protection against current abnormalities: Current of the skinpilot can be detected, and its abnormal output can be controlled.

Wear sensing: The status of connection between the FLOSSOM flexible conductive film fabric and the skinpilot can be detected so that current is on when the connection is correct and stable.

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Connection identification

Connection between the magnetic connection ends on two sides of the product, and contact between the FLOSSOM flexible conductive film fabric and skin of the skinpilot can be detected. In the case of successful connection and contact, the skinpilot will give two beeps, and the level indicator will breathe in green, suggesting successful connection. At this time, the power/level selection button can be pressed again to output energy of level 1 to start skin care.

In the case of failed connection and contact, the level indicator will breathe in red. At this time, level selection and energy output will fail.

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Timing

Upon start of skin care, the product will begin to output a current and time. In the case of failed connection or contact as mentioned above, the product will pause its current output and timing. After being properly worn and pressed on the power/level selection button, it will continue with current output and timing to continue timing for care skin.

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Automatic power-off

If the skinpilot has not been operated in 10 minutes after power-on, it will be powered off automatically.

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Skin care control on applet

Current level adjustment: In skin care, the current level can be controlled on the applet: The current level for each facial area can be adjusted for personalized use of the product. The adjustment data of all the facial areas can be imported into the skinpilot through only one button and used offline directly in future skin care.

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Adaptability test


Before the first use of the product, please do a skin test to avoid skin problems: apply a little essence onto the inner side of the arm. In the case of discomfort, stop using the product.

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Charging method

Before the first use of the product, please fully charge the skinpilot as below.

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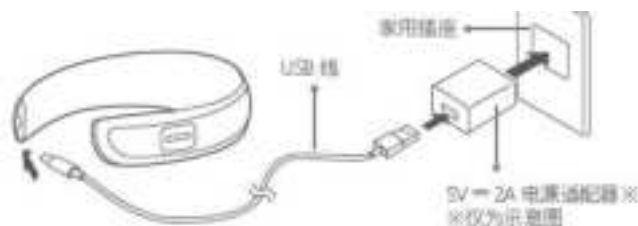Notes:

1. Please charge the skinpilot with the provided Type-C cable. The power adapter for charging should be prepared by user.
2. Before charging the skinpilot, please dry it (including its charging port) and the power adapter.

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
Charging duration and battery life

Charging duration: About 60 minutes (depending on the charging conditions)

Battery life: About 60 minutes (depending on the operating conditions and the level selected)

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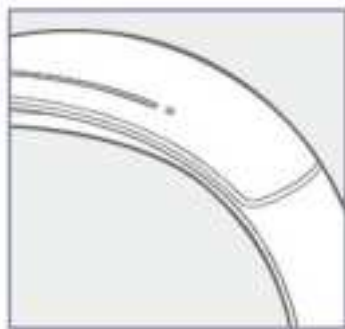When the skinpilot is used:

- The battery level indicator remains lit in white: High battery level
- The battery level indicator flashes in red: The product has a low battery level and needs to be charged.

When the skinpilot is charged:

- The battery level indicator remains lit in red: The battery is being charged.

The battery level indicator remains lit in green: The battery is fully charged.

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Usage

Recommended frequency of use

It is suggested to use the product twice a week. The frequency of use can be properly reduced according to the skin condition.

Control of skin care on applet

The skinpilot can be connected to the following mobile applet through Bluetooth.

- ① Before the first use of the product, bind it first.
- ② Scan the following code to enter the FLOSSOM applet.

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③ Bind the product.

- Bind the product as instructed in the applet.
- The product only needs to be bound for a single time. Control of skin care can be directly done on the applet after that.

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④ Connect the product to the applet.

- Enter the FLOSSOM applet and select the product.
- Power on the product, enable the Bluetooth function on the mobile phone, connect the product and do skin care.
- Some functions of the FLOSSOM applet may be updated.

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STEP 1: Remove the makeup and prepare for use of the product.

Before use of the product, please remove the makeup on the face and around the eyes and clean the skin.

Remove the contact lens and the metallic accessories (such as earrings, necklace, bracelets and watch) (if any).

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STEP 2: Power on and wear the skinpilot.

Press and hold the power/level selection button. Two beeps suggest successful power-on. After that, the product will be on standby, with its level indicator breathing in white.

As in the figure, take the product and wear it on the head like a hair band, making sure that its left and right sides are behind the ears and that the power/level selection button is on the right hand side.

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STEP 3: Fully wet the film fabric.

Pour the FLOSSOM-brand essence into the flexible conductive film fabric to wet it evenly for about 1 minute, making sure any position with the conductive material is fully wetted.

The remaining essence (if any) can be poured onto your palm and applied on the face gently.

 **Note:**

Keep the essence from the connectors on two sides of the flexible conductive film fabric. Apply it on the flexible conductive film fabric in contact with the face only.

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STEP 4: Wear the flexible conductive film fabric.

Wear the flexible conductive film fabric onto the face and adjust it to fit the face shape. Make sure its part with the printed flexible conductive material is in good contact with the face and do not fold it. The cut positions of the film fabric's part for the chin can be properly folded to fit the face shape.

Connect the magnetic connection ends on two sides of the film fabric magnetically to two sides of the skinpilot. When there are two beeps and the level indicator is breathing in green, it suggests successful connection and readiness for skin care.

△Note:

When connecting the two magnetic connection ends, make sure there's no hair between them and the skinpilot to ensure good conduction.

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STEP 5: Adjust current level and do skin care.

The current level can be adjusted as below on the skinpilot or the applet.

Adjustment on the skinpilot: Press the power/level selection button on the right side gently. The skinpilot will then output the current of level 1. The skinpilot will output the current of level 2 if it is gently pressed for the second time, or the current of level 3 if it is gently pressed for the third time. The output currents rotate in the sequence of level 1, level 2 and level 3.

Adjustment on the applet: Click the skin care page on the applet to set the current levels for different facial areas. After that, the applet will control the current output by the skinpilot, and the level indicator on the skinpilot will breathe in blue. The power/level selection button on the

skinpilot is disabled in this mode of adjustment till disconnection from the applet.

呼吸灯闪烁	The indicator is flashing.
亮灯	The indicator is lit.
电量指示灯	Power indicator
档位指示灯	Level indicator
蓝牙指示灯	Bluetooth indicator
开关/档位键	Power/level section button
一档	Level 1
二档	Level 2
三档	Level 3
小程序控制输出 呼吸灯为蓝色	Current output control on applet The indicator is breathing in blue.

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STEP 6: End the skin care.

Upon end of the skin care, there will be two beeps, the skinpilot will be powered off automatically, and the power indicator and the level indicator will both go out, suggesting that skin care has ended.

Take off the magnetic connection ends on two sides of the film fabric and remove the flexible conductive film fabric and the skinpilot in sequence. After that, face washing and further skin care can be done if you want.

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STEP 7: Clean and maintain the product.

Cautions in maintenance:

- After each use, please wipe the skinpilot with clean, soft cloth or tissues. Do not flush or soak it in water.
- Before maintenance, please power off the product to avoid electric shocks, short circuit and fires.
- Do not wipe the product with any solvent or medicine, such as alkaline cleaner, alcohol, thinner, gasoline or nail polish remover, or any abradant.
- The flexible conductive film fabric is a disposable product. Do not use it repeatedly.

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Cautions and safety warnings for use

Please use the product as instructed to avoid adverse effects.

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Use

- Please use the product for the purpose designated in this manual only.
- When in use, do not fold the printed flexible conductive material.
- When in use, do not pull the flexible conductive film fabric violently, otherwise there may be poor contact, rupture, etc.
- The flexible conductive film fabric is a disposable product. Do not use it repeatedly.
- Do not use the product while taking a bath.
- Do not throw the product into fire or use it close to any inflammable (such as alcohol, thinner, gasoline or nail polish remover sprayer).
- Do not use or charge the product when it is hot or being heated.
- Do not use the product when it smells abnormally, gives off smoke or is abnormally hot.
- Do not use the product when it is damaged.

- Do not disassemble, repair or refit the product.
- Do not use the product for a commercial purchase or another purpose needing long-time use.
- Avoid falls and collisions of the product.
- Do not scratch the flexible conductive film fabric with any hard object.
- Do not use the product on pets.
- Stop the product when you feel uncomfortable in use.
- Using the product repeatedly on a single site may result in hazards for the skin there.
- There may be adverse effects if this manual is not followed in use.
- The product contains a radio transmission module with a CMIIT ID of 2022DP11232.

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Charging and power supply

- Please keep the power adapter from water. In the case of accidental contact with water, dry it before use.
- Do not damage or refit the USB cable.
- After charging, unplug the power adapter with dry hands and keep it properly.
- In the case of any abnormality, power off the product immediately.

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Storage and safekeeping

- Do not wipe the product with any solvent such as alcohol, thinner, gasoline or nail polish remover.
- Do not store the product with essence on.
- Do not put the product close to flame or any heater, under sunlight or in any other hot or humid place.
- Please keep the product in a place unreachable to children.
- Before long-time storage of the product, please unplug the power adapter from the power socket and then unplug the USB cable from the skinipilot.
- Before storing the skinipilot, wipe it and put it into the storage bag.
- Please keep the product in a clean, well-ventilated place.

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Troubleshooting

Fault	Possible cause/recommended solution
Failure of power-on	The skinpilot has a low battery level. Please charge it.
	Please press the power/level selection button and hold it for 2s.
	The skinpilot has failed. Please contact the customer service of the authorized store where the product was purchased or call our customer service hotline (400-848-0056) directly.
Failure of charging	Please use the original USB cable and a 5V/2A power adapter prepared by yourself.
	The contact between the skinpilot and the power adapter, or that between

	the power adapter and the power socket is poor.
	The power socket is not electrified. Please check whether it is electrified and normal.
Failure of level adjustment	The contact between the magnetic connection ends of the flexible conductive film fabric and the skinpilot is poor.
	When the FLOSSOM applet is used, the highest level is imported into the skinpilot. In this case, the level of the skinpilot is set as the highest one, and it is impossible to select a higher level. Please enter the skin care page of the FLOSSOM applet and reset the level of the skinpilot.
	The level of the skinpilot cannot be adjusted through the power/level selection button on the skinpilot when the level indicator on the skinpilot is breathing in blue when using the FLOSSOM applet to control current output. To adjust the level on the skinpilot, please disconnect from the FLOSSOM applet first.
	The skinpilot has failed. Please contact the customer service of the authorized store where the product was purchased or call our customer service hotline (400-848-0056) directly.
Only minutes of use even after a full charge	After skin care, the product can will power off automatically.
	The skinpilot powers off automatically because of low battery level. Please charge it before continued use.
	The battery will expire soon. Please contact the customer service of the authorized store where the product was purchased or call our customer service hotline (400-848-0056) directly.
Strong current stimulation	Please select a lower level.
	If the current stimulation in certain area is strong, select a lower current level on the FLOSSOM applet.
Insufficient or no current stimulation	Select a higher level.
	Make the flexible conductive film fabric fit the skin well.
	Before use, please apply mist spray or toner onto the flexible conductive film fabric.
	It is normal because the muscles in different facial areas may have different sensibility for the current stimulation.
Red level indicator in use	The flexible conductive film fabric does not fit the face well, or its magnetic connection ends have a poor contact with the skinpilot. Please adjust for good contact or fit. When the level indicator is breathing in green, press the power/level selection button for continued skin care.
Failure of connection of the magnetic connection ends on the two sides to the skinpilot	The skinpilot is worn reversely. Please take it off and wear it on the head like a hair band, making sure that its left and right sides are behind the ears and that the power/level selection button is on the right hand side. Refer to the figure on the right part in STEP2 on page 6 in this manual.
	The skinpilot has failed. Please contact the customer service of the authorized store where the product was purchased or call our customer service hotline (400-848-0056) directly.

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Q & A

Q: Can the product be used in a bathroom?

A: The product cannot be used in a bathroom (or a bathtub).

Q: Can the product be used in sleep?

A: The product cannot be used in sleep.

Q: Why do left half and right half of the face have different feelings in use?

A: Left half and right half of the face are usually asymmetric. Please adjust the current levels for different facial areas on the FLOSSOM applet.

Q: Why is there toothache in use?

A: It may happen to users receiving dental treatment, with dental implant or tooth filling. In the case of toothache, please consult a doctor immediately.

Q: Can the product be used immediately after sunburn or shaving?

A: Using the product immediately after sunburn or shaving may result in inflammation. In the case of any skin problem, please stop using the product immediately and consult a doctor.

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Disposal of built-in rechargeable lithium battery

How to take out the battery:

1. Before disposal the product, please take out its battery.
2. Before taking out the battery, please power off the product for safety.
3. For any help, please contact the customer service of our flagship store.

How to dispose the used battery:

Send the used battery to a used battery recycling organization for professional treatment.

⚠️ Notes:

- The operation instructions above apply to disposal of the skinpilot only. The built-in lithium battery should be disposed according to the local laws and regulations about recycling.
- The operation instructions above do not apply to normal repair of the product.
- If the skinpilot fails due to unauthorized disassembling, we will not provide the normal after-sales service for it.

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- ① Open the bottom shell along the buckling directions on the sides with help of a small screwdriver and loosen the buckles.

底壳: Bottom shell

- ② Take out the sponge.

海绵: Sponge

- ③ Remove the screws.

螺钉: Screws

- ④ Take out the battery pressing plate.

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- ⑤ Separate the parts on the top shell.

- ⑥ Open the mainboard on the right, pull out the battery terminal and take out the battery.

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
Hazardous substances or elements

	Hazardous Substances/Elements					
Component	Pb	Hg	Cd	Cr(VI)	PBB	PBDE
USB cable	○	○	○	○	○	○
PCB	○	○	○	○	○	○
Battery	○	○	○	○	○	○
Plastic parts	○	○	○	○	○	○
Hardware	x	○	○	○	○	○
○: The content of a hazardous substance in all homogeneous materials of the component is less than the limit in GB/T26572. x: The content of a hazardous substance in any homogeneous material of the component is more than the limit in GB/T26572. (However, the component with a “x” symbol meets the RoHS requirement and is exempted.)						

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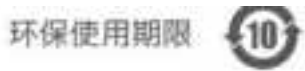Environmental life

This sign refers to an environmental life of 10 years for an electronic product. In this period, hazardous or toxic substances in an electronic product used normally will not leak or mutate and use of the electronic product will not cause serious environmental pollution, personal injuries or property losses

Please read this manual.

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

RF Exposure Information

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator and your body.

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After-sales information

After-sale service and warranty service

- The product (except for its consumables) has a warranty service of two years following its date of purchase if it is used normally.
- If the product needs to be repaired out of its warranty period, its user will need to bear the transport cost.
- For any fault out of the warranty period, the actual repair cost will be charged.
- The repairs in the warranty period may result in expenses. Refer to the content under “Circumstances not covered by free repair”.
- The warranty service applies to the Chinese mainland only.

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In any of the following circumstances, we will not provide free repair for the product.

- The product fails or is damaged due to improper use, repair or refit;
- The product fails or is damaged due to toppling over or fall in use;
- The product fails or is damaged due to a fire, earthquake, flood, or natural disaster, public disaster, or abnormal voltage;
- The product fails or is damaged due to an abnormality of any electric appliance connected to it;
- Battery, any consumable or any other part of the product has been replaced;
- The product fails or is damaged due use in a non-household condition (commercial application needing long-time use);
- The product fails or is damaged due to use on a vehicle or ship;

- The product fails or is damaged due to transport, moving or drop after being purchased.

The battery is a rechargeable battery is a consumable and is not covered by the warranty service. For its replacement, please contact the customer service of our flagship store or call our customer service hotline (400-848-0056) directly.

Any circumstance not mentioned on the warranty card

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Contact information for after-sales service and warranty service

For a product purchased from our flagship store or any store or dealer authorized by us, you can contact the customer service of any store above or call our customer service hotline (400-848-0056) directly and mail it to the address designated by us to be repaired for free. If the product cannot be repaired for free because it's given to another person as a gift, please contact the customer service of or our flagship store or call our customer service hotline (400-848-0056) directly.

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Contact information

Retailer: Flossom (Shanghai) Beauty Technology Co., Ltd.

Address: 8F, Building 2, 1919 Baziqiao Road, Fengxian District, Shanghai

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Warranty card

Warranty Card	
Product name	FLOSSOM Wearable Intelligent Beauty Apparatus
Model	WE200-M
Customer	
Contact information	
Channel/store the product is purchased from	
Reason of repair	
Date of filling	