



www.mortara.com

REF 9515-210-50-ENG Rev A1

Surveyor S2 Telemetry Device

USER MANUAL

Manufactured by Mortara Instrument, Inc., Milwaukee, Wisconsin U.S.A.



CAUTION: *United States federal law restricts this device to sale by or on the order of a physician.*



Copyright © 2018
by Mortara Instrument, Inc.
7865 N. 86th Street
Milwaukee, Wisconsin 53224
E-mail: techsupport@mortara.com

This document contains confidential information that belongs to Mortara Instrument, Inc. No part of this document may be transmitted, reproduced, used, or disclosed outside of the receiving organization without the express written consent of Mortara Instrument, Inc. V1.0.0 (2018-02)

Mortara is a registered trademark of Mortara Instrument, Inc. Surveyor™ is a trademark of Mortara Instrument, Inc. All other trademarks and registered trademarks are the property of their respective owners.

Product License

This product and its accompanying software are subject to Mortara Instrument's End User License Agreement. By using this product and its accompanying software, you agree to the terms and conditions specified therein.

Open-Source Software

Mortara Instrument uses a variety of software and components created by the open-source community. The use of open-source software and components is essential to creating robust, high quality products, as these have been under constant scrutiny and review by the open-source community. Upon explicit written request, Mortara Instrument will make available to customers the source code for open-source components used in this product within three years after the purchase of the product. Mortara Instrument may charge a fee for delivering such source code. Permission to use, copy, modify, and/or distribute this software for any purpose with or without fee is hereby granted, provided that the above copyright notice and this permission notice appear in all copies.

Additional Notices

Portions of this product incorporate software and components that are copyrighted by their respective owners including:

Cypress WICED Studio

*Copyright (c) 2017
All rights reserved.*

SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING (BUT NOT LIMITED TO) THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT OWNER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION), HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY; WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE); ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

TABLE OF CONTENTS

1.	GENERAL STATEMENTS.....	7
	Technical Support and Service	7
2.	NOTICES.....	8
	Manufacturer's Responsibility.....	8
	Responsibility of the Customer	8
	Equipment Identification	8
	Copyright and Trademark Notices.....	8
	Other Important Information	8
3.	WARRANTY INFORMATION.....	9
	Your Mortara Warranty	9
4.	USER SAFETY INFORMATION	10
	Safety Regulations.....	10
	Warnings.....	10
	Cautions.....	14
5.	EQUIPMENT SYMBOLS AND MARKINGS.....	15
6.	GENERAL CARE.....	16
	Precautions for S2.....	16
	Precautions for Li-Ion Battery Charger and Cord.....	16
	Inspection Prior to Clinical Use.....	16
	Cleaning and Disinfecting	17
	Maintenance	19
	Battery life	20
	Product life	20
	Decommissioning and Disposal	20
7.	ELECTROMAGNETIC COMPATIBILITY (EMC)	21
	Guidance and Manufacturer's Declaration: Electromagnetic Emissions.....	22
	Guidance and Manufacturer's Declaration: Electromagnetic Immunity	22
	Guidance and Manufacturer's Declaration: Electromagnetic Immunity	23
	Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and the Equipment.....	24
	USA Radio Regulations.....	25
8.	INTRODUCTION.....	26
	General Information	26
	Indications For Use.....	26
	System Overview.....	27

Operating controls and indicators	28
Carrying Pouches.....	29
Battery Charger.....	30
9. UNPACKING AND SET UP	31
Checking Contents.....	31
Battery Installation	32
10. PATIENT PREPARATION FOR QUALITY ECG	33
Quality ECG Data Acquisition	33
Skin Preparation	33
Electrode Placement	34
Electrode Locations for 4-Wire Cable	35
Pacemaker Patients	35
Checking ECG Electrode and Lead Wire Signal Quality	36
11. OPERATION.....	37
Powering On the S2.....	37
Screen Time-Out & Reactivation	37
Powering Off the S2	37
Main Menu Screen.....	38
Waveform Screen	46
Parameters Screen.....	47
Demographics Screen	48
12. CONFIGURING THE S2	49
Host Settings.....	51
Network Settings	51
Language Settings.....	53
Device Diagnostics	53
S2 Version	53
13. ECG MONITORING.....	54
Power-On.....	54
Patient Prep.....	54
Patient Cable Hookup.....	54
Confirm Good Signal.....	54
Patient Sessions	55
15. PRODUCT SPECIFICATIONS.....	56
General	56
Environmental	56
Power Requirements & Battery	57
ECG Specifications.....	57
Wireless Network Specifications.....	58
16. TROUBLESHOOTING	59

Power and Battery	59
Display and Touchscreen.....	59
ECG Trace	59
Network Transmission.....	60
17. REORDERING ACCESSORIES & CONSUMABLES	62
Power accessories	62
Mounts / holders	62
ECG accessories	62
18. APPLIED STANDARDS.....	63

1. GENERAL STATEMENTS

Technical Support and Service

Mortara Instrument, Inc.

Global Headquarters

7865 North 86th Street Milwaukee,
WI 53224 U.S.A.

Tel: 414.354.1600
Tel: 800.231.7437
Fax: 414.354.4760
Internet: <http://www.mortara.com>

Mortara Instrument Europe

Customer Care:

IDA Business Park – Johnstown
Navan, County Meath Ireland
Tel.: +353 46 90 67790
Email: europe@welchallyn.com

Service Center:

Centaurusweg 123, 5015 TC
Tilburg, The Netherlands
RMA
Tel: +353 46 9067790
eme.servicecenter@welchallyn.com

Tech. Support

Tel: +353 46 9067790
eme.techsupport@welchallyn.com

Mortara Instrument Australia

Head Office

Suite 4.01, 2-4 Lyonpark Road
Macquarie Park, Sydney
NSW 2113 Australia

Tel: 1800 650 083
Fax: +61 2 9562 0982

Mortara Instrument France

15, rue Léon Gaumont
ZAC de la Pentecôte
44700 ORVAULT
France

Tel: +33 253355530
+33 157324994
Fax: +33 253358539

Mortara Instrument Germany

Hofgartenstraße 16
72379 Hechingen
Deutschland

Tel.: +49 (0)7471984114-0
Fax: +49 (0)747984114-90
E-Mail: info@welchallyn.com

Mortara Instrument Italy

Via Cimarosa, 103/105
40033 Casalecchio di Reno (Bologna)
Italia

Tel: +39 051 2987811
Fax: +39 051 6133582
E-mail: clienti.mortarait@welchallyn.com

Mortara Instrument Netherlands

“Amerika” Gebouw– 7e verdieping
Hoogoorddreef 15
1101 BA Amsterdam
Nederland

Tel.: 020 206 1360
E-mail: Infonl@welchallyn.com

Mortara Instrument UK

Clinitron House, Excelsior Road
Ashby de la Zouch
Leicester LE65 1JG

Tel: 0207 365 6780
Fax: 0207 365 9694

2. NOTICES

Manufacturer's Responsibility

Mortara Instrument, Inc. (Mortara) is responsible for the effects on safety and performance of the Surveyor™ S2 telemetry device



WARNING: Only Mortara authorized service providers should perform servicing of the S2 to ensure that the correct maintenance and calibration procedures are followed and that the S2 returns to proper operation. Refer to Section 1 for a list of technical support and service providers.

- The telemetry device is used in accordance with the instructions for use.
- The telemetry device is correctly maintained according to the standards authorized by Mortara using original spare parts.
- The telemetry device is used with original accessories and supplies that are in compliance with the standard specifications described in this manual.
- The electrical installation of the relevant room complies with the requirements of appropriate regulations.

Responsibility of the Customer

The user of this telemetry device is responsible for ensuring the implementation of a satisfactory maintenance schedule. Failure to do so may cause undue failure and possible health hazards. This manual must be kept in a safe place to prevent its deterioration and/or alteration. The user and Mortara authorized personnel must have access to this manual at any time. The user of this telemetry device must periodically check the accessories, their functionality and integrity.

Equipment Identification

Mortara equipment is identified by a serial and reference number on the back of the telemetry device. Care should be taken so that these numbers are not defaced.

Copyright and Trademark Notices

This document contains information that is protected by copyright. All rights are reserved. No part of this document may be photocopied, reproduced, or translated into another language without prior written consent of Mortara.

Other Important Information

The information in this document is subject to change without notice.

Mortara makes no warranty of any kind with regard to this material including, but not limited to, implied warranties of merchantability and fitness for a particular purpose. Mortara assumes no responsibility for any errors or omissions that may appear in this document. Mortara makes no commitment to update or to keep current the information contained in this document.

3. WARRANTY INFORMATION

Your Mortara Warranty

MORTARA INSTRUMENT, INC. (hereafter referred to as “Mortara”) warrants that components within Mortara products (hereafter referred to as “Product/s”) will be free from defects in workmanship and materials for the number of years specified on documentation accompanying the product, or previously agreed to by the purchaser and Mortara, or if not otherwise noted, for a period of twelve (12) months from the date of shipment.

Consumable, disposable or single use products such as are warranted to be free from defects in workmanship and materials for a period of 90 days from the date of shipment or the date of first use, whichever is sooner.

Reusable product such as, but not limited to, BATTERIES, PATIENT CABLES, LEAD WIRES, MAGNETIC STORAGE MEDIUMS, CARRY CASES or MOUNTS, are warranted to be free from defects in workmanship and materials for a period of 90 days. This warranty does not apply to damage to the Product/s caused by any or all of the following circumstances or conditions:

- a) Freight damage;
- b) Supplies, accessories and internal parts NOT approved by Mortara;
- c) Misapplication, misuse, abuse, and/or failure to follow the Product/s instruction sheets and/or information guides;
- d) Accident;
- e) A disaster affecting the Product/s;
- f) Alterations and/or modifications to the Product/s not authorized by Mortara;
- g) Other events outside of Mortara’s reasonable control or not arising under normal operating conditions.

THE REMEDY UNDER THIS WARRANTY IS LIMITED TO THE REPAIR OR REPLACEMENT WITHOUT CHARGE FOR LABOR OR MATERIALS, OR ANY PRODUCT/S FOUND UPON EXAMINATION BY MORTARA TO HAVE BEEN DEFECTIVE. This remedy shall be conditioned upon receipt of notice by Mortara of any alleged defects promptly after discovery thereof within the warranty period. Mortara’s obligations under the foregoing warranty will further be conditioned upon the assumption by the purchaser of the Product/s (i) of all carrier charges with respect to any Product/s returned to Mortara’s principal place or any other place as specifically designated by Mortara or an authorized distributor or representative of Mortara, and (ii) all risk of loss in transit. It is expressly agreed that the liability of Mortara is limited and that Mortara does not function as an insurer. A purchaser of a Product/s, by its acceptance and purchase thereof, acknowledges and agrees that Mortara is not liable for loss, harm, or damage due directly or indirectly to an occurrence or consequence there from relating to the Product/s. If Mortara should be found liable to anyone under any theory (except the expressed warranty set forth herein) for loss, harm, or damage, the liability of Mortara shall be limited to the lesser of the actual loss, harm, or damage, or the original purchase price of the Product/s when sold.

EXCEPT AS SET FORTH HEREIN WITH RESPECT TO REIMBURSEMENT OF LABOR CHARGES, A PURCHASER’S SOLE EXCLUSIVE REMEDY AGAINST MORTARA FOR CLAIMS RELATING TO THE PRODUCT/S FOR ANY AND ALL LOSSES AND DAMAGES RESULTING FROM ANY CAUSE SHALL BE THE REPAIR OR REPLACEMENT OF DEFECTIVE PRODUCT/S TO THE EXTENT THAT THE DEFECT IS NOTICED AND MORTARA IS NOTIFIED WITHIN THE WARRANTY PERIOD. IN NO EVENT, INCLUDING THE CLAIM FOR NEGLIGENCE, SHALL MORTARA BE LIABLE FOR INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES, OR FOR ANY OTHER LOSS, DAMAGE, OR EXPENSE OF ANY KIND, INCLUDING LOSS OF PROFITS, WHETHER UNDER TORT, NEGLIGENCE OR STRICT LIABILITY THEORIES OF LAW, OR OTHERWISE. THIS WARRANTY IS EXPRESSLY IN LIEU OF ANY OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND THE WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE.

4. USER SAFETY INFORMATION



WARNING: Means there is the possibility of personal injury to you or others.



CAUTION: Means there is the possibility of damage to the telemetry device.

NOTE: Provides information to further assist in the use of the telemetry device.

***NOTE:** This manual may contain screen shots and pictures. Any screen shots and pictures are provided for reference only and are not intended to convey actual operating techniques. Consult the actual screen in the host language for specific wording.*

Safety Regulations

- The Surveyor S2 (henceforth referred to as either Surveyor S2 or S2) is a medical telemetry device.
- The S2 accessories are  labeled, according to applicable standards.
- The S2 with all accessories that have a physical or logical connection with it, forms part of a Medical Electrical System.
- The S2 complies with various safety and performance regulations as mentioned in this manual (Applied Standards).



Warnings

- This manual gives important information about the use and safety of this telemetry device. Deviating from operating procedures, misuse or misapplication of the telemetry device, or ignoring specifications and recommendations could result in increased risk of harm to users, patients and bystanders, or damage to the telemetry device.
- Users are expected to be licensed clinical professionals knowledgeable about medical procedures and patient care, and adequately trained in the use of this telemetry device. The S2 telemetry device captures and transmits data reflecting a patient's physiological condition to a Monitoring Node that when reviewed by a trained physician or clinician can be useful in determining a diagnosis.
- Before attempting to use this device for clinical applications, the operator must read and understand the contents of the user manual and other accompanying documents. Inadequate knowledge or training could result in increased risk of harm to users, patients and bystanders, or damage to the telemetry device. Contact Mortara for additional training options.
- Operation of the equipment beyond its specified ranges, or beyond normal physiological conditions of human subjects, may cause inaccurate results.
- A possible explosion hazard exists. Do not use the device in the presence of a flammable anesthetic mixture. Do not mount any part of the device closer than 25 cm from outlets of flammable gases, including oxygen.
- For proper operation and the safety of users or patients and bystanders, equipment and accessories must be connected only as described in this manual.

- Repairs and modification must be made by authorized and trained technical personnel. Unauthorized modifications and repairs will void the S2 warranty and may pose a danger to patients and users.
- If additional devices beyond the S2 are connected to the patient, leakage currents could add up and should be accounted for.
- The S2, as all medical equipment or systems, requires special precautions regarding EMC, and should be installed and put into service according to the EMC information provided in the installation procedure to obtain a sufficient degree of immunity as well as not to create disturbance to other equipment. Refer to the specific EMC instructions in this manual.
- The quality of the signal produced by the device may be adversely affected by the use of other medical equipment, including but not limited to defibrillators, electrosurgery equipment and ultrasound machines. Do not use the S2 system in the presence of imaging equipment such as magnetic resonance imaging (MRI) and tomography systems. Simultaneous operation may damage the device or lead to erroneous results.
- There is no known safety hazard if other equipment, such as pacemakers or other stimulators, is used simultaneously with the device; however, disturbance to the signal may occur.
- Portable and mobile RF communications equipment may affect medical electrical equipment or systems as well as the S2 and its accessories. Do not operate the S2 near sources of high frequency emissions (e.g. microwaves). Unauthorized wireless devices, such as personal access points and WLAN hotspots including those available on personal smartphones, may also interfere with the operation of the system.
- A telemetry device in conjunction with a Monitoring Node is not intended to replace clinical assessments. It is important that a qualified individual regularly supervise the patient.
- The S2 is restricted to use on one patient at a time.

Power Supply Warnings

- Only use the recommended batteries. Use of alternate batteries may damage the device or cause other hazards. Do not use a battery that is physically damaged.
- The S2 is a battery operated device that transmits data reflecting a patient's physiological condition to a receiving device. During operation failure, data transmission and LCD information will cease to occur. In the case of battery depletion, replace the batteries on the device to resume function. In mission critical conditions, it is advisable to have a backup device available.
- When charging the Li Ion battery, ensure that the electrical installation complies with local safety requirements for the use environment.
- Regularly check all cables for damage and proper connection. Do not use equipment with a damaged cable.
- The S2 contains an internal battery. The following precautions should be taken regarding the internal battery:
 - Do not immerse in water.
 - Do not heat or throw in fire.
 - Do not leave in conditions over 60 °C.
 - Do not crush or drop.
 - Only use the approved batteries.
 - Follow the instructions in the disposal section of this manual when taken out of service.
- The S2 rechargeable battery must be initially fully charged prior to use.

- When the S2 initially powers on, the screen will illuminate if the battery is installed properly and charged. Remove the S2 from service and contact Mortara if the screen does not activate when a new or fully charged battery is initially installed.
- Batteries are known to leak their contents when stored for an extended period of time in unused equipment. Always remove the battery after completing operating the telemetry device. Always place rechargeable batteries in the battery charger when not in use. This ensures that the battery is recharged for the next time the telemetry device is operated.
- There is a potential pinch hazard when applying the battery compartment cover to the device that could result in minor injury. Care should be taken to avoid entrapment of fingers when performing this operation.

Accessories, Cables, and External Connections Warnings

- The S2 is designed to meet applicable specifications when using Mortara-approved patient cables and accessories. Use of non-approved cables and accessories may result in reduced performance or electromagnetic interference, and may pose possible patient and user safety concerns.
- Do not use excessive force on any of the connection cables and handle all accessories with care.
- Conductive parts of the ECG patient cables, electrodes, and associated connections of type CF applied parts, including the neutral conductor of the patient cable and electrode should not come into contact with other conductive parts including earth ground.
- To avoid the possibility of serious injury or death during patient defibrillation, do not come into contact with telemetry device or patient cables.
- Accessories may be provided with separate user manuals. Read these manuals thoroughly and refer to them for specific functions. It is recommended to keep all manuals together.
- To avoid potential for spread of disease or infection, single-use components and accessories (e.g., electrodes, , pouches, etc.) must not be reused. To maintain safety and effectiveness, ECG electrodes and sensors must not be used beyond their expiration date or useful life.
- All accessories including cables, connectors and other patient-applied parts supplied with the S2 do NOT contain any latex. If the patient develops an allergic reaction or rash, immediately remove the accessory and inform Mortara.

Defibrillation and Electrosurgery Warnings

- The S2 has not been designed for use together with electrosurgery equipment.
- The S2 is defibrillator protected in compliance with IEC 60601-2-27 standards if used with Mortara-approved patient cables. Defibrillation while using non-approved patient cables may damage the device beyond repair and cause a safety hazard to the patient.

ECG Warnings

- Excessive patient movement could interfere with the operation of the system.
- Proper clinical procedure must be employed to prep the electrode and sensor sites, and to monitor the patient for excessive skin irritation, inflammation, or other adverse reactions. Electrodes and other sensors that are intended for short-term use should be removed from the patient promptly following use.

- When using the 4-wire ECG cables, it is not possible to acquire a 12-lead ECG with the S2.
- The S2 allocates no less than 10 mm (1mV) of vertical space to each ECG waveform displayed.
- Do not use sensors and/or cables showing signs of physical damage, as these may produce erroneous measurements.



Cautions


- This device must be installed as part of a system in conjunction with the MONITORING NODE and in accordance to guidance and minimum characteristics per requirements provided by Mortara for deployment of the system on the hospital/clinic's IT network. Refer to those requirements as well as Manufacturer Disclosure Statement for Medical Device Security (MDS2) statements provided by Mortara before deploying and using the system.
- The device and patient cable should be cleaned between each use. Cleaning must be performed with the system turned off and battery removed. Let all parts dry well before turning the power back on.
- Prevent liquids from penetrating the system, components, or the monitor. Do not spray the system with liquid cleaning agents. Do not allow the system, components or accessories to become in contact with unknown substances which may compromise its mechanical or electrical integrity. If liquids have penetrated the system, open by authorized personnel for inspection and let dry completely.
- Do not attempt to clean the telemetry device or patient cables by submersing into a liquid, autoclaving, or steam cleaning as this may damage equipment or reduce its usable life. Wipe the exterior surfaces with a warm water and mild detergent solution and then dry with a clean cloth. Use of unspecified cleaning/disinfecting agents, failure to follow recommended procedures, or contact with unspecified materials could result in increased risk of harm to users, patients and bystanders, or damage to the telemetry device.
- No user-serviceable parts inside. Damaged or suspected inoperative equipment must be immediately removed from use and must be checked/repared by authorized service personnel prior to continued use.
- The S2 accommodates rechargeable batteries. If the rechargeable battery appears to become defective, refer to Mortara Technical Support.
- Do not pull or stretch patient cables as this could result in mechanical and/or electrical failures. Patient cables should be stored off of the floor away from bedrails and wheels to avoid cable damage. Roll the patient cables into a loose loop prior to storage.
- When necessary, dispose of the telemetry device, its components and accessories (e.g., batteries, cables, electrodes), and/or packing materials in accordance with local regulations.
- Check that all operating and environment conditions such as ambient temperature meet the device's specifications. Allow the device to stabilize within its intended operating environment for a minimum of two hours prior to use. Do not use the device outside the operating and environment conditions specified (section 15).
- Do not exert excessive pressure on the touchscreen LCD or use a sharp or hard object with it. Excessive pressure may permanently damage the display.
- The device is UL classified:



UL-classified device in the USA and Canada.

- Upon request, Mortara can supply a service manual that includes test instructions as well as a list of spare parts that must be used with the S2.

5. EQUIPMENT SYMBOLS AND MARKINGS

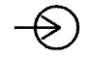 Indicates compliance to applicable European Union directives	 Do not dispose as unsorted municipal waste. Per European Union Directive 2002/96, requires separate handling for waste disposal according to national requirements
	 Defibrillator-proof type CF applied part
 Warning/Caution	 Consult accompanying documents
 Catalog number for relevant Mortara part	 Serial number
 This end up	 Keep away from sunlight
 Fragile, handle with care	 Keep dry
 Temperature limitation	 Direct Current
 Power On/Wake Screen	 This icon indicates the status of the WLAN connection
 Latex free Does not contain natural rubber latex.	ECG  Patient Cable Input

6. GENERAL CARE

Precautions for S2

- Power off and remove the battery from the S2 before inspecting or cleaning.
- Protect the S2 from water and other liquids.
- Never immerse the S2 in water or other fluids.
- Do not drop the S2 or subject to shock and/or vibration.
- Do not use organic solvents, ammonia-based solutions, or abrasive cleaning agents that may damage equipment surfaces.

Precautions for Li-Ion Battery Charger and Cord

- Remove the AC power from a Li-Ion Battery Charger before inspecting or cleaning.
- Protect the AC power cord and Li-Ion Battery Chargers from water or other liquids.
- Never immerse the AC power cords or Li-Ion Battery Chargers in water or other fluids.
- Do not subject the Li-Ion Battery Chargers to shock and/or vibration.
- Do not use organic solvents, ammonia-based solutions, or abrasive cleaning agents that may damage equipment surfaces.

Inspection Prior to Clinical Use

Inspect your equipment prior to clinical operation. Do not use the equipment and contact an authorized service representative for servicing if there are concerns about integrity of the system.

- Verify that all cables and connectors are securely seated.
- Check the case and chassis of the S2, AC power cord, and Li-Ion Battery Charger for any visible damage.
- Inspect the S2 touchscreen and controls for proper function and appearance.
- Inspect the Li-Ion Battery Charger indicator LEDs for proper operation and indication of battery charging and AC connection.
- Inspect patient accessories for any visual damage.
- S2 patient input connector - Verify the pins on the patient input connector are all present and are not bent or damaged in any way. The recessed area for the patient connector should be free from debris and clean. Use compressed air to remove any debris that has entered into the connector area.
- S2 display – Verify there are no deep scratches or physical damage to the S2 telemetry device's display. Inspect the display bezel to ensure it is firmly adhered to the device housing. Contact Mortara Technical Support if the display or display bezel require replacement.
- S2 battery door – Verify the S2's battery door for proper opening and closing. Inspect the plastic door assembly for signs of excessive wear or cracking, including the door seal to prevent fluid ingress. Replace the battery door assembly if necessary.

- S2 Battery Compartment – Inspect the S2’s battery spring contacts and the battery door latching mechanism for signs of excessive wear. If the battery compartment has been damaged, contact Mortara Technical Support for assistance.
- Battery Charger bay – Inspect the Li-Ion Battery Charger’s contacts and mechanisms for signs of excessive wear.
- Rechargeable Li-Ion Battery – Follow the instructions and note the cautions labeled on the rechargeable battery. Contact Mortara for replacement.
- Device Labeling – Inspect the device labeling for signs of wear and legibility. If the labeling is no longer clear and legible, contact Mortara Technical Support for assistance.

Cleaning and Disinfecting

This section describes the procedures for cleaning and disinfecting the S2, its sensors and accessories.



WARNING: Follow these instructions to clean and disinfect the S2 and its accessories. Improper cleaning may cause damage that is not immediately apparent, leading to possible safety hazards, device malfunction, and/or spread of infectious agents between persons.

Disinfecting agents

The S2 is compatible with the following disinfectants:

- Clorox Healthcare® Bleach Germicidal Wipes (use according to instructions on product label), or
- a soft, lint-free cloth dampened with a solution of sodium hypochlorite (10% household bleach and water solution) minimum 1:500 dilution (minimum 100 ppm free chlorine) and maximum 1:10 dilution as recommended by the APIC Guidelines for Selection and Use of Disinfectants.

Cleaning and disinfecting the S2 monitor



WARNING: Do not immerse the S2 in water or any other fluids. The S2 is not designed to be immersed in liquid and doing so may result in liquid entering the device leading to possible safety hazards and/or device malfunction.



CAUTION: Do not steam autoclave, gas sterilize, or irradiate the S2 as these may result in damage to the device.



WARNING: Ensure the battery door is securely in place when cleaning the S2 to avoid risk of liquid entering into the device which may lead to a possible hazard and/or device malfunction.

To clean the S2:

1. Switch off the S2. If on, press and hold the Push-Button button on the right side of the device to turn on the display screen. Press the configuration key on the touch screen, then select [Shutdown Device].
2. Disconnect patient cables from the S2.

3. Thoroughly wipe the surface of the S2 with a clean, lint-free cloth dampened with water for general cleaning, or use one of the above recommended agents for disinfection.



WARNING: Do not oversaturate the cleaning cloth. Liquid pooling on the device may enter into the device possibly leading to a safety hazard and/or device malfunction.

4. Dry the device with a clean, soft, dry, lint-free cloth.

Cleaning sensors and other accessories

The S2 is compatible with a number of sensors and accessories, each with unique cleaning needs. Follow the cleaning instructions provided in the directions for use shipped with those items.



WARNING: Do not reuse sensors or accessories indicated as single-patient use; as this may facilitate the spread of infectious agents between persons.



WARNING: Always clean and/or disinfect reusable sensors and accessories between patients to reduce the risk of spreading infectious agents between persons.

Maintenance

The following table shows the recommended maintenance procedures for the S2, patient accessories, and Li-Ion Battery Charger. The S2 should be serviced once a year by a Mortara authorized service technician; however, it is good practice to periodically ensure the telemetry device is in proper working order. This can be performed by a clinician or biomed at the hospital or healthcare delivery organization familiar with the S2 telemetry device, ECG signal acquisition, as well as general maintenance/calibration of biomedical equipment.

To accomplish these steps in their entirety and verify the correct operation of the system, appropriate patient simulators or other equipment may be required. Refer to the service manual for further details.

Functionality	Procedure
Mechanical Integrity	Check for cracks, abrasive edges and other signs of damage.
ECG Telemetry	- Connect ECG leads to Patient Simulator. - Start a new monitoring session. - Verify that waveforms for all leads are properly shown on the LCD.
Li-Ion Battery Charger	- Check for cracks, abrasive edges and other signs of damage. - Check that all connectors and AC cord length is unbroken and smooth along its length. - Verify proper LED indicators during battery charging.
ECG Cables	Approved Cleaning Agents - Enzymatic detergent such as ENZOL (US) or CIDEZYME (outside the US). - Distilled water. - Disinfectant solution (such as CIDEX OPA, or a 10% solution of household bleach (5.25% sodium hypochlorite) in distilled water). - Soft, lint-free cloths and/or soft-bristled brushes. - Protective gloves and eyewear. Procedure - Remove the battery from the telemetry device. - Put on gloves and protective eyewear. - Prepare the detergent according to the manufacturer's instructions, and also the disinfectant solution, in separate containers. - Apply detergent to product using a soft, lint-free cloth. If material is dried on, allow to sit for 1 minute. Do not immerse cable ends or lead wires in liquid as it can cause corrosion. - Wipe smooth surfaces with the cloth. - Use a soft-bristle brush on visibly soiled areas and irregular surfaces. - Remove detergent from product using cloth dampened in distilled water. - Repeat as necessary. - Apply disinfectant solution on affected area using a soft cloth. Allow product to sit for 5 minutes. - Wipe excess solution and clean product again with cloth dampened in distilled water. - Allow 2 hours for drying.



WARNING: Only Mortara authorized service providers should perform servicing of the S2 to ensure that the correct maintenance and calibration procedures are followed and that the S2 is returned to proper operation. Refer to Section 1 for a list of technical support and service providers.



WARNING: Do not use the S2 telemetry device and/or its accessories and parts beyond their product life or past their expiration date, since performance of these may be degraded leading to inaccurate or misleading results and/or materials may break down over time releasing substances that may be harmful.

Battery life

The S2 may be powered only using a (1) Rechargeable Li-Ion Battery Pack.. As the battery pack ages, its ability to store charge is gradually reduced. If the battery pack no longer provides the run time needed, it should be replaced. (See table below in Chapter 15 – Power Requirements & Battery for operation time)



WARNING: Remove the batteries if the S2 will not be used for an extended period to avoid possible leakage of harmful substances from the battery.



WARNING: Use only APPROVED BATTERIES and BATTERY CHARGER as listed in the Accessories section. Use of unapproved batteries may cause a hazard and/or damage the S2, and will void the warranty.

NOTE: Excessive wireless network traffic, signal strength, and network dropouts may affect the battery life.

Product life

The S2 has a defined product life of 5 years excluding accessories, cables and batteries. As required, product service, accessories and spare parts are available through Mortara or its authorized partners. Using the telemetry device or its accessories and components beyond their defined life may lead to damage to the equipment or a safety hazard to the user.

Decommissioning and Disposal

Dispose of the telemetry device, its components and accessories (e.g., batteries, cables, electrodes), and/or packing materials in accordance with local regulations. Do NOT incinerate the battery.

7. ELECTROMAGNETIC COMPATABILITY (EMC)

When using the telemetry device, assess the electromagnetic environment affected by surrounding devices.

An electronic device may either generate or receive electromagnetic interference. Testing for electromagnetic compatibility (EMC) has been performed on the telemetry device according to the applicable international standards.

The telemetry device should not be used adjacent to or stacked with other equipment. If the telemetry device is used in this manner, verify the telemetry device operates in an acceptable manner in the configuration in which it will be used.

Fixed, portable, and mobile radio frequency communications equipment may affect the performance of medical equipment. See the appropriate EMC table for recommended separation distances between the radio equipment and the telemetry device.

The use of accessories, transducers, and cables other than those specified by Mortara may result in increased emissions or decreased immunity of the equipment.

Guidance and Manufacturer's Declaration: Electromagnetic Emissions

The equipment is intended for use in the electromagnetic environment specified in the table below. The customer or the user of the equipment should ensure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment: Guidance
RF Emissions CISPR 11	Group 2	The S2 must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
RF Emissions CISPR 11	Class A	
Harmonic Emissions IEC 61000-3-2	Not Applicable	
Voltage Fluctuations/ Flicker Emissions IEC 61000-3-3	Not Applicable	

Guidance and Manufacturer's Declaration: Electromagnetic Immunity

The equipment is intended for use in the electromagnetic environment specified in the table below. The customer or the user of the equipment should ensure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment: Guidance
Electrostatic discharge (ESD) EN 61000-4-2	+/- 8 kV contact +/- 15 kV air	+/- 8 kV contact +/- 15 kV air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst EN 61000-4-4	+/- 2 kV for power supply lines +/- 1 kV for input/output lines	+/- 2 kV for power supply lines +/- 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	+/- 1 kV differential mode +/- 2 kV common mode	+/- 1 kV differential mode +/- 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage fluctuations and Interruptions	<5% UT for 0.5 cycles 40% UT for 5 cycles 70% UT for 25 cycles <5% UT for 5s	<5% UT for 0.5 cycles 40% UT for 5 cycles 70% UT for 25 cycles <5% UT for 5s	Note that monitoring is interrupted at the level "< 5% UT for 5s", but equipment remains safe (as specified in EN 60601-1-2).
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: UT is the AC Mains voltage prior to application of the test level.

Guidance and Manufacturer's Declaration: Electromagnetic Immunity

The equipment is intended for use in the electromagnetic environment specified in the table below. The customer or the user of the equipment should ensure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment: Guidance
Conducted RF EN 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms 150 kHz to 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the equipment, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2 \sqrt{P}$ $d = 1.2 \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m 80 MHz to 2.7 GHz	

a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radios, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the equipment is used exceeds the applicable RF compliance level above, the equipment should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the equipment.

b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [3] V/m.

Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and the Equipment

The equipment is intended for use in the electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the equipment can help to prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the equipment as recommended in the table below, according to the maximum output power of the communications equipment.

Rated Maximum Output Power of Transmitter (W)	Separation Distance According to Frequency of Transmitter (m)		
	150 KHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = 1.2 \sqrt{P}$	$d = 1.2 \sqrt{P}$	$d = 2.3 \sqrt{P}$
0.01	0.12 m	0.12 m	0.23 m
0.1	0.38 m	0.38 m	0.73 m
1	1.2 m	1.2 m	2.3 m
10	3.8 m	3.8 m	7.3 m
100	12.0 m	12.0 m	23.0 m

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

USA and Canadian Radio Regulations

USA (FCC)

This device is equipped with Transmitter Module FCC ID: O7P-341 .

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- This device may not cause harmful interference.
- This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy. If not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the distance between the equipment and the receiver
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help

The user may find the following booklet prepared by the Federal Communications Commission helpful:

The Interference Handbook

This booklet is available from the U.S. Government Printing Office, Washington, D.C.
20402. Stock No. 004-000-0034504.

Mortara is not responsible for any radio or television interference caused by unauthorized modification of the devices included with this Mortara product, or the substitution or attachment of connecting cables and equipment other than specified by Mortara.

The correction of interference caused by such unauthorized modification, substitution, or attachment will be the responsibility of the user.

This device with transmitter module has been tested to SAR and complies with FCC exposure requirements for portable devices. SAR testing has been done at a distance of 0 mm from the body.

CANADA (ISED)

This device complies with Industry Canada's license-exempt RSSs. Operation is subject to the following two conditions: (1) This device may not cause interference; and (2) This device must accept any interference, including interference that may cause undesired operation of the device. Le présent appareil est conforme aux CNR d'industrie Canada applicable aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes: (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage, est susceptible d'en compromettre le fonctionnement.

8. INTRODUCTION

General Information

This user manual provides information for users of the Mortara Surveyor S2 telemetry device and its related components. It is written for clinical professionals expected to have a working knowledge of medical procedures and terminology as required for monitoring cardiac patients.

The Surveyor S2 telemetry device is a small, lightweight device designed to acquire an ECG and transmit this data to a MONITORING NODE.

This user manual explains:

- Relevant safety cautions, warnings and notices.
- How to configure, use and maintain the S2 telemetry device as well as its components and accessories.
- How to change the device configuration to adapt to the hospital IT network environment.
- How to properly acquire and transmit patient ECG signals to a MONITORING NODE.

Please also refer to the MONITORING NODE user manual for additional information.



WARNING: *The Surveyor S2 telemetry device by itself is not intended to be used as a vital signs physiological monitor.*



CAUTION: *The Surveyor S2 is has not been designed for direct cardiac applications involving direct contact with a patient's heart.*



CAUTION: *The Surveyor S2 has not been designed to be used together with electrosurgery equipment.*

NOTE: *This manual may contain renderings of various display screens. Any screen images are provided for reference only and are not intended to convey actual operating techniques.*

Indications For Use

- The Surveyor S2 is indicated for use in adult, adolescents, and children patient populations. The Mortara Surveyor S2 facilitates the monitoring of ECG signals.
- The Surveyor S2 is a prescription device intended to be used by knowledgeable healthcare professionals within a healthcare facility or clinical pharmacology unit.
- The Surveyor S2 is indicated for use in a clinical setting by a physician, or by trained personnel acting on the orders of a licensed physician. It is not intended as a sole means of diagnosis.
- The Surveyor S2 is indicated for use as a radiofrequency physiological signal transceiver, receiving and delivering real-time acquisition and transmission of simultaneous electrocardiographic data, while allowing the patient to be ambulatory within the range of the antenna network.
- The Surveyor S2 is indicated for use to acquire and output electrocardiographic data.

System Overview

The S2 telemetry device provides a means to acquire and transmit simultaneous ECG data to a MONITORING NODE while allowing the patient to be ambulatory within the coverage range of the WLAN network.

The S2 uses one (1) rechargeable battery pack.

The following equipment is necessary to use the S2:

- One (1) Rechargeable Li-Ion Battery Pack
- Applicable patient cable, lead wires, and electrodes.
- Carrying pouches.
- Battery charger.
- Customer provided 802.11 a/b/g/n 2.4 GHz (a/g/n is recommended) or 5 GHz wireless access points and network infrastructure.
- MONITORING NODE monitoring system.

Available configurations

The S2 is offered in the following, factory-set parameter configurations:

- 4-wire ECG


INTRODUCTION

4-wire ECG

When configured with the 4-wire ECG port, the S2 provides cardiac telemetry using a shielded 4-wire ECG cable.

Operating controls and indicators

The operating controls and indicators of the S2 are shown in the following diagram:

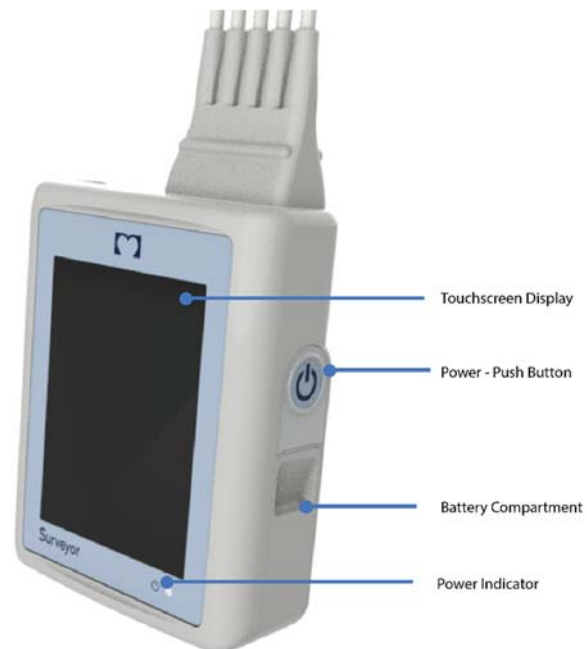

Display Screen

Screen allows for display of icons, menus, waveforms, and other device, patient, and status information.

The S2 uses a touchscreen user interface for operator interaction which allows control over the device via selection of menu items and icons as well as entry of relevant configuration and patient information by pressing on the screen.

Power Push-Button

By pressing and holding this button once, it turns the S2 device on. The Waveforms Screen should display by default.

After 60 seconds of inactivity the device will enter a low power state where the display will turn off. In this state, the display screen will be off but the functionality (i.e. ECG telemetry) will continue. To turn the display back on, hold the Power Push-Button.

Power Indicator

This bi-color LED will indicate these conditions in the following priority :

1. Flash green rapidly when: Device On, not Connected (Online or Offline)
2. Flash yellow when: Device On, any ECG lead in fail
3. Flash green slowly when: Device On and Connected

Battery Compartment

This section of the device stores the Li-Ion battery.

The S2 device is designed to operate on a single rechargeable Li-Ion battery.

The battery is replaceable without the operation of a tool.

Carrying Pouches


The pouches are designed to fit the contours of the S2 and accommodate extended wearing of the device while the patient is ambulatory or stationary. The transparent film allows viewing of the screen.



CAUTION: Pouches are designed for single patient use, and should NOT be reused.



CAUTION: While the adhesive pouches are made of biocompatible material and contain no latex, they should NOT be applied directly to the patient's skin.

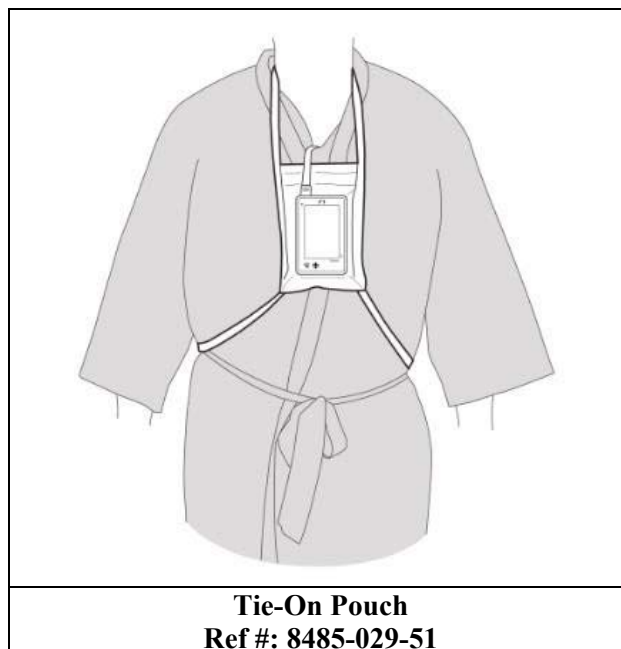

Battery Charger

The following Li-Ion Battery Charger RRC-SCC1120 – Single Cell Charger 1120 can be used to charge the rechargeable Li-Ion Batteries.



WARNING: To avoid the risk of electric shock, battery chargers must be connected to a supply main with protective earth.

To charge rechargeable Li-Ion batteries, load a battery into an open battery bay. The battery fits snugly in its charging bay in only one orientation, It snaps in place when slid down onto the charging contacts on the battery bay.



The LED on the battery charger indicate the status of the battery charging: Orange light indicates that the inserted battery is of the correct type and is currently being charged; Green light indicates that the battery is charged and can be removed for use; Red blinking indicates battery detection phase; Red light indicates no battery inserted, battery over/under temperature, charger over temperature battery over voltage, battery charge timer time-out-error or input voltage too low.

9. UNPACKING AND SET UP

Checking Contents

Your S2 is shipped with the following components:

- Surveyor S2 Telemetry device
- Li-ion rechargeable battery
- DVD containing the S2 User Manuals

The following optional S2 accessories are available separately:

- Pouch Options:
 - Disposable Tie-On Pouch
- Battery Charger

Battery Installation

NOTE: *If using the S2 rechargeable Li-ion battery, ensure it is fully charged prior to first use.*

The battery compartment is accessible via the removable battery door on the right side of the device.

1. Remove the battery door by firmly pushing down along the groove of the door and remove.
2. Load the battery into the battery compartment:

Align the battery so that once it is slid into place it adheres with the contacts

- a. Insert the battery into the compartment and slide to lock it mechanically in place.
3. Replace the battery door. As shown below, align the wedges of the door with the grooves of the S2 then push it into place until it is sealed completely.



10. PATIENT PREPARATION FOR QUALITY ECG

Quality ECG Data Acquisition

Obtaining quality ECG data is important in continuous ECG monitoring. A quality ECG signal depends largely on the patient prep and electrode placement. Good contact between the electrodes and the patient's skin and correct placement of the electrode can help ensure obtaining quality ECG data.

A good quality ECG contains:

- Clearly discernible P waves, QRS complexes, and T waves.
- Steady, even, crisp baseline.
- Absent of respiratory variability, artifact, noise, and other interference.

A good quality ECG may enhance the performance of the Monitoring Node and may lessen false erroneous alarm notifications.

A poor quality ECG may be caused by many factors:

- Poor site preparation.
- Poor electrode application or failure to refresh electrodes regularly.
- Patient movement.
- Interference by other equipment in the room.
- Poor quality ECG becomes synonymous with artifact and interference in the ECG waveforms.

A poor quality ECG may manifest in several ways:

- Fast baseline artifact.
- Erratic baseline.
- Sharp “spikes.”
- Rolling, wandering waveforms as seen with patient breathing patterns.
- Difficult to discern P waves or atrial fib waves from noise.
- Inability to discern P waves, QRS complexes, T waves.

Artifact and interference in the ECG waveforms may be caused by using accessories, lead wires, and ECG cables other than those specified to work with the S2. Always use accessories, lead wires, ECG cables, and other accessories specified to work with the S2.

Skin Preparation

In continuous ECG monitoring, the goal of skin preparation is to minimize the contact resistance between the patient's skin and the ECG electrode. Follow the facility's standard of care when preparing the patient's skin for ECG electrode placement and ECG monitoring.

To prepare the patient's skin for electrode placement:

1. Explain the procedure to the patient.
2. Maintain patient privacy during skin prep and electrode placement.
3. Locate the correct anatomical landmarks for electrode placement.
4. Clip or shave excess hair in the areas marked for electrode placement.
5. Remove residual skin oils, creams, and lotions by gently abrading the skin with a small gauze pad.

NOTE: With elderly or frail patients take care to not abrade the skin causing discomfort or bruising. Clinical discretion should always be used in patient preparation.

Electrode Placement


To apply electrodes:

1. Use pre-gelled, Ag/AgCl disposable electrodes.
 - a. Do not use electrodes after their expiration date, or if the gel has dried out.
 - Store electrodes in an air tight container.
 - Electrodes dry out if not stored properly leading to loss of adhesion and conductivity.
 - b. Always use the same electrodes. Do not mix electrode brands or types. Using different types of electrodes may cause baseline artifact and noise in the ECG tracing.
2. Apply the electrodes in the following manner:
 - a. Attach the electrode to the ECG lead wires prior to attaching the electrode to the patient's chest.
 - b. Place the electrode in the properly prepared, correct location by using a circular motion on the electrode adhesive area.
 - c. Gently press the electrode adhesive to the patient's skin until the entire outer surface of the electrode is adhered to the patient's chest.
 - d. Once the electrode adhesive is attached, gently press on the gel area to ensure proper gel to chest contact. Avoid dislodging the gel as the displaced gel can increase baseline artifact and noise in the ECG tracing.
 - e. Test for firm electrode contact by slightly tugging on the electrode to check for adhesion among the entire electrode surface. If the electrode moves freely, change the electrode. If the electrode does not move easily, a good adhesive contact has been obtained.

Refer to the Electrode Location section below for further details on correct anatomical landmarks for electrode placement.


Best Practice Recommendation: Change electrodes as per hospital standard of care, or at least every 24 hours to enhance patient skin care and the acquisition of quality ECG data. Clinical discretion should always be used in patient preparation.

Electrode Locations for 4-Wire Cable

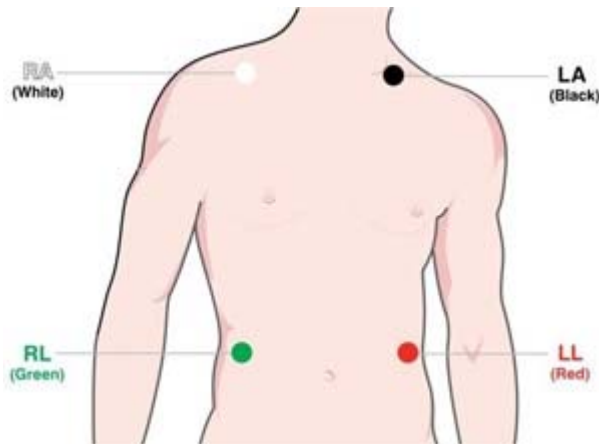

4-Wire Lead Placement (AHA)

- Place the RA (white) electrode under the patient's right clavicle, at the midclavicular line within the rib cage frame.
- Place the LA (black) electrode under the patient's left clavicle, at the midclavicular line within the rib cage frame.
- Place the LL (red) electrode on the patient's lower left abdomen within the rib cage frame.
- Place the RL (green) electrode on the patient's lower right abdomen within the rib cage frame.

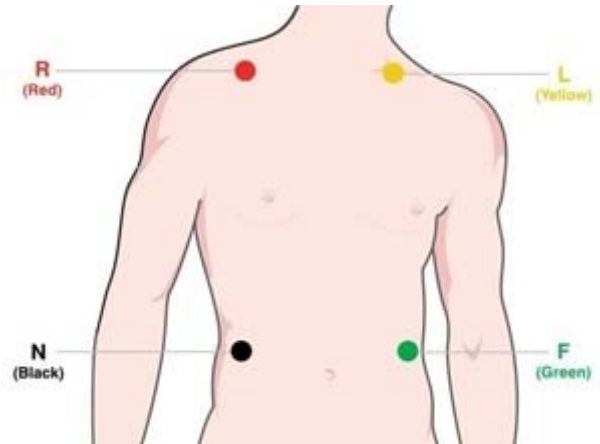

4-Wire Lead Placement (IEC)


- Place the R (red) electrode under the patient's right clavicle, at the midclavicular line within the rib cage frame.
- Place the L (yellow) electrode under the patient's left clavicle, at the midclavicular line within the rib cage frame.
- Place the F (green) electrode on the patient's lower left abdomen within the rib cage frame.
- Place the N (black) electrode on the patient's lower right abdomen within the rib cage frame.

Pacemaker Patients


Pacemaker patients may require a modified electrode placement based on the physical location of the patient's pacemaker generator device. Do not place an ECG electrode directly over the pacemaker generator as this may lead to artifact and noise on the ECG tracings.

Best Practice Recommendation: Place the electrode patches 3 – 5 from the pacemaker generator area. For example, if the pacemaker generator is located in the left subclavian area, relocate the Left Arm electrode closer in towards the center of the chest.

Use the following lead placement for monitoring a pacemaker patient using 4-wire cable:

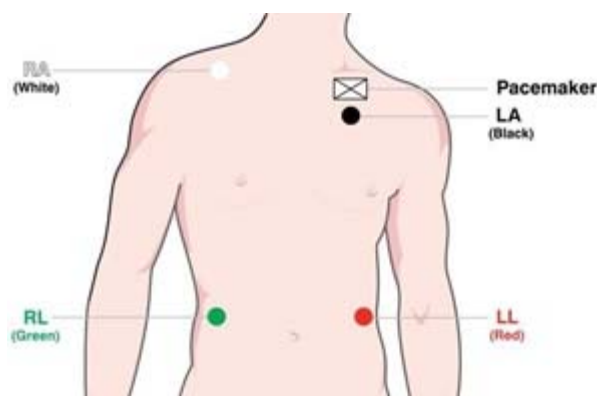

4-Wire Lead Placement for Pacemaker Patients (AHA)

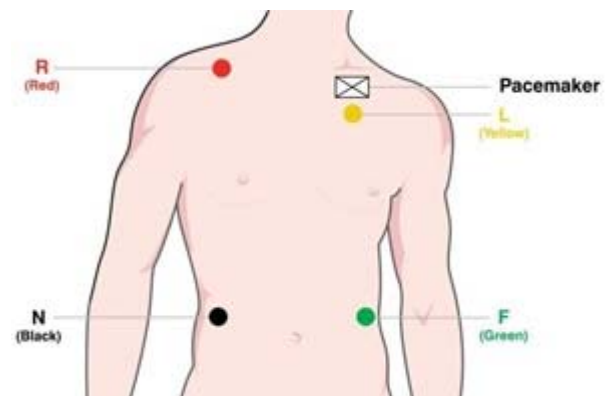

4-Wire Lead Placement for Pacemaker Patients (IEC)

Checking ECG Electrode and Lead Wire Signal Quality

Once the patient has been properly prepared, the electrodes have been attached in the correct anatomical location, and the patient ECG cable is connected to the S2, the S2's ECG screen displays the patient's ECG tracings.

Check to ensure the ECG tracing is free of artifact and noise with a clean ECG baseline as patient condition permits.

Refer to the Operation chapter of this user manual to check the quality and waveforms for each lead. If the ECG contains artifact or noise, review the steps for proper electrode site preparation and placement and repeat these steps to obtain proper ECG signals.

11. OPERATION

Powering On the S2

- To power ON the S2, place batteries in the unit as described in Chapter 9.
- Press and hold the Power Push-Button on the right side
- The power indicator will be lit green
- The display will flash blue momentarily and then display the Waveforms screen

Screen Time-Out & Reactivation

- To preserve power and to protect against accidental activation of menus, the display will turn off after 60 seconds of inactivity.
- To reactivate the display, press and hold the power Push-button until the screen is activated.

Powering Off the S2

To power OFF the S2, you can either remove the batteries or access the shutdown menu from the settings menu*. In the latter case, the shutdown menu also stops any active monitoring of telemetry data from the S2 on the MONITORING NODE.



WARNING: *Powering off the transmitter may terminate the monitoring of the patient and cease all alarms. Ensure that the patient is properly discharged or otherwise cared for prior to powering off the S2.*

Main Menu Screen

The following diagram shows the general layout of the S2's main functional screens:

Status Bar



MAIN MENU

Configure Device

Shutdown Device



The top portion of the Main Screen is the Status Bar. The remainder of the Main Screen shows one of the following functional displays:

- **Configure Device Option:** upon tapping on “Configure Device” active area, the device shall switch to Passcode Screen before proceeding to the Main Configuration Screen.
- **Shutdown Device Option:** upon tapping on “Shutdown Device” active area, the device shall switch to Passcode Screen before proceeding to Shutdown the device.

Navigation Bar buttons appear along the bottom of the screen, as appropriate to each display screen.

Status Bar Area

The top portion of the LCD screen shows status information regarding the current state of the device. It displays the (1) Time of Day, (2) WLAN and Connection Status, in addition to the (3) battery status as described below:

(1) Time The following information is shown: Displays the time of day in 24H format without seconds, as received from MONITORING NODE.	
(2) WLAN Connection Status The following information is shown: - When connected to WLAN, the status bar shall include an indication of the level of the signal. - When not connected to WLAN, the status bar shall display an indication in place of Signal Level and Communication Status	
(3) Battery Capacity Status The status bar shall display the remaining battery capacity in percentage.	

Time of Day

The current time of day is displayed on the top of the status area. The current date and time is synchronized between the S2 and the MONITORING NODE. If the connection between the MONITORING NODE is disrupted or the S2 battery removed the time information is resynchronized on the next connection to the MONITORING NODE.

WLAN Connection Status

The battery level indicator in the Status Bar Area displays one of the following symbols depending on the status of the WLAN Connection




The device was successful in establishing a connection to the wireless network.



The device was not successful in establishing a connection to the wireless network.

Navigation Bar

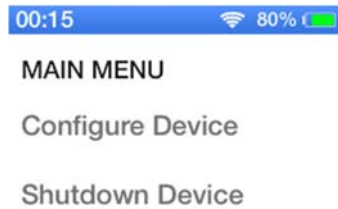


Navigation Bar Area



The bottom portion of the LCD screen contains the bar for switching between Monitoring Screens. It contains the (1) Main Menu Screen, (2) Waveform Screen, (3) Parameters Screen and (4) Demographics Screen as described below:

(1) Main Menu Button • Upon tapping on the Main Menu Screen Button, it will bring you to the Main Menu Screen • Here you will be able to shut down the device or configure it	
(2) Waveform Button • Upon tapping on the Waveform Button, it will bring you to the Waveform Screen • Here you will be able to view the ECG waveform	
(3) Parameters Button • Upon tapping on the Parameters Button, it will bring you to the Parameters Screen • Here you will be able to view the Heart Rate of the patient at a glance	
(4) Demographics Button • Upon tapping on the Demographics Button, it will bring you to the Demographics Screen • Here you will be able to view the patient's information	

Main Configuration Screen

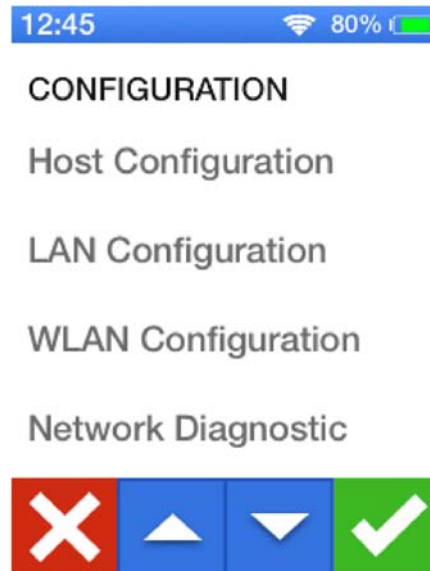

Located on the bottom of the Main Configuration Screen is the Configuration Navigation Bar. In this Configuration Navigation Bar, it contains buttons that will (1) Cancel, (2) Up and Down Buttons and (3) Confirm the actions made during the time of entering the Main Configuration Screen.

(1) Cancel Button

- Upon tapping on the Cancel Button, it shall discard all changes and switch back to the Main Menu Screen



(2) Up and Down Buttons

- Upon tapping on these buttons, it shall scroll the current area content in a per page fashion, in case not all elements can fit in a single pagearea.



(3) Confirm Button

- Upon tapping on the Confirm Button, it shall store and apply all changes through an automatic reboot. If there are no outstanding changes then it will simply return to the Main Menu Screen.



In order to get to the Main Configuration Screen, you have to input the correct passcode when prompted on the Passcode Screen pictured below.



The image shows a mobile device screen with a blue status bar at the top displaying '00:15' and '80%' battery. Below the status bar, the word 'Passcode' is centered. A vertical cursor is positioned below 'Passcode'. Below the cursor is a numeric keypad with four rows and four columns. The buttons are: Row 1: 1, 2, 3, <; Row 2: 4, 5, 6, abc; Row 3: 7, 8, 9, ?!#; Row 4: ESC, 0, ., OK.

1	2	3	<
4	5	6	abc
7	8	9	?!#
ESC	0	.	OK

- The default passcode for the S2 device is **7865**.
- The **ESC** button is pressed to cancel data entry and go back to the parent screen
- The **OK** button is pressed to confirm data entry to go back to the parent screen or proceed to the next screen
- The < (backspace button) is pressed to delete the last character in the entry

Upon entering the correct Passcode, it will proceed to the Main Configuration Screen. Here, you will be able to configure the S2 device in these respective fields, (1) Host Configuration, (2) LAN Configuration, (3) WLAN Configuration and view information such as (4) Network diagnostic, (5) Device diagnostic, and (6) Version.

(1) **Host Configuration**


- Upon tapping on the Host Configuration area, it will proceed to a Host Configuration Screen.
- In this view, you will be able to edit the following:
 - Central IP
 - Bed ID
 - Unit ID
 - Base Port
 - Port Number

12:45 80%

Central IP
172.17.1.20


Bed ID
0

Unit ID
1

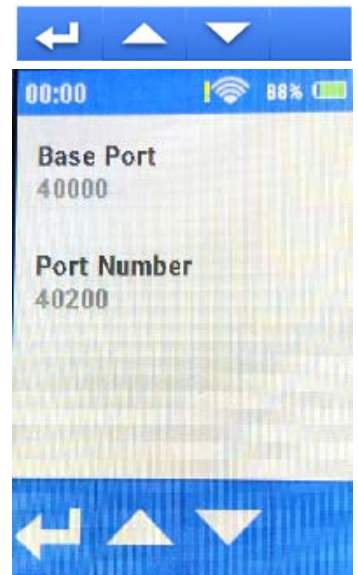

(2) **LAN Configuration**


- Upon tapping on the LAN Configuration area, it will proceed to a LAN Configuration Screen
- In this view, you will be able to edit the following:
 - o Enable DHCP
 - If this setting is set to true then the following settings are NOT able to be edited
 - IP
 - Subnet
 - Gateway
 - If this setting is set to false then the following are able to be edited
 - IP
 - Subnet
 - Gateway

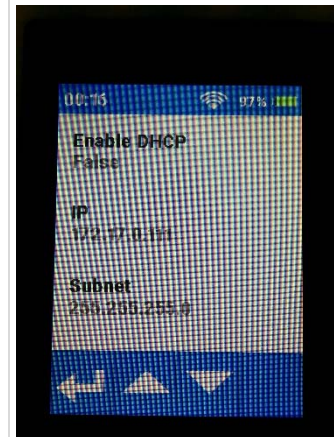

Pictured here is an example of when the "Enable DHCP" is enabled. The other fields are then editable.


(3) WLAN Configurations

- Upon tapping on the WLAN Configuration area, it will proceed to a WLAN Configuration Screen
- In this view, you will be able to edit the following:
 - o ESSID
 - o Passkey
 - o Band
 - This is a toggle to allow the user to switch between 2.4Gz and 5Gz

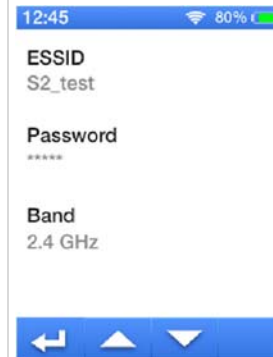


(4) Network Diagnostics

- Upon tapping on the Network Diagnostics area, it will proceed to a Network Diagnostics Screen
- In this view, you will be able to view the following information:
 - o Network Status
 - o IP Address
 - o AP MAC
 - o Station MAC
 - o Signal Strength
 - o Signal Noise
 - o Quality
 - o Channel
 - o Error
 - o Current ESSID

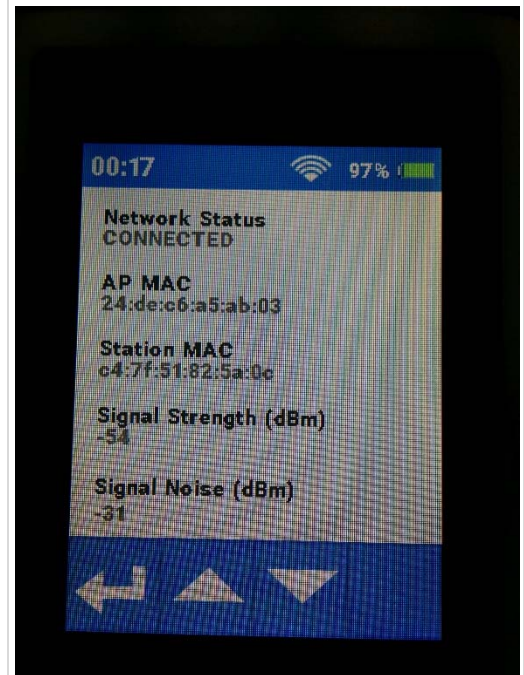


(5) Device Diagnostics

- Upon tapping on the Device Diagnostics area, it will proceed to a Device Diagnostics Screen
- In this view, you will be able to view if an error occurred and be able to clear it

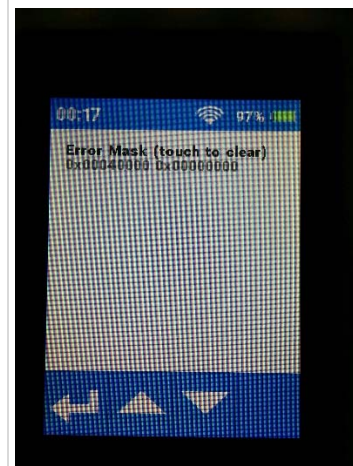

(6) **Version**


- Upon tapping on the Version area, it will proceed to a Version Screen
- In this view, you will be able to view the current version of the device

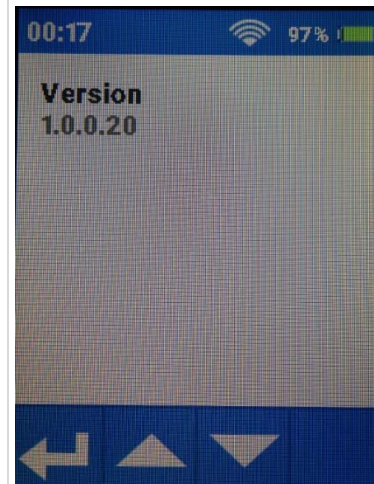

Waveform Screen

When initially powered on, the S2 will default to the Waveform Screen. It begins monitoring right away.

Main Waveforms Screen

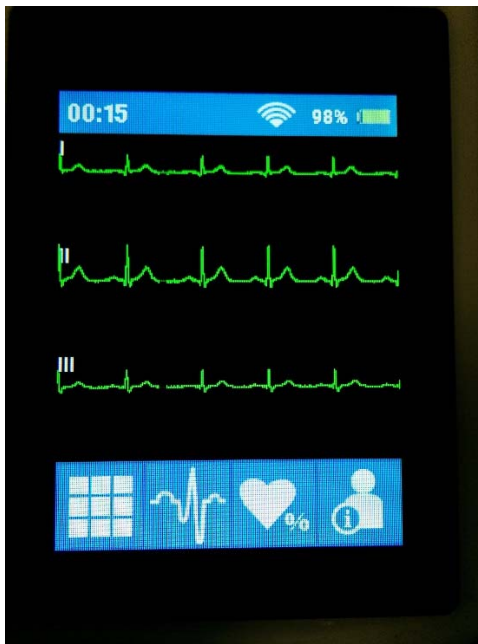


In the monitoring screen, the S2 displays the on-demand ECG waveforms as follows:

- 1) The S2 is capable of displaying derived ECG leads depending on the type of ECG cable attached
 - When a 4-wire cable is in use, the S2 is capable of displaying Leads I, II, and III.
- 2) The ECG waveform(s) is green in color.
- 3) The ECG waveform(s) will be drawn as a square wave when in lead fail.
- 4) The S2 displays a minimum 3 seconds of the ECG waveform(s)
- 5) ECG waveform(s) displayed is high pass filtered

Parameters Screen

When the S2 is in an active session, it is possible to display the patient's heart rate by pressing the Parameters Icon in the  Navigation Bar

Parameters Screen

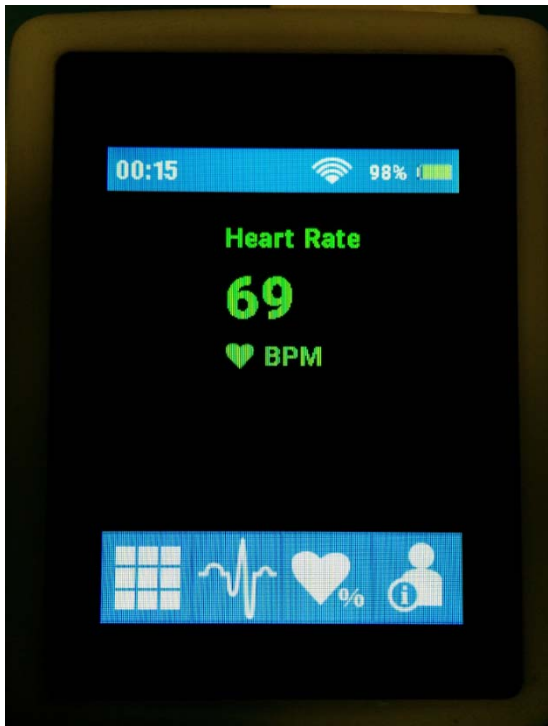

This will show the patient's current heart rate as reported by the MONITORING NODE.

***NOTE:** By design, no other parameters or alarms are shown on the S2's display. Refer to the MONITORING NODE user manual for relevant parameters, alarms and other information available there.*

Parameters Icons

1. The Heart Rate will display “???” when there is an **INVALID** parameter
2. The Heart Rate will display “?” when the parameter has **EXPIRED**
3. The Heart Rate will display “---” when the parameter is **DISABLED**
4. The Heart Rate as received from the MONITORING NODE will be displayed together with “bpm” unit of measure.

To return to the main screen, press the Main Menu Icon 

Demographics Screen

When the S2 is in an active monitoring session, it is possible to receive information from the MONITORING NODE about the demographics of the patient.

Parameters Screen

The screenshot shows a mobile application interface. At the top, the status bar displays the time 00:15, a Wi-Fi signal icon, and a battery level of 98%. Below the status bar, a dark header bar contains the text 'Slot Name: Bed S2T' and 'Running'. The main content area is titled 'PATIENT INFO' and displays a table of patient demographics. At the bottom, there is a blue navigation bar with four icons: a 3x3 grid, a heart rate line, a heart with a percentage sign, and a person icon with an information 'i' symbol.

Last name	First name
Smith	John
Patient ID	Sex
ID0001	Male
Date of Birth	Age
09/27/1976	41 Years

The following information will be displayed on the Main Area for sessions assigned to the patient and are acquired from the MONITORING NODE

1. First Name
2. Last Name
3. Patient ID
4. Date of Birth
5. Age
6. Sex

Demographics information will be replaced by “**Anonymous**” if no patient is assigned to the monitoring session.

The Status of the monitoring session (Running or Stopped) will be indicated as well as the bed/slot name.

The Age will be displayed in one of the following units of measurement:

1. Years
2. Months
3. Weeks
4. Days

Demographics not received for more than 50 seconds will be considered **EXPIRED** and all the fields will be replaced by “?”.

12. CONFIGURING THE S2


This chapter describes how to configure the S2 including changing the settings during the operation as well as general configuration for setup of the S2.

Select the configuration area settings when inside the Main Menu Screen.

NOTE: *This particular screen is a password-protected function and only available through a passcode. The S2 is factory shipped with passcode 7865. Refer to your facility's administrator for additional information.*


Upon entry of the passcode, the menus for this screen will be available. If the wrong passcode is entered, the “Unlock failed” message is displayed and the correct passcode must be entered to enter this screen. Alternatively, this operation may be cancelled.

Passcode Entry Screen

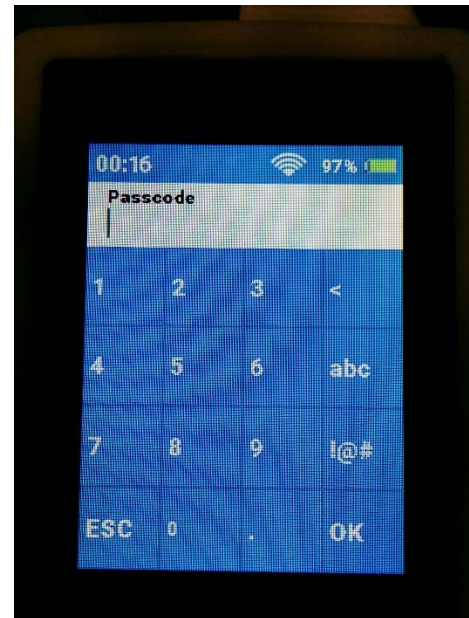

Upon successful entry of the passcode, the S2 is unlocked and the configuration settings screen is displayed as shown.


Configuration Settings Screen

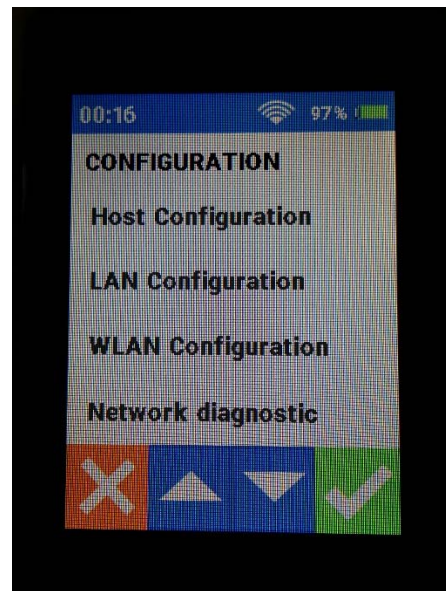

Host Settings

The settings for the MONITORING NODE including its IP address, and the unit and bed ID for this particular S2, as well as base IP port and IP port number are defined here. If you are uncertain about these settings, please consult Mortara or your IT Network Administrator.

Host Settings Screen

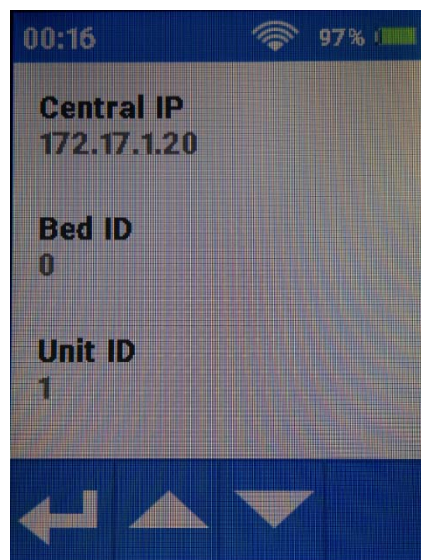


Network Settings


The network settings of the S2 are configured as shown below.

Network Settings Screen

LAN

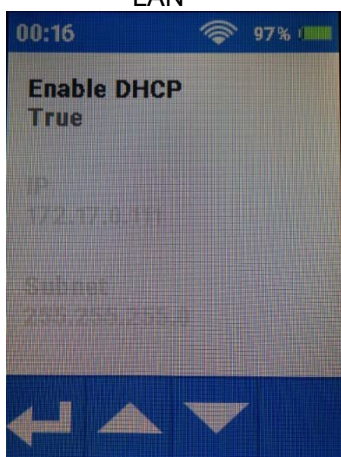

WLAN

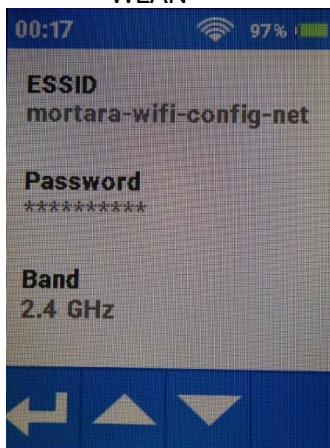

Available settings include:

- LAN
 - IP Address Method: Whether dynamically assigned by the network through a DHCP service or manually assigned here to a static address. If “False” is selected, the IP address, subnet mask and gateway address must be defined.
- WLAN
 - SSID, Band, & Passkey: The SSID and Passkey of the WLAN network that will be used for telemetry transmission.

If you are uncertain about these settings, please consult Mortara or your IT Network Administrator.


Language Settings

The S2's user interface supports only the English language as of now.

Device Diagnostics

This device diagnostics screen provides an easy access method to determine the current state of the device. It quickly assesses if the device has a current error as well as clear previous errors.


Device Diagnostics Screen

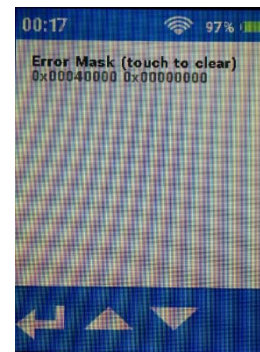

S2 Version

This screen displays the current version of the S2 software. In case of problems or issues, it may be necessary to share this information with Mortara.

S2 Version Screen

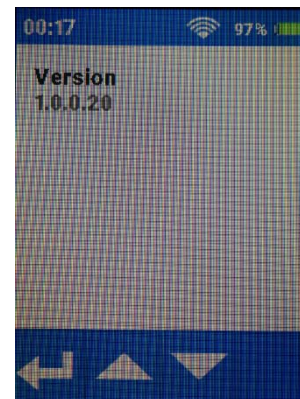

13. ECG MONITORING

To perform ECG acquisition on a patient using the S2 and MONITORING NODE:

Power-On

1. Insert freshly charged rechargeable batteries in the device.
2. Press the Power On button to power up the S2.

Patient Prep

3. Prep the patient and place electrodes per clinical guidelines and provided instructions.

Patient Cable Hookup

4. Attach the patient cable to the electrodes on the one end, and to the S2 on the other end.

ECG MONITORING

Confirm Good Signal

4. Confirm that all leads are OK. Select the waveforms icon () to

view and confirm good signal on all leads.



Patient Sessions

Sessions are assigned to a slot, managed, and ended at the MONITORING NODE.

15. PRODUCT SPECIFICATIONS

General

Instrument Type	Digital telemetry device with ECG.
Dimensions	3.5" x 2.8" x 1.1" (89 × 71 × 28 mm)
Input Channels	Continuous 3 lead signal acquisition and transmission with 4-wire cables.
ECG Leads Transmitted	4-wire: I, II and III
Device Classification	Type CF, battery operated
Special Functions	ECG display, lead fail, battery notification
Defibrillator Protection	IEC 60601-2-27
Function Keys	Control button; touchscreen menu navigation
Display type	High definition, antiglare Color TFT-LCD with LED backlight and capacitive touch panel
Display size	Display size: 77.3mm x 56.1mm Display Active Area: 57.6mm x 43.2mm Display Visual Area: 58.6mm x 44.2mm
Display resolution	240 x 320 pixels

Environmental

Temperature	Operating temperature: +10° to +40° C (+50° to +104° F) Storage temperature: -40° to +70° C (-40° to +158° F)
Humidity	Operating humidity: 10% to 95% RH, non-condensing Storage humidity: 10% to 95% RH, non-condensing
Altitude	Operating & Storage: 62 kPa to 106 kPa corresponds to an approx. altitude range of -260 m to 5000 m (-850 ft to 16404 ft)
Cooling	Passive (no fan)
Weight	< 200g (7.1 oz.) with rechargeable battery pack

Power Requirements & Battery

Rechargeable Battery Type	Rechargeable Lithium-ion battery pack
Rechargeable Battery Run-Time	32 hours (when configured with ECG only)
Battery Charging Time	4 hours

ECG Specifications

ECG	With 4-wire ECG
Simultaneous Leads Available	4-wire: I, II and III
Acquisition Rate:	40,000 samples/s. for pacemaker detection, 500 samples/s for transmission
Resolution:	.117 μ V originally, reduced to 2.5 μ V for transmission
Dynamic Range	$\pm$ 350mV
CMRR:	According to applicable standards (IEC 60601-2-27)
Max. Auxiliary Patient Current:	< 10 μ A
Frequency Response:	0.05 to 150 Hz
Input Impedance:	>2.5 M Ω at 10 Hz
Electrodes:	Compatible with Ag/AgCl disposable electrodes
HR Report from the MONITORING NODE	Less than one second. Heart rate is only available when the S2 is actively communicating with the MONITORING NODE.
Recovery from Defibrillation Discharge	Less than 5 seconds, compliant with IEC 60601-2-27.

Wireless Network Specifications

Wireless protocol	802.11 a/b/g/n 2.4 GHz or 5 GHz Disable 802.11b support on WLAN infrastructure. If 802.11b support is required, do not specify any 802.11b data rate as mandatory.
Frequency Range	2.400 GHz to 2.483 GHz; 5.15 GHz to 5.35 GHz, 5.725 GHz to 5.825 GHz
802.11a/g supported rates	6, 9, 12, 18, 24, 36, 48, 54 Mbps
802.11n supported rates	MCS0 through MCS7 6.5, 13, 19.5, 26, 39, 52, 58.5, 65 Mbps
Transmit Power	15 dBm Configure WLAN APs maximum transmitting power to 15 dBm
Channels	1, 6, 11 with 2.4GHz 36, 40, 44, 48, 149,153,157, 161 with 5GHz
Wireless security protocols	WPA2-PSK-AES, passphrase
WMM / QoS	Yes
SSID Requirements	Dedicated SSID is required. Hidden SSID not supported. SSID can be up to 24 ASCII characters long, per 802.11. SSID prioritization and minimum bandwidth configuration on WLAN infrastructure is recommended. Period Session Timeout or similar policies should be disabled or excluded in the SSID dedicated to the S2: these policies might cause data gaps of a few seconds whenever the AP terminates the session.
IP address assignment	Static or DHCP
IP protocols and ports	UDP: one port per device, according to the system configuration (Base port number 30000). Actual ports will be communicated upon definition of the system configuration; ICMP
Maximum devices per Access Point	16
Data Rate (payload)	Outbound: < 80 kbps Inbound: < 24 kbps
Minimum RSSI / SNR	-70 dBm in the coverage area
Minimum SNR	20 dB in the coverage area
Target total channel occupancy	< 50%

16. TROUBLESHOOTING

The following table provides guidance for investigating issues that may occur during operation of the S2. Contact Mortara Technical Support and Service by calling 1-888-MORTARA (US & Canada) or +1-414-354-1600 (Worldwide) or via e-mail at techsupport@mortara.com for further assistance.

Power and Battery

Symptom	Possible Causes	Suggested Resolution
The S2 is not working and display does not light up.	Internal system failure.	Power cycle the S2 with fresh batteries and try again. If problems persist, stop using the S2 and contact Mortara Technical Support.
	Battery will not hold charge.	Replace the battery according to instructions in the General Care & Maintenance chapter.

Display and Touchscreen

Symptom	Possible Causes	Suggested Resolution
The touchscreen is not working properly.	Touchscreen failure.	Power cycle the S2 and try again. If problems persist, stop using the S2 and contact Mortara Technical Support.
The display is not working properly.	Display failure.	Power cycle the S2 and try again. If problems persist, stop using the S2 and contact Mortara Technical Support.

ECG Trace

Symptom	Possible Causes	Suggested Resolution
ECG signal is noisy.	Excessive patient movement.	Confirm electrode site preparation; confirm correct ECG positioning, move electrodes if needed.
	Electrical noise from auxiliary equipment.	Isolate the patient from auxiliary equipment, if possible.
	Bad electrode contact.	Check to make sure electrodes are still securely attached to the patient, and reattach if necessary. Remember the importance of good skin preparation techniques.

Network Transmission

Symptom	Possible Causes	Suggested Resolution
1. X symbol on WLAN signal indicator. 2. "AP MAC: Not-Associated" in the WLAN diagnostic screen. 3. Invalid IP address if configured for DHCP. 4. Inability to ping in WLAN diagnostic screen.	Incorrect WLAN configuration (e.g., SSID, Security protocol, Passkey, Passkey type).	Check WLAN network settings; contact your IT Network Administrator to correct issues.
	WLAN infrastructure issues (e.g. device placed on black list).	Check WLAN network settings; contact your IT Network Administrator to correct issues.
1. Normal symbol on WLAN signal indicator. 2. Network info correctly populated in WLAN diagnostic screen. 3. Invalid IP address if configured for DHCP. 4. Inability to ping in WLAN diagnostic screen.	Incorrect LAN configuration (Method, IP, Subnet, Gateway).	Check LAN network settings; contact your IT Network Administrator to correct issues.
	LAN infrastructure issues (e.g. DHCP server is not providing an IP address).	Check LAN IPv4 network settings; contact your IT Network Administrator to correct issues.

TROUBLESHOOTING

Symptom	Possible Causes	Suggested Resolution
Any of the following: 1. Intermittent traces on MONITORING NODE, 2. X or low signal symbol on WLAN signal indicator. 3. Intermittent "AP MAC: Not-Associated" in WLAN diagnostic screen, low signal, high noise, low quality indexes. 4. Intermittent Valid IP address if configured for DHCP. 5. Inconsistent Ping results in WLAN diagnostic screen.	Poor WLAN network coverage.	Check IT Network Administrator to correct issues.

17. REORDERING ACCESSORIES & CONSUMABLES

Use the following Mortara part numbers to obtain spare parts or to reorder accessories:

Power accessories

Description	Part Numbers
LI-ION BATTERY PACK	4800-020
BATTERY CHARGER	4102-003

Mounts / holders

TIE-ON DISPOSABLE POUCH - BOX 100**	8485-029-51
-------------------------------------	-------------

ECG accessories

ECG LEAD SET 4 WIRE CLIP ENDS AHA GRAY	9293-059-80
ECG LEAD SET 4 WIRE CLIP ENDS IEC GRAY	9293-059-81
ECG LEAD SET 4 WIRE SNAP ENDS AHA GRAY	9293-059-82
ECG LEAD SET 4 WIRE SNAP ENDS IEC GRAY	9293-059-83
HOOKUP KIT MONITORING 10E SINGLE**	9294-009-50
ELECTRODES, 10/POUCH **	9300-032-52
ELECTRODES, 300/CASE (10/pouch, 30 pouches/case) **	9300-032-50

Other

S2 USER MANUAL	9515-210-50-CD
----------------	----------------

* Technical Service Installation Required.

** This item is intended for single patient use. Items warranted to be free of defects in workmanship and materials for a period of 90 days or first use, whichever comes first.

To order additional supplies, contact your Mortara Customer Service Representative.

18. APPLIED STANDARDS

IEC 60601-1:2005/A1:2012-Ed.3.1, Medical electrical equipment – Part 1: General requirements for basic safety and essential performance

European Union: CENELEC EN60601-1:2006 Medical Electrical Equipment Part 1: General Requirements for Basic Safety and Essential Performance.

USA: AAMI ES60601-1:2005 Medical electrical equipment Part 1: General Requirements for Basic Safety and Essential Performance.

Canada: CSA C22.2#60601-1:2008 Ed:3 Medical Electrical Equipment – Part 1: General requirements for basic Safety and essential performance.

IEC 60601-1-2:2014 Medical Electrical Equipment – Part 1-2: General requirements for safety – Collateral Standard: Electromagnetic Compatibility.

IEC 60601-2-27:2011-08 Edition 3.0, Medical Electrical Equipment. Part 2-27: Particular requirements for the basic safety and essential performance of electrocardiographic monitoring equipment.

(Notes: Applies to 3/5-wire and 10-wire acquisition modules; when used in conjunction with a MONITORING NODE.)

ETSI EN 301 489-1:2017-02 v2.1.1, Electromagnetic compatibility and Radio Spectrum Matters, Part 1: Common technical requirements.

ETSI EN 301 489-17:2017-02 v3.1.1, Electromagnetic compatibility and Radio Spectrum Matters, Part 17: Specific conditions for Broadband Data Transmissions Systems.

ETSI EN 300 328:2016-11 v2.1.1, Electromagnetic compatibility and Radio Spectrum Matters, Wideband transmission systems; Data transmission equipment operating in the 2.4 GHz ISM band and using wide band modulation techniques.

ISO 10993-1:1997-12-15, 2nd Edition, ISO 10993-5:2009/(R)2014, 2nd Edition, ISO 10993-10:2002-09-01, 2nd Edition ISO 10993-1:2009/Cor.1:2010, Biological evaluation of medical devices – Part 1: Evaluation and testing within a risk management process

AAMI/ANSI/ISO 10993-5:2009/(R)2014, Biological evaluation of medical devices – Part 5: Tests for in vitro cytotoxicity.

ISO 10993-10:2010, Biological evaluation of medical devices – Part 10: Tests for irritation and skin sensitization

Title 47 CFR Part 15, US Federal Communications Commission Part 15 Rules for Intentional Radiators.

2014/53/EU: European Council Directive on the harmonization of laws relating to radio equipment.. (RED Directive)

2002/96/EC:2003, Waste Electrical and Electronic Equipment (WEEE).

2002/95/EC:2003, Restriction of the use of Certain Hazardous Substances in Electrical and Electronic Equipment (RoHS).

IEC 62366:2014-Ed.1.1, Medical devices – Application of usability engineering to medical devices

93/42/EEC: Council Directive of 14 June 1993 concerning medical devices. (Medical Device Directive)