


Oclean Smart Sonic Electric Toothbrush

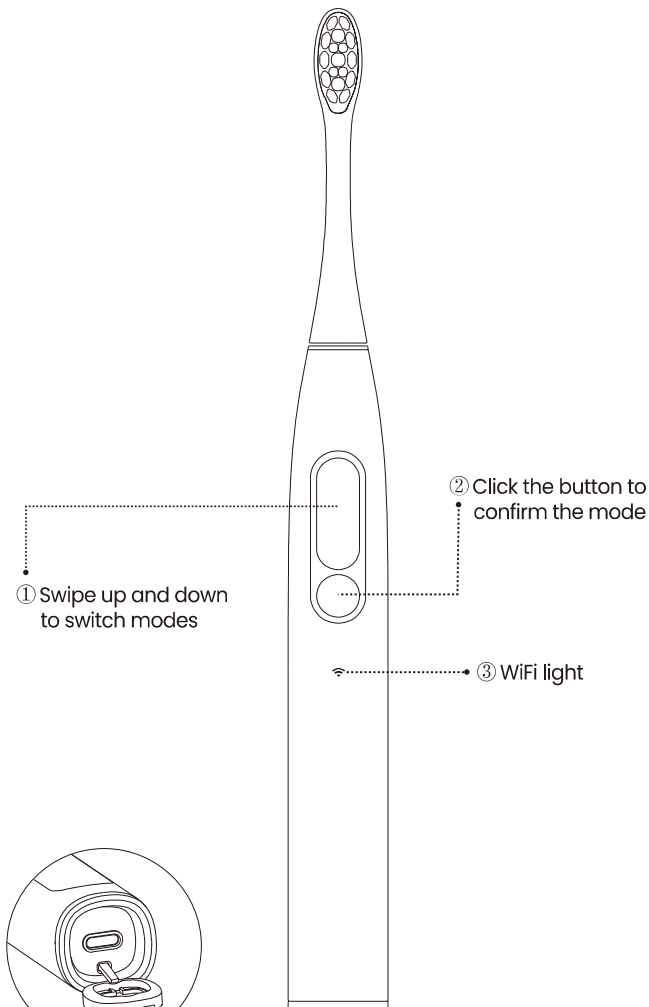

Detailed component description of the machine

1

English

Thank you for using Oclean electric toothbrush

Important safety tips

- Children aged 8 years and over and persons with reduced physical, sensory or intellectual capabilities or lack of experience and knowledge may use this appliance as long as they are supervised or instructed on safe use of the appliance and understand the hazards. Children are not allowed to play with this appliance. Children are not allowed to perform cleaning and user maintenance without supervision.
- This product can only be used with Oclean special brush heads, and we will not be responsible for any consequences caused by non-certified third-party brush heads.
- This product is only for personal cleaning and is not recommended for medical use.
- This product is for indoor use only. Do not place it in direct sunlight. It is recommended to place it in a dry and ventilated place at room temperature.
- Please keep the body dry before charging, wipe it dry with a soft cloth, and clean the charging port with a dry cloth regularly to prevent stains from accumulating.
- The appliance can only be powered at the safety extra-low voltage corresponding to the mark on the appliance.
- Do not expose the device or battery to excessively high temperatures. The normal temperature range for battery charging is 0-45 degrees.
- Be careful not to let metal objects touch the terminals of battery-powered devices or batteries to avoid the risk of short circuit.

Brush Head

When replacing the brush head, please pull it out or insert it vertically, and do not rotate it left or right.
Dentists recommend replacing your brush head every 2 to 3 months.
People with oral sensitivity, such as those in the following situations, should pay attention:
People who have just undergone oral periodontal surgery or are currently undergoing treatment.
People who have other symptoms in the mouth, including loose dentures or persistent heavy bleeding from the gums.
Other people with oral sensitivity, etc.
It is not recommended to use this product or consult a professional dentist before use.
This product is not suitable for use by disabled persons, mentally ill persons or persons lacking relevant experience and knowledge [including children (0-12 years old)] unless used under the supervision and guidance of their guardians.

Usage Guide

- Lift the screen to wake up the interface, slide the screen up and down to switch the brushing mode, and click the button to start brushing. During brushing, press and hold the button for 2s to cycle through the modes.
- 5 brushing modes:
"Sensitive Care" - specially designed for people with sensitive teeth, it gently cleans while reducing irritation to sensitive areas.
"Quick Clean" - for efficient daily cleaning, effectively removing food residues and dirt.
"Deep Clean" - goes deep between teeth and hard-to-reach areas to powerfully remove stubborn stains and bacteria.
"Whitening Polishing" - Intense cleaning makes the tooth surface smoother and whiter.
"Personal customization" - users can intelligently formulate oral care plans on the APP according to their own gum and oral conditions, such as: brushing intensity and brushing time settings, etc. After completing the settings, the toothbrush will synchronize the currently set mode.
Therefore, you can choose different modes according to your preferences and the dentist's advice.
- When the APP is configuring WiFi, the WiFi light will breathe. When WiFi is connecting, the white light will breathe. When WiFi is successfully connected, the white light will be on.
- Restore factory settings: Press and hold the button for 8 seconds to restore the toothbrush to factory settings.



Google Play Store & iOS APP Store

Install the APP: After scanning the QR code to download the "Oclean Care+" APP, please search and connect the toothbrush in the APP.
Note: The connection in the mobile phone's "System Settings->Bluetooth Settings" is invalid.
APP connection (Homepage -> Add device -> Select WiFi -> Enter password -> Connecting -> Connected successfully)

Healthy brushing

- Move the toothbrush slowly with the bristles at a 45-degree angle to the gum line;
- During brushing, the vibration of the toothbrush will pause every 30 seconds, prompting you to change the brushing position.
- If you do not pause during use, the device will automatically shut down after the brushing mode ends.

Charging method:

Wipe the bottom of the toothbrush dry before charging, and open the Type C port at the bottom of the machine to charge. While charging, the screen interface will display the charging status at the same time.

Low battery reminder:

When the toothbrush is running low on battery, the screen will remind you that it needs to be charged.



- Do NOT disassemble, impact, squeeze or throw into fire.
- If the product is severely swollen, do not continue to use it.
- Do NOT expose electrical appliances or batteries to excessively high temperatures.
- Do NOT use the battery after it has been immersed in water.
- Be aware of the danger of short-circuiting the terminals of battery-powered appliances or batteries by metal objects.

2

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID :2AN50-Y2089

CE WARNING

- the appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction
 - children being supervised not to play with the appliance
 - Signal words WARNING, CAUTION, DANGER in uppercase having a height as specified: Uppercase letter of the text explaining the signal word not smaller than 1,6 mm: Moulded in, engraved, or stamped markings either raised above or have a depth below the surface of at least 0,25 mm, unless
 - This appliance contains batteries that are non-replaceable
1. Cleaning and user maintenance shall not be made by children without supervision.
2. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

Before disposing of the product, the battery must be removed from the product according to the following steps:

- Open the bottom waterproof cover and remove the screw plug.
- Use a screwdriver to unscrew the bottom screws and open the bottom cover.
- Use a screwdriver to open the clips on both sides, take out the movement, use pliers to cut off the battery nickel sheet pins from B+ and B-, and take out the battery.

- Before discarding the device, the battery must be removed from the device;
- When removing the battery, the appliance must be powered off;
- Batteries should be disposed of safely.

Product Specifications

Product name	Oclean X Elite WiFi Smart Sonic Electric Toothbrush
Product model	Y2089
Rated input	5V=0.5A*
Rated input power	2.5W
Battery capacity	1100mAh
Compatible devices	compatible with Android 5.0, IOS 9.0 and above, support Bluetooth 4.2
Frequency range	2400-2483.5MHz
Waterproof rating	IPX7
Charging method	Type-C

*When charging an adapter with an output of 5V=0.5A (the maximum power of the output is less than 15W) and in compliance with the requirements of EC 60335-1 must be used for charging

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer Address: 28G, Building 3, Dachong Business Center (Phase 3), No.18, Dachong 1st Road, Dachong Community, Yuehai Street, Nanshan District, Shenzhen City, Guangdong Province, China



EU Authorized Representative
Agency Name: Mbrush UAB

European Address: Krokuvos str. 13,
09314 Vilnius, Lithuania

Email: service@mobat.eu
Post Code: LT44237

Made in China

FDA Registration Number: 3014246581



3