图号:E3-EN

材质:128g铜版纸 封面需要过哑油固色处理 成品尺寸:105×143mm

颜色: PANTONE3252C PANTONE8C Black

装订方式: 骑马钉

产品名称统一: Eye Massager 页码目录更新 蓝牙名称变更 20201020 增加海外认证用警告用语 20210413



Thank you for choosing our Eye Massager. Please read this manual carefully before use in order to use the product correctly and take full advantage of its excellent performance. We believe it will bring convenience and comfort to your life. Please keep this manual for your future reference. If any problems occur when using the product, please contact SKG after-sales service outlet or our company. We will be glad to help you. Thank you!

•	4
Troubleshooting	4
Precautions ·····	5
Technical parameters	7
Electrical schematic diagram	7
Product diagram ·····	8

This product is a household product. Not intended for commercial use.

As we continually improve our products, we shall reserve the right to make technical changes without prior notice.

Preparation

Please check whether the battery has run out. If the battery is too low, which leads to failure to boot, please charge the product before use. After the product is fully charged, make sure to remove the power adapter to avoid long-time charging which reduces product life

Charging

- 1. When the battery is low, the massager will prompt "Please charge if the battery is low." Please adjust the product to proper angle, expose charging interface and insert the USB cable. 2. During charging, the red indicator light is on. After fully charged,
- the green indicator light is on. 3. During charging process, the product doesn't support starting up. (the product initiates self-protection procedure to stop working automatically when it's charged).
- 4. The massager adopts built-in chargeable battery with charging time of 3-5 hours.
- 5. During usage, if the voice prompts low battery alarm, please stop use and reuse after full charged.
- Multiple charging modes

	Computer USB interface	
Charging modes	Mobile power	
	Standard power plug 5V 1A(recommended)	

Usage time

The default working time is 15 minutes. To prolong usage time, reboot the massager.

How to use

- 1.Long press (1) for two seconds to turn on the massager with the voice prompts "HI, SKG". Short press to switch modes. The product has five modes:
- Active mode: air pressure + vibration + hot compress + music Soothing mode: air pressure + hot compress + music
- Eye care mode: air pressure + vibration + music Sleep mode: hot compress + music Relaxation mode: air pressure + music
- After 15 minutes working the massager will turn off with the voice prompts "Thanks for using, goodbye". 2. The default volume for booting is level 3. Short press ◀ to
- switch to 3--4--5--off--1--2--3... so on and so forth. 3.By default, the music and the Bluetooth is on when booting.
- The music is on and the Bluetooth is on by default when booting. Short press ♪ to switch music. Long press to pause
- 4.The Bluetooth name of the eye massager is "SKG-E3-EN". For first time use, please connect the audio equipment with Bluetooth. Afterwards, Bluetooth will be matched automatically when the audio equipment is on. When the Bluetooth is connected, the voice prompts "Bluetooth is connected". When the Bluetooth is disconnected, the voice prompts "Bluetooth is disconnected"

5.Before wearing, pleas remove contact lens and hair pins, and

adjust elastic bands to the most comfortable position.

- Use clean water and wet cloth to wipe the dirty surface. Don't use alcohol and other volatile liquid to wipe. The product cannot be dissembled and put in water for washing. • For long-term idleness, cover the product to avoid dust.
- During cleaning process, prevent water or other liquids from entering the machine and avoid short circuit.

Failure	Possible cause	Treatment method	
Unable to boot	Low power The product is being charged	Charge the massager fully Use after charging	
Unable to charge or no charging indication	The USB cable, massager and power adapter are not well-connected	Re-connect after unplugging	
Shut down after short-term usage after being fully charged	End of battery life	Contact SKG customer service for replacement	

- ∧ Notes
- Make sure to carefully read and follow the operation guidelines of the user's manual before using the product. The company is not held liable for any consequence arising from incorrect usage or failure to abide by the user's manual.
- Don't use the product after wearing make-ups as the product may mess your make-up and cosmetics may pollute the product surface. In case the product surface is polluted, wipe it with clean water.

- Don't use the eye massager during driving.
- Don't use the eye massager at places subject to excessively
- Don't use the product if you are patients suffering from stroke,
- Patients who had eye surgery and heart disease are not allowed
- The product is personal care massage product and cannot replace
- Children, pregnant women, women in menstrual period are not
- high-temperature environment.

Marning

- high temperature or excessive humidity.
- high blood pressure and low blood pressure.
- Don't use the product if you have eye symptoms such as
- glaucoma, cataract, amotio retinae, pinkeye. Avoid using the eye massager after drinking.
- to use the product.
- Make sure you are not wearing contact lens and have all hair pins and accessories removed.
- medicines and medical devices.
- allowed to use eye massager.
- As the product contains lithium battery, don't place it in



Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received,

This device complies with part 15 of the FCC Rules.

Users subject to sensitive skins, facial allergic dermatitis,

• Immediately stop using the product in case of discomfort

Don't dissemble the product without authorization.

Don't massage parts other than your eyes with the product.

the guidance of doctors.

during usage.

should be careful when using it.

plastic materials in faces are advised to use the product under

• The product can be heated. Users who are not sensitive to heat

including interference that may cause undesired operation. Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void

the user authority to operate the equipment. Notes: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the f



ollowing measures:





Date of manufacture: _

Inspector No.:

Children should be supervised to ensure that they do not play with the appliance. The appliace must only be supplied at safety extra low voltage

--Increase the separation between the equipment and receiver.

--Consult the dealer or an experienced radio/TV technician for help. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a

--Connect the equipment into an outlet on a circuit different

--Reorient or relocate the receiving antenna.


from that to which the receiver is connected.

person responsible for their safety.

corresponding to the marking on the appliance. Do not expose the appliance or battery to excessive temperatures. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.

must be careful when using the appliance.



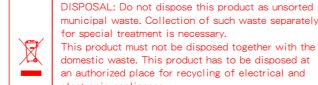

municipal waste. Collection of such waste separately for special treatment is necessary. This product must not be disposed together with the domestic waste. This product has to be disposed at an authorized place for recycling of electrical and By collecting and recycling waste, you help save

natural resources, and make sure the product is disposed in an environmental friendly and healthy This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This device and its antenna(s) must not be co-located or operating in conjunction with any other antenna or transmitter.

rear parameters				
roduct name	Eye Massager	Product code	1205060008	
em code	36833	Rated input voltage	5V==	
1odel	E3-EN	Rated input power	5W	
attery capacity	1100mAh			
tandards	IEC 60335-1	IEC 60335-2-32	2	



