
HubiCare

Heating and Cooling with 3D Air Compression

EYE MASSAGER

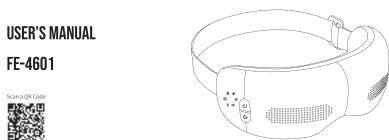


Statement of Responsibility

- 01. Our company reserves the right to modify any information in this User's Guide without prior notice and without any responsibility.
- 02. Unless required by applicable laws, our company does not provide any express or implied warranty for the contents of this User's Guide, including but limited to the warranty of sale or use for a specific purpose.
- 03. Our company will not bear any losses arising from personal and property injuries or other social hazards caused by your illegal operation of this product.

Dear User:

Preface

You are welcome to use the products of our company product. In order to fully experience our products, you must carefully read and save this user guide.

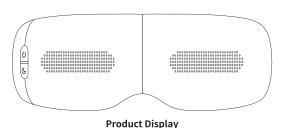

Attention

⚠ Hazard warning:

Do not use this product under any of the following circumstances

- a. For those who wear contact lenses, remove the contact lenses before use.
- b. Patients in the rehabilitation stage after eye surgery.
- C. Patients who are receiving treatment and feel abnormal.
- d. Patients with hypertension, high intraocular pressure, arteriosclerosis or cerebral thrombosis.

7.34°C & 38°C two levels heat selection

8. 22°C cold feeling selection

10. Low power tips

11. Type-C jack

9. Two button 7 modes selection

- e. Patients with high fever, malignant tumor and osteoporosis.
- F. Patients with infectious diseases.

⚠ be careful

Attention

- 01. Before using or cleaning the product, be sure to pull the plug out of the socket in advance.
- 02. Please keep the product away from water, fire and corrosive environment. 03. If the product accidentally falls into the water during charging, please pull out the plug
- immediately. Do not directly touch the product falling into the water. 04. Please do not use this product in humid, high temperature and high electromagnetic
- environment.
- 05. This product contains lithium battery. Please do not disassemble, modify or throw it into fire to avoid danger.
- 06. This product has the function of constant temperature heating. If it is used for a long time or is sensitive to temperature, it may cause low-temperature burns.
- 07. Do not use this product continuously for a long time, or it may cause fatigue and discomfort in the massage area.
- 08. This product is not applicable to any patient or person without operation ability. 09. Do not use this product with children or persons with limited ability without adult care.
- 10. This product is a household electronic product, which is not used as a medical device and is not applicable to diagnosis, treatment and other commercial purposes.

Packaging Instructions

Before using the product, please check the integrity of the internal accessories. In case of any missing or damaged parts, please contact the seller and replace them with new ones.

The list of items in the package is as follows:

01. Eye massager * 1

02. Type-C charging cable * 1

03. User Manual * 1

Preparation before use

Please scan the OR code for

user manual in preferred languages

Please confirm whether the product is powered. If the power is too low, and the machine function is turned off or cannot be started, please charge it before use.

- 1. Connect the USB charging cable to the charging hole of the product.
- 2. The red indicator blinks during charging, and fully lights on blue after charging.
- 3. During the charging process, the product does not support power on. If the product is powered on for charging, the product will automatically shut down.
- 4. The host adopts built-in rechargeable battery, and the charging time is 3.5-4 hours.
- 5. In the process of use, voice prompts the low battery alarm, please pause the use, and use it after full charge.

Preparation before use

Please confirm whether the product is powered. If the power is too low, and the machine function is turned off or cannot be started, please charge it before use.

Multiple Charging Methods


RF Exposure Information

Computer USB interface/mobile power/5V==1A power adapter

Note: The default working time of power on is 10 minutes.

If you need to extend the service time, you need to turn on the host again.

Operation Guide

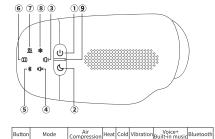

- U Relax Mode V V V

- (Auto Mode / Relax Mode / Relieve Mode / Energy Mode)
- Power-up with "Low Level Heat" with White light
- Short Press: Adjust 3 sleep modes
- (Soothing Mode / Sleep Mode / Cod Mode)
- Double-press quickly: Go to next music
- 4. Speaker Mute 5 Bluetooth
- Air pressure function

1. Power Button U

- Long Press: ON/OFF

- 9. Battery: Low Power<15% Red LED blinking
- When Full Charging Blue LED Lighting

Product Function

- 1. Voicebroadcast
- 3. Bluetooth function selection
- 4. Can adjust to next song function
- 5. Speaker mute selection
- 6. Multiple Vibration Function

Transportation and storage specifications

Atmospheric pressure: 86kpa~106kpa Ambient temperature: - 10°C~65°C

Relative humidity: <80% RH, no corrosive gas and well ventilated room.

Transportation: This product is equipped with simple shock proof facilities, which is suitable for air, railway, highway and ship transportation, and should avoid rain and snow splashing, inversion and collision. If the storage period exceeds 6 months, open the box for inspection and power on for half an hour. The room shall be well ventilated to avoid the corrosion of smn corrosive gas in strong sunlight.

Product Specifications

Product: Heating And Cooling With 3D Air Compression Eye Massager

Model: FE-4601

Rated Power: Input 5V 1A / Output DC 3.7V

Battery: 3.7V, 2500mAh

Manufacturer: Xiamen Healthpal Electronic Co,.Ltd (CHINA)

Charging: USB C-type

Weight: 430g

Size: L210 x W115 x H80(mm)

Distributed by HuBDIC USA

Carson, CA 90246 United States of America Customer service center: 424-488-2241

Troubleshooting and maintenance

1. If any fault occurs during operation, please conduct the following inspection.

2. If the basic fault of the product cannot be solved after troubleshooting, please call the after-sales customer service personnel for consultation. Do not disassemble this product by

	Basic Fault	Possible Causes	Basic elimination methods
	Unable To Start	No power	Use After Charging
	Air Pressure Cannot Be Used	Enterprise Transferred To Closed Status	Adjust The Air Pressure Mode To Another Mode
	NO COLD	Working With No Cold Or The Use Time Is Too Long	Please manually adjust to the ice coating mode or shut down the product for half an hour before use
	NO HEAT	Working With No HEAT	Manually Adjust To The Heating Mode
	Unable To Charge	USB Charging Cable Is Not Inserted Properly	Insert USB charging Cable Corretly

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- · Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two

conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The device has been evaluated to meet general RF exposure requirement.

The device can be used in portable exposure condition, compliance with exposure requirements.

- 2. Built-in 4 inner music
- Short Press: Adjust 4 Modes
- Double-press guickly:
- Adjust to "High Level Heat" with Blue light
- 2. Sleep Button 🕓
- Long Press: Speaker mute
- 3 Vibration
- 8. Cold
- When Charging Blue LED blinking