

TENS NECK MASSAGER

SHARPER IMAGE®



RECHARGEABLE BATTERY WARNINGS:

- This product contains batteries that are non-replaceable.
- There is a risk of fire and personal injury if battery is punctured, damaged, or misused.
- Do not dismantle, crush, incinerate or short-circuit battery.
- Never expose battery to extreme temperatures or moisture.
- Keep away from flammable materials and direct sunlight.
- Only charge using the provided charger, cable, and/or adapters.
- Rechargeable battery is only to be charged by an adult.
- Adult supervision is required while charging the battery.
- Allow battery to cool to room temperature before charging.
- Charge on a non-flammable surface and keep away from flammable products.
- Discontinue charging if battery becomes hot, swells, smokes, or gives off a strong odor.
- Battery chargers used with the product are to be regularly examined for damage to the cord, plug, enclosure, and other parts. In the event of such damage, chargers must not be used until the damage has been repaired.
- Recycle or dispose of battery according to federal, state, and local laws.
- Keep these instructions for future reference.

Keep this instruction manual as it contains important information for future reference.

Due to continuous product improvements, the picture on the manual may differ slightly from the actual product.

SHARPER IMAGE® NAME AND LOGOS (INCLUDING THE 77 LOGO AND TOMORROW'S TOMORROW) ARE TRADEMARKS OR REGISTERED TRADEMARKS. DISTRIBUTED BY AND © 2021 MERCHSOURCE, LLC, IRVINE, CA 92618. ALL RIGHTS RESERVED. FOR CUSTOMER SUPPORT, PLEASE CALL 1-800-374-2744.

Item # 1014706

Printed in China

CONGRATULATIONS ON YOUR PURCHASE OF THE TENS NECK MASSAGER

WARNING: Contains button or coin cell battery. Hazardous if swallowed—see instructions.

- WARNING:** This product contains a Button or Coin Cell Battery. A swallowed Button or Coin Cell Battery can cause internal chemical burns in as little as two hours and lead to death. Dispose of used batteries immediately. Keep new and used batteries away from children. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.
- Do not short-circuit the battery supply terminals
 - Do not dispose of batteries in fire, batteries may explode or leak.
 - Recycle or dispose of batteries according to federal, state, and local laws.
 - Keep these instructions for future reference.

CARE AND TIPS

- Power down before cleaning. Clean with a 75% alcohol wipe.
- Do not use detergent or soap when cleaning.
- Do not submerge in water.
- To access the remote battery, use a coin to remove the battery cover.

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, this is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

CONTENTS

- TENS Neck Massager
- Wireless Remote

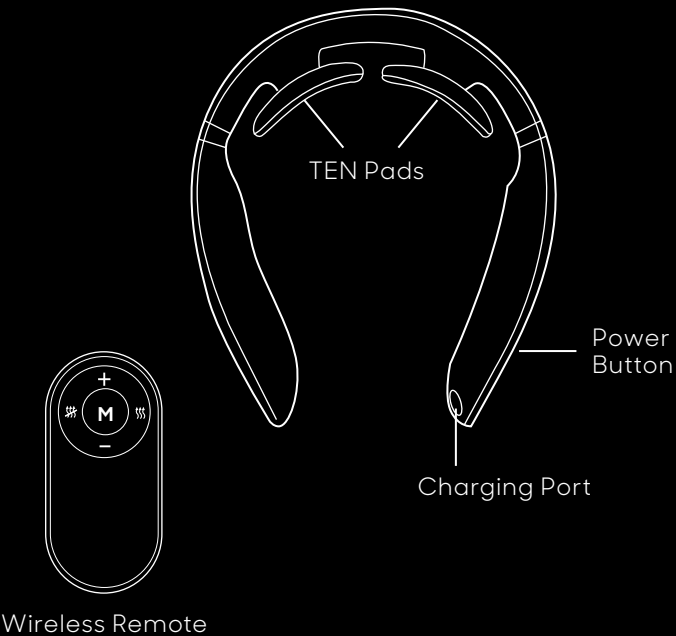


FEATURES

- TENS massage with heat mode
- 3 massage modes and 15 strength levels
- Wireless remote with magnetic mounting

GETTING STARTED

1. Charge using the provided USB cable and unplug when fully charged.
2. Secure the rubber charging port cover.
3. Remove any electronic items or jewelry you are wearing.
4. Using a clean, damp towel, moisten the sides and back of your neck.
5. Place the massager around your neck. Make sure the TENS pads are in full contact with your skin.
6. Turn on the massager by pressing the power button. A tone will be heard when massager is turns on.
7. Find your comfort level. Press the [+] button to increase massage intensity and the [-] button to decrease massage intensity.
Note: Massager will turn off if left at the lowest setting.
8. Press the [M] button to toggle through the 3 massage modes.
Note: changing modes resets massage intensity to lowest level.
9. Press the [Speaker Icon] button to turn heat on and press the [Speaker Icon] button to turn heat off.
10. Press and hold the [-] button turn off sound cues. Press and hold the [+] button to turn on sound cues.
11. Use for up to 15 minutes at a time but no more than 30 minutes per day.

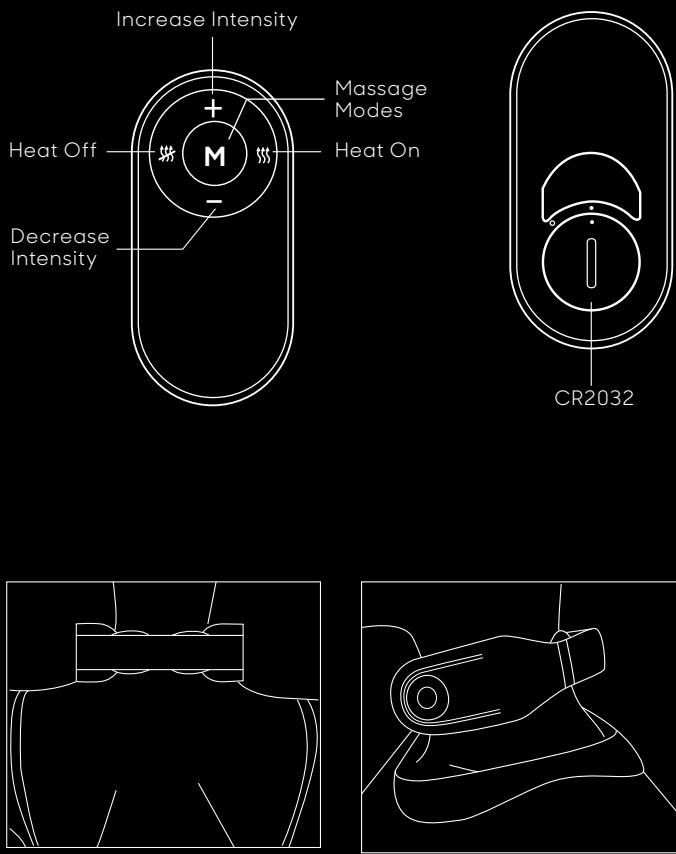

Warning: Use massager only as indicated. Do not use massager for any other body areas.