

Welcome

Thank you for purchasing our product. We recommend you use it with the supporting app to learn more about teeth brushing.

Important Safety Tips

It is recommended that this toothbrush be used only by a person over 12 years old.

- It is normal to feel slightly uncomfortable when the toothbrush is used for the first time.

- The toothbrush can only be used with Oclean special heads. We will not bear any responsibility due to the use of any third-party accessory.


This toothbrush is intended for personal brushing only. It is not recommended that the toothbrush be shared with others.

- Avoid exposing the toothbrush to direct sunlight.
- Please power off the charging holder when it is not in use and before cleaning it.

Unrecommended Use

Persons who have just undergone oral surgery for treatment. It is recommended that the consumer who has mouth disease, such as seriously hemorrhaged gums, consult a dentist before using this toothbrush.

- Persons who have mental instability, those who are disabled and those who lack relevant experience unless the product is used under supervision.

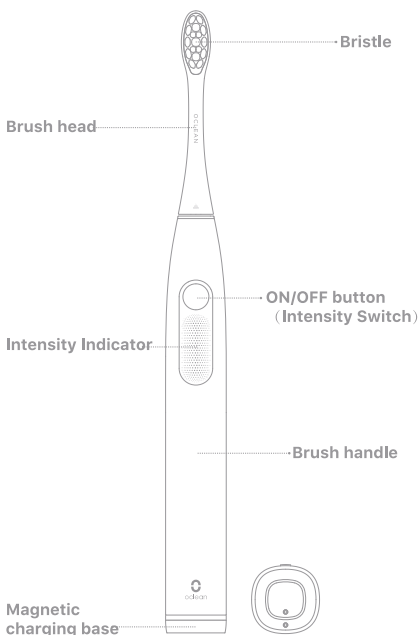

Instructions

• First-Time Use

The toothbrush can be awakened by pressing the ON/OFF button for the first time on the brand-new toothbrush. (The Bluetooth function will be activated and remain on.) After installing the Oclean App, please enter the app to complete the connection and settings of the toothbrush.

• Install the App and Toothbrush Settings

You can use the app to customize your brushing programs. It will then complete the process of tooth cleaning, tooth whitening and gum massaging automatically in one brushing session. There is no need to open the Oclean app in the daily brushing process. Instead, simply open the app every few days in order to sync the brushing data.

• Direct Use Without App Settings

Without the toothbrush settings in Oclean App, the toothbrush will default to run for 2 minutes as a "standard cleaning plan", in which every 30 seconds the toothbrush will remind the user of the need to change the brush area by vibration, and will automatically stop running after 2 minutes of working.

• Switch Intensity position during operation

Long press (2 seconds) to change the intensity. Different numbers of lights represent the high, medium and low intensity.

• Change Brushing Area

When the toothbrush is working, there will be a reminder with vibration to remind the user to change the brushing area every 30 seconds.

• Intensity

When the toothbrush is working, the roughly strength range can be judged according to the length of the LED indicator. When only one light is on, means the product is running in a lower intensity, and when all three lights are on, means it's running with a higher intensity.



Low Intensity
Cleaning Level 1-11
Whitening Level 1
Massage Level 1
Sensitive Level 1



Medium Intensity
Cleaning Level 12-24
Whitening Level 2
Massage Level 2



High Intensity
Cleaning Level 25-32
Whitening Level 3
Massage Level 3

• Low Battery

The indicator flashes when battery percentage is lower than 20%, needs to be charged to ensure normal use.

• Charging Percentage

When the toothbrush charging, the charging indicator shows a breathing flashing state. The charging percentage can be judged according to the length of the indicator. The red indicator flashes when battery percentage is lower than 20%. One white indicator flashes when battery percentage is around 20-40%. Two white indicator flashes when battery percentage is around 40-70%. Three white indicator flashes when battery percentage is around 70-100%. Three white indicator is on constantly when battery is full.

• Brushing Plan and App Settings

You can either choose or customize a brushing program through the supporting Oclean app. You can also upgrade the firmware to continuously update its performance.

How to Charge

Place the toothbrush on the charging base, the intensity indicator starts to show the charging progress. It can be fully charged in about 2 hours. (The recommended parameters of charger are between 5V/1A-5V/2A).

Important Charging Tips

The toothbrush runs with a fast-charge mechanism. It is recommended that the toothbrush be placed on the charging base only when it is necessary. Please don't put the toothbrush on the charging base for a long time, as it may cause the built-in battery to decay prematurely so that the user experience is affected.




1. No disassembly, impact, extrusion or throwing into fire
2. Please stop using the toothbrush if the battery has expanded
3. Please don't put the battery in high temperature environment
4. Do not use the battery if water entered the device.


Recommended Tooth Brushing Method

• Tooth Brushing Time


Dental professionals recommend that each brushing process should last two to three minutes.

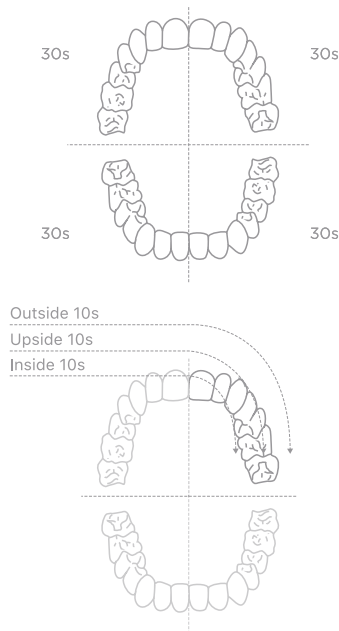

Recommended Positioning for Different Areas

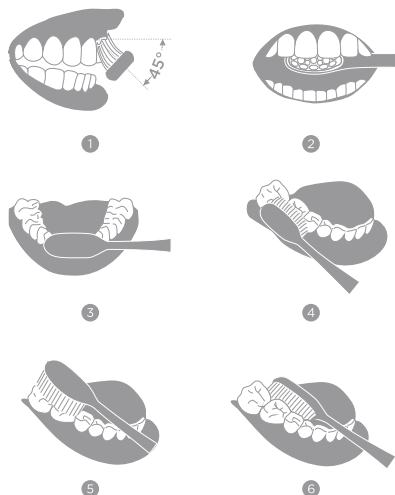


Cleaning and Maintenance

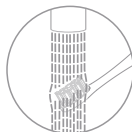

Cleaning Brush Head

Clean the brush head after each use and air-dry it naturally.

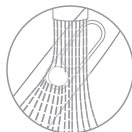


Cleaning the Body

Clean the body with a soft cloth.



Cleaning Charger

Clean the dirty parts of the charger with a soft cloth. Be sure to remove the battery plug from the socket prior to cleaning the charger, and use it again only when it is fully dry.

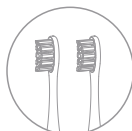

Changing the Head

It is recommended to change the brush head every two to three months.




Cleaning the Handle Base

Always keep the bottom clean and dry before charging. Clean the dirty parts of the handle with a soft cloth.

• Steps




Insert the brush head carefully onto the body of the toothbrush. A small gap should be left between the head and the main body so the brush head can function properly.

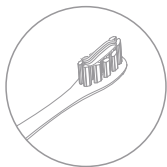

To prevent toothpaste from splashing, smear a small amount of toothpaste on the head before using it. (A sonic electric toothbrush conserves more toothpaste than a manual toothbrush.)




Soak the toothbrush head and put it into your mouth, then turn on the toothbrush to start brushing.

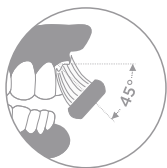

Slowly move the bristles at 45° to the gum line to brush your teeth properly.



Hold the toothbrush handle and slowly move it in a circular motion. Do not press too hard.

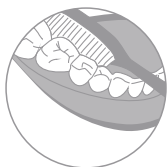

Please rinse with clean water after brushing.

• Specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Max Torsion	220gf.cm
Frequency	2.4GHz
Waterproof	IPX7
Power	1.8W
Charger	USB charger
Input	5V \pm 1A
Battery Life*	30 days
Battery Capacity	800 mAh
Charging Time	2 hours

* Calculated by 2 minutes/time, twice/day.

• FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

• IC WARNING

This device contains licence-exempt transmitter (s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

(1) This device may not cause interference.
(2) This device must accept any interference, including interference that may cause undesired operation of the device.

Cet appareil contient des émetteurs / récepteurs exemptés de licence conformes aux RSS (RSS) d'Innovation, Sciences et Développement économique Canada. Le fonctionnement est soumis aux deux conditions suivantes :

(1) Cet appareil ne doit pas causer d'interférences.

(2) Cet appareil doit accepter toutes les interférences, y compris celles susceptibles de provoquer un fonctionnement indésirable de l'appareil.

CAN ICES-3(B)/NMB-3(B)

Facebook Page



www.facebook.com/ocleanofficial

Follow the
Oclean Official Facebook Page

Download App



www.oclean.com/download

If the product you have purchased is found to have a fault or defect, please contact us.

North America

Call us at - (866)OCLEAN9 or (866) 625-3269

Email: support.us@oclean.com

Other Region

Email: service@oclean.com

Manufacturer: Shenzhen Yunding Information Technology Co., Ltd.

Manufacturer's address: #32G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Production address: #6-6-101, Yanshan Avenue, Yanchuan Community, Yanluo Street, Bao'an District, Shenzhen



Hereby, Shenzhen Yunding Information Technology Co., Ltd. declares that the radio equipment type Y3 is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address:
<http://www.oclean.com/en/service.aspx>



R 201-190289



For more regulatory information and product certification/compliance logos related to Oclean X, please go to Settings->About->Regulatory on your toothbrush.