
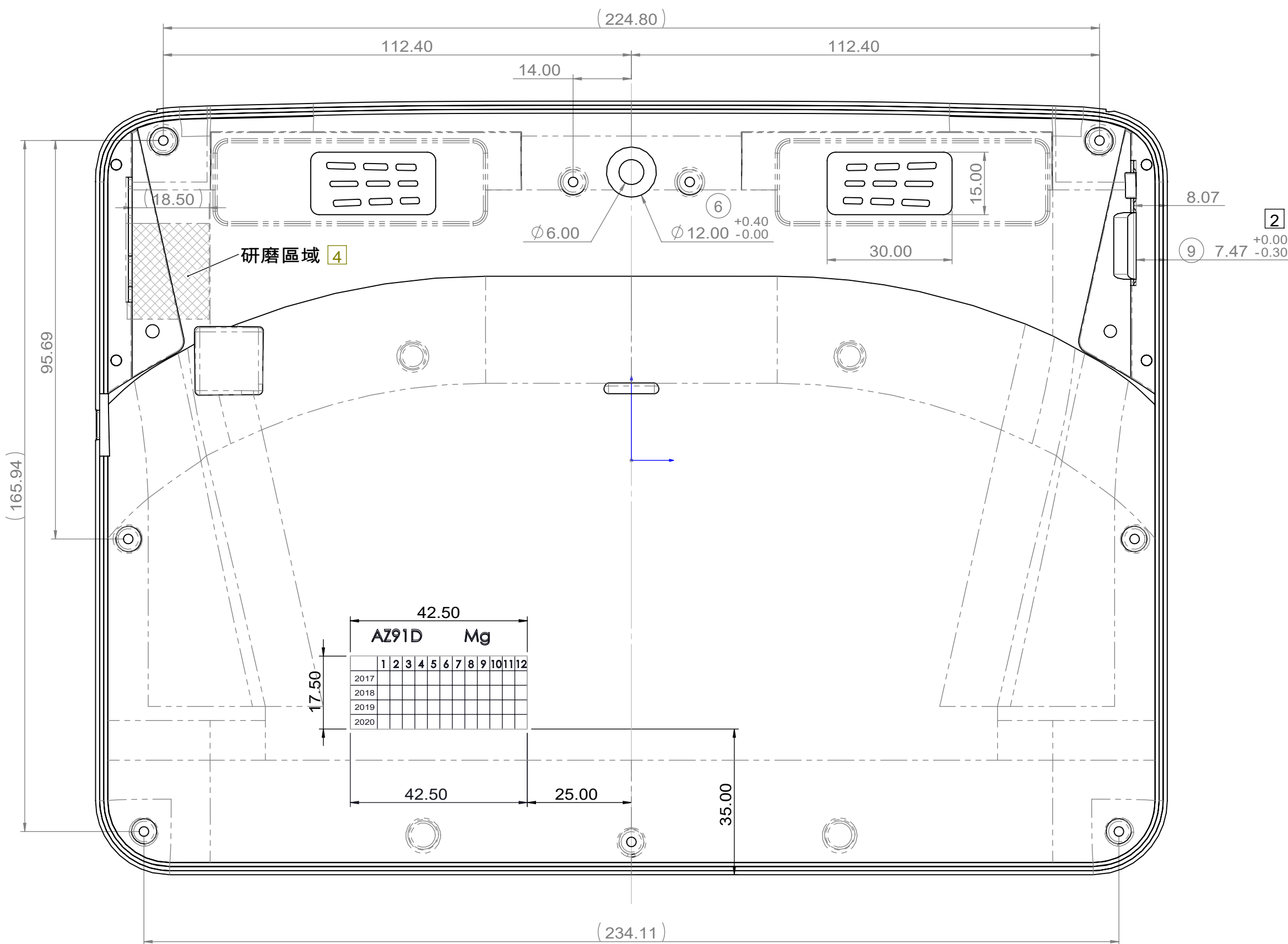
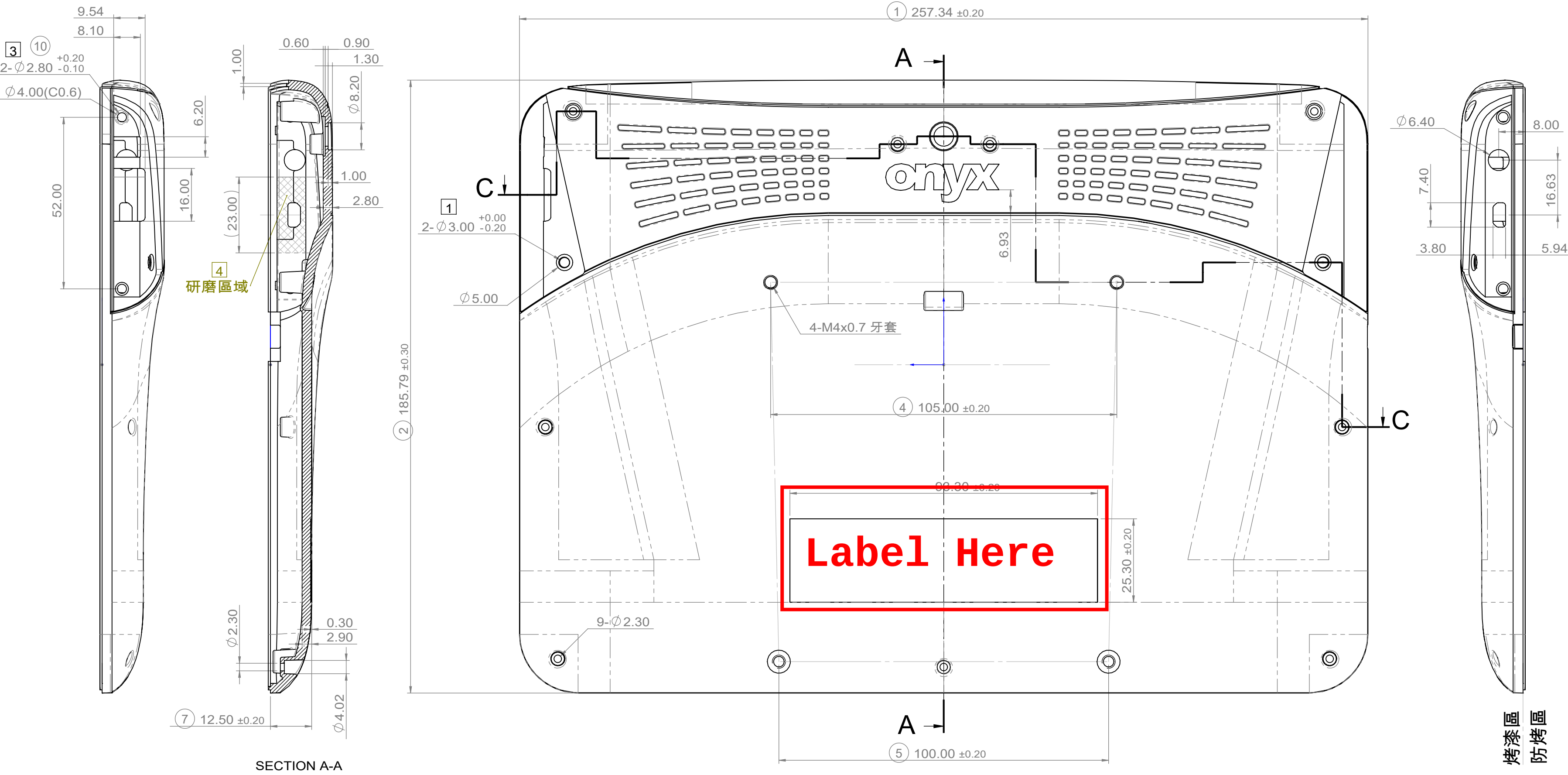
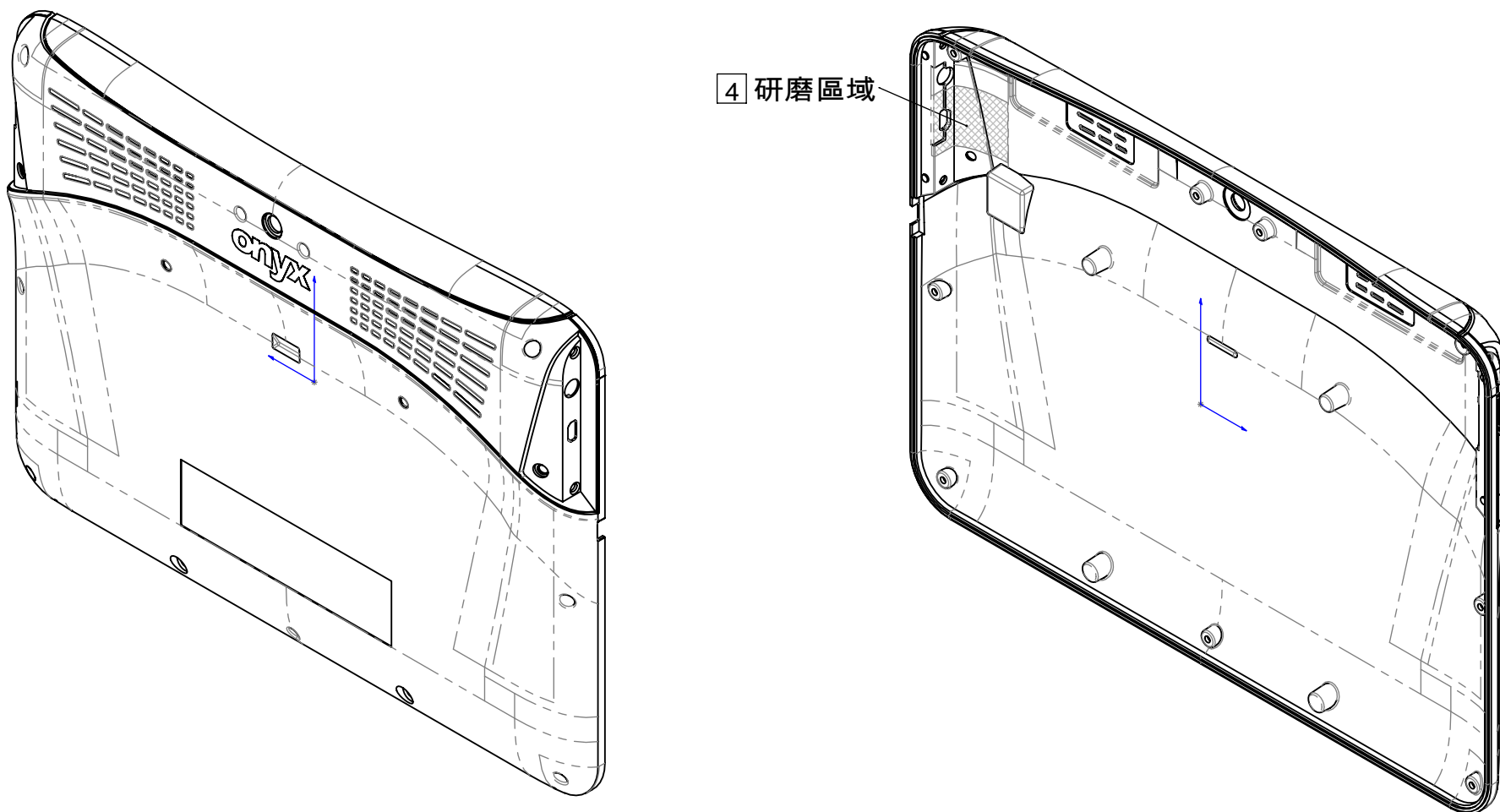


NOTE:

- 1.No sharp edges/ chips,no sink mark, injection mark, flow line,flash and welding line should be visible on appearance surface.
- 2.Warping should be 0.3mm per 100mm ; part max 0.5mm.
- 3.The geometry dimensions of enclosure is based on 3d model
- 4.Surface finish:Don't have shrink marks and scratches
- 5.This drawing contains critical to function dimension only, refer to the 3d-cad file for complete part information.
- 6.Parts should meet RoHS.
- 7.The interior of the material should be added to the date code and the recycling symbol.

- 1 孔徑由2.2-->3.0mm(2017/07/06)
- 2 增加尺寸7.47+0.0/-0.3mm(2018/09/21)
- 3 孔徑增加公差由2.8-->2.8+0.2/-0.1mm(2019/07/06)
- 4 去除皮膜(2022/06/16)

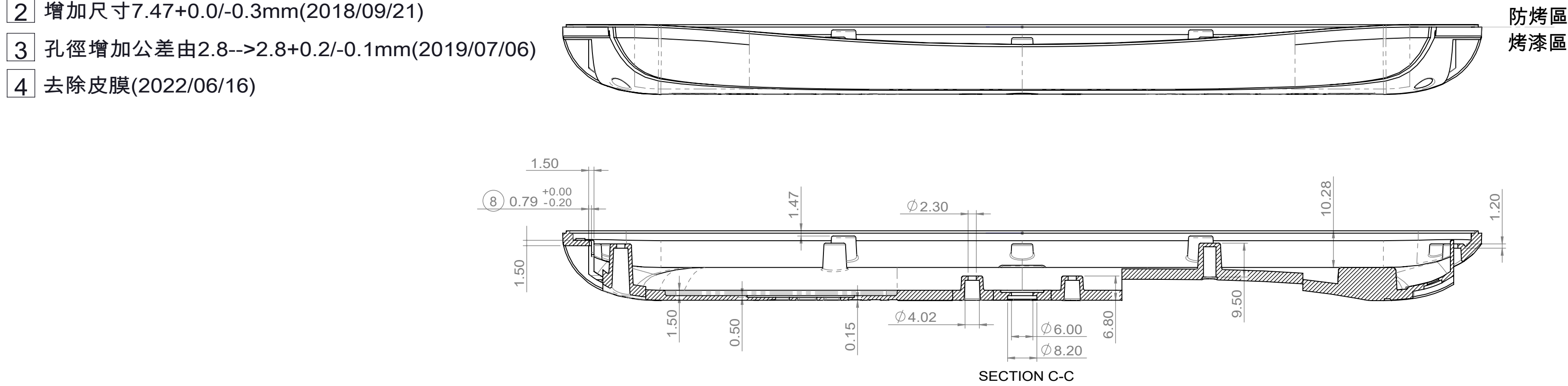

Paint color : PMS Cool Gray 10C H00

Print color : PMS# 7468C

GENERAL TOL. ≥0 ±0.1 ≥10 ±0.2 ≥50 ±0.3 ≥200 ±0.5 ≥500 ±0.8 ANGLE TOL. ±0.5°	MATERIAL: Mg-Al		PART No. M94M101030		onyx Healthcare Inc. 			
	SPEC.: AZ91D		DWG No. M94M101030		醫揚科技股份有限公司			
	FINISH: PMS Cool Gray 10C H00		APPROVED	CHECKED	DESIGNED	NAME Rear Cover for Piccolo		
	UNIT: MM		Ernest Johnny		MODEL No. MD101			
SCALE: 1/1					2017 / 04 / 18	2017 / 04 / 18	REV. 0.4	SHEET: 1 / 1