

Medtronic MEDTRONIC CONFIDENTIAL – THIS DOCUMENT IS THE PROPERTY OF MEDTRONIC, INC., AND MUST BE ACCOUNTED FOR. INFORMATION HEREON IS CONFIDENTIAL. DO NOT REPRODUCE IT, REVEAL IT TO UNAUTHORIZED PERSONS, OR SEND IT OUTSIDE MEDTRONIC WITHOUT PROPER AUTHORIZATION.		AUTHOR
		Luke Munroe
REVISIONS		
REV	DESCRIPTION	DATE
A	Initial release	10 May 2021

SPEC NUMBER	REV	SHEET
B05631	A	1 of 5

1. This specification is a representation of the final Label Design Output generated from a combination of label parts within the associated PDM BOM.
2. These label images are the combinations of the Transform template file, LDS (Label Model Data) record, sample trigger data and representative views of the label die lines.
3. The Serial Number, Use-By Date, Lot Number, Date Of Manufacture, PIN Number, GTIN, GTIN Version (if used), Product Number (REF) and Barcodes are shown for example only. These fields are variable and are to be supplied by Manufacturing at the time of printing. Not all fields will be used for all labels.
4. The Label Set is to be printed using a 600 DPI LaserJet printer or equivalent. The individual labels and PIC are to be printed using a 300 DPI Zebra thermal transfer printer.
5. The Trigger Data shown on the Label Set Image might differ from the printed label.
6. Spacing and text wrapping shown in the image may differ slightly from actual printed label. Content of printed version is accurate.
7. The label printed in Manufacturing must be at an equal or later Model Data and Template revision as shown on the image.
8. This specification is controlled by Neuromodulation Package Engineering.
9. This specification has been assessed for key design characteristics which are defined as
 - Critical KC_C - Identified characteristic is classified as Critical and impacts safety.
 - Major KC_M - Identified characteristic is classified as Major and impacts safety and/or performance.

Medtronic MEDTRONIC, MINNEAPOLIS, MINNESOTA USA	TITLE LABEL IMAGE SPECIFICATION NEURO 97800 CE MARK				
	IMAGE SPEC TEMPLATE: A86694_c				
	DATE 10-May-2021	SCALE N/A	ENGINEERING Luke Munroe		
	MODEL 97800	Size B	SPEC NUMBER B05631	REV A	SHEET 2 of 5

 **Medtronic**

97800

InterStim™ X

SN WWW000000W

Implantation date




 **FOR YOUR PATIENT FILES**

 **MR Conditional**

 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD
 Medtronic 97800		 SN WWW000000W YYYY-MM-DD



InterStim™ X

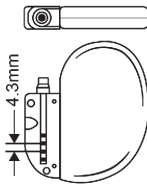


MR
Conditional



Sacral
Neuromodulation

REF **97800**

SN **WWW000000W**

 **YYYY-MM-DD**

 **YYYY-MM-DD**

PIN **W000WWW00W**

MD




www.medtronic.com/manuals

STERILE EO

CE0123

Rx Only







Package contents

Implantable Neurostimulator

Accessories

Product documentation

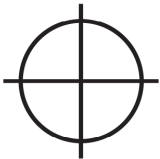

Medtronic Inc., 710 Medtronic Parkway
 Minneapolis, MN 55432 USA

Medtronic B.V., Earl Bakkenstraat 10
 6422 PJ Heerlen, The Netherlands


www.medtronic.com

Manufactured in:
 Juncos, Puerto Rico, USA

Country of Origin:
 USA

EC REP

MO14347C001 B - M020161C001 A

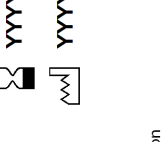	InterStim™ X	REF 97800 SN WWW000000W YYY-MM-DD YYY-MM-DD YYY-MM-DD YYY-MM-DD
 	 Sacral Neuromodulation	 
	  Implantable Package contents Accessories	
 Rx Only	 Medtronic, Inc., 710 Medtronic Parkway Minneapolis, MN 55432 USA www.medtronic.com	M014347C001 B - M020161C001 A
STERILE PACK LABEL		M014347C001 B - M020161C001 A

Medtronic

InterStim™ X

REF 97800

SN WWW000000W

YYYY-MM-DD

M014347C001 B -
M020161C001 A

(01)00000000000000(17)YYMMDD(21)WWW000000W

 **Medtronic**

M014347C001 B - M020161C001 A

REF 97800

FCC ID: LF597800

This product complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation








TemplateNumber:	M014347C001_B
ImageSpecID:	B05631
ImageSpecRev:	
SpecName:	B05090_Txt
ImageSpecTemplate:	A8694_c
LabelDefNo:	B05631-007
OrientationLabel:	
ImageSpecTitle:	LABEL IMAGE SPECIFICATION NEURO 97800
Engineering:	CP MARK Luke Munroe
PrinterName:	Prt_Bmp
FM:	0000000000
MatrixID:	9876543210
QuantityOfLabels:	1
DocumentNumber:	M020161C
DocumentRevision:	A
GTIN:	000000000000000
UPNCCC:	01234567891234
URD:	YYNNMMDD
ManufactureDate:	YYMMDD
LotSerial:	YYNNMMDD
LotSerialNumber:	NNNN00000000W
AddressCode:	007
QualificationStatus:	90
i	
BusinessUnit:	NRI
LaserPrinter:	Prt_Bmp
DotPrinter:	Prt_Bmp
ProductNumber:	97800
TemplateType:	NRI Label Templates
ArchiveTrigger:	210510 145131-9664-B05631-M014347C001

UNIVERSAL LABEL SET (DEVICES) Rev. : E

* A 1 4 6 3 2 0 0 1 *	LOT NUMBER
-----------------------	------------

LASER

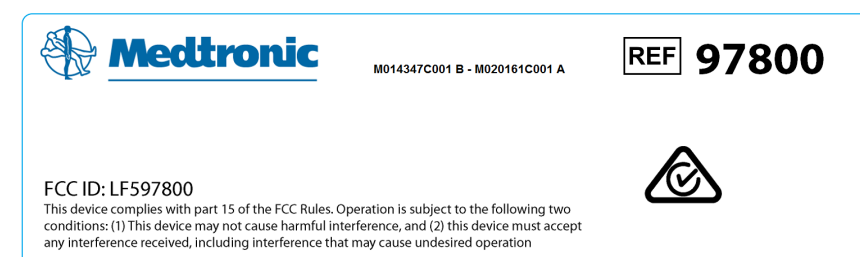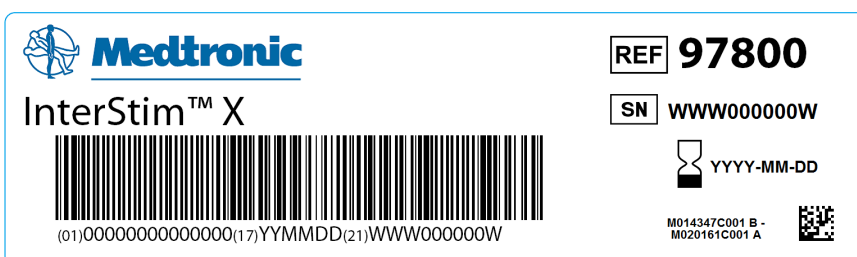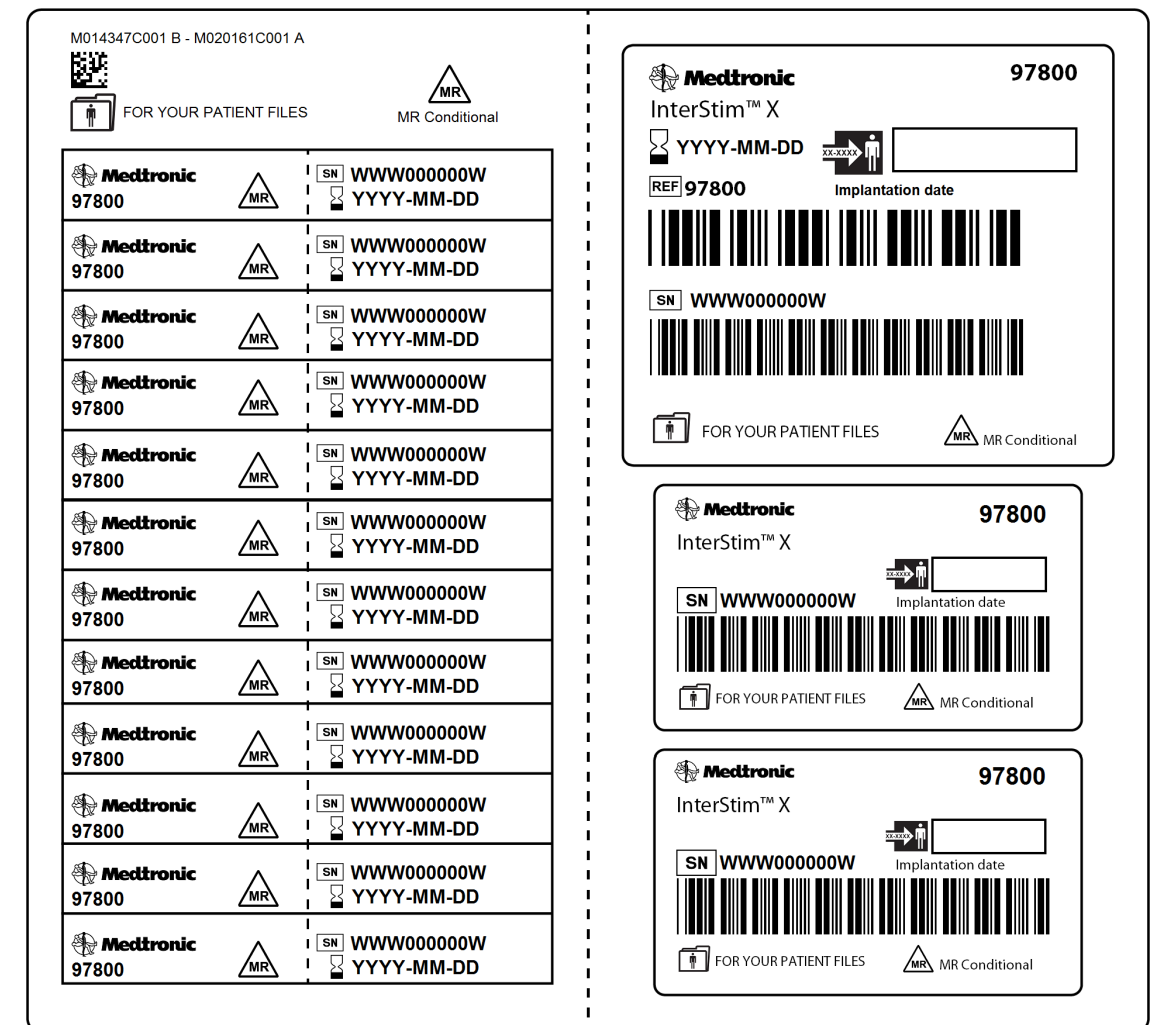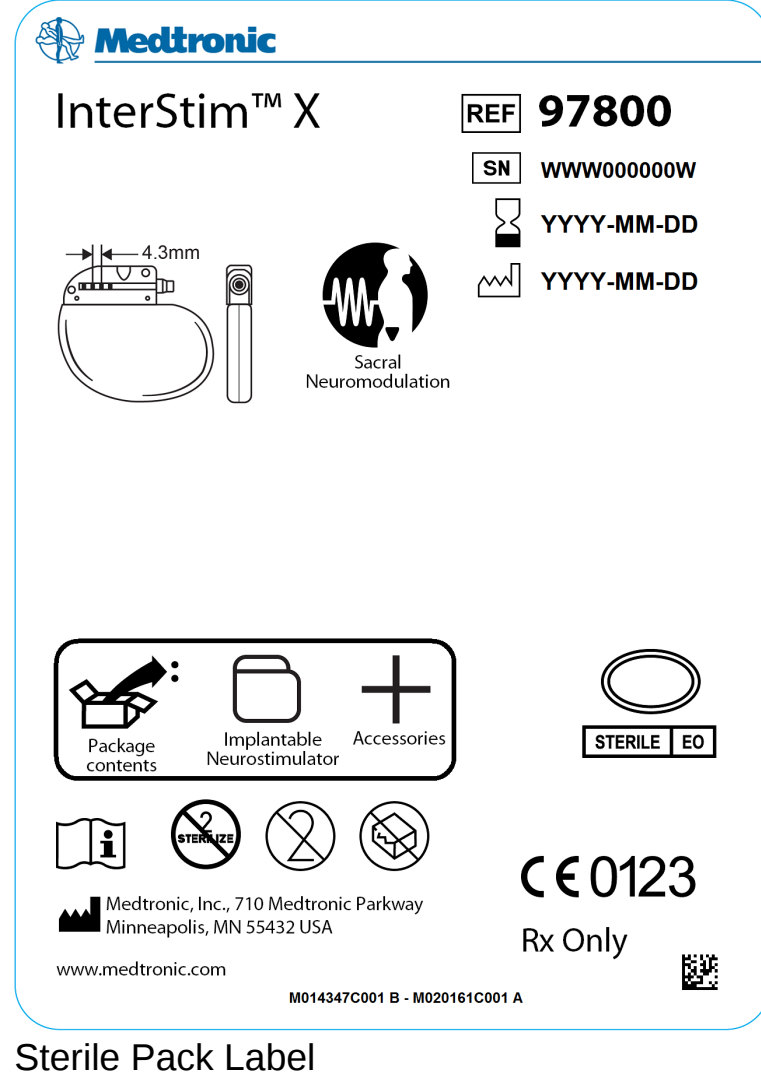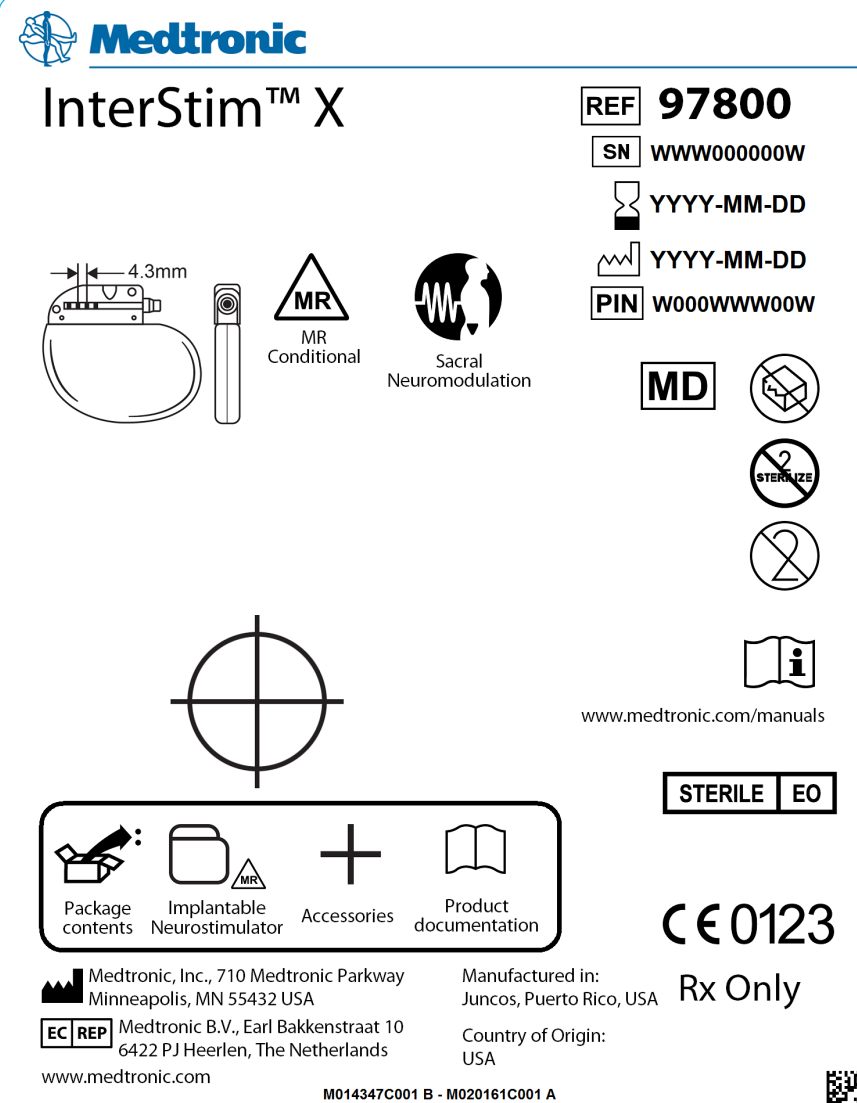

ZEBRA

NOTE: IMAGES NOT TO SCALE

SPEC NUMBER	REV	SHEET
B05631	A	4 of 5

en Stimulator / **fr** Stimulateur / **de** Stimulator / **es** Estimador /
nl Stimulator / **it** Stimolatore / **no** Stimulator / **fi** Stimulaattori /
sv Stimulator / **el** Διεγέρτης / **da** Stimulator / **pt** Estimador /
pl Stymulator / **cs** Stimulátor / **hu** Stimulátor / **sk** Stimulátor /
ro Stimulator / **bg** Стимулатор / **sl** Stimulator / **hr** Stimulator

MD 97800 InterStim™ X

UDI



UDI-DI: (01)00000000000000

SN

WWW000000W



Medtronic, Inc., 710 Medtronic Parkway
Minneapolis, MN 55432 USA, www.medtronic.com

M014347C001 B
M020161C001 A

PIC