
APPENDIX 5 : Dosimetric E-Field Probe Calibration
(ET3DV6,S/N: 1684)

Schmid & Partner Engineering AG

Zeughausstrasse 43, 8004 Zurich, Switzerland, Phone +41 1 245 97 00, Fax +41 1 245 97 79

IMPORTANT NOTICE

USAGE OF PROBES IN ORGANIC SOLVENTS

Diethylene Glycol Monobuthy Ether (the basis for HSL1800 and M1800 liquids), as many other organic solvents, is a very effective softener for synthetic materials. These solvents can cause irreparable damage to certain SPEAG products, except those which are explicitly declared as compliant with organic solvents.

Compatible Probes:

- ET3DV6
- ET3DV6R
- ES3DV2
- ER3DV6
- H3DV6

Important Note for ET3DV6 Probes:

The ET3DV6 probes shall not be exposed to solvents longer than necessary for the measurements and shall be cleaned daily after use with warm water and stored dry.

**Schmid & Partner
Engineering AG**

Zeughausstrasse 43, CH-8004 Zurich
Tel. +41 1 245 97 00, Fax +41 1 245 97 79

Schmid & Partner Engineering AG

Calibration Certificate

Dosimetric E-Field Probe

Type:

ET3DV6

Serial Number:

1684

Place of Calibration:

Zurich

Date of Calibration:

November 20, 2002

Calibration Interval:

12 months

Schmid & Partner Engineering AG hereby certifies, that this device has been calibrated on the date indicated above. The calibration was performed in accordance with specifications and procedures of Schmid & Partner Engineering AG.

Wherever applicable, the standards used in the calibration process are traceable to international standards. In all other cases the standards of the Laboratory for EMF and Microwave Electronics at the Swiss Federal Institute of Technology (ETH) in Zurich, Switzerland have been applied.

Calibrated by:

D. Vetter

Approved by:

Thomas Vetter

Probe ET3DV6

SN:1684

Manufactured:	April 3, 2002
Last calibration:	May 10, 2002
Recalibrated:	November 20, 2002

Calibrated for DASY Systems

(Note: non-compatible with DASY2 system!)

DASY - Parameters of Probe: ET3DV6 SN:1684

Sensitivity in Free Space

NormX	1.55 $\mu\text{V}/(\text{V}/\text{m})^2$
NormY	1.54 $\mu\text{V}/(\text{V}/\text{m})^2$
NormZ	1.58 $\mu\text{V}/(\text{V}/\text{m})^2$

Diode Compression

DCP X	96	mV
DCP Y	96	mV
DCP Z	96	mV

Sensitivity in Tissue Simulating Liquid

Head	900 MHz	$\epsilon_r = 41.5 \pm 5\%$	$\sigma = 0.97 \pm 5\% \text{ mho/m}$
Head	835 MHz	$\epsilon_r = 41.5 \pm 5\%$	$\sigma = 0.90 \pm 5\% \text{ mho/m}$
ConvF X	6.8 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	6.8 $\pm 9.5\%$ (k=2)	Alpha	0.31
ConvF Z	6.8 $\pm 9.5\%$ (k=2)	Depth	2.51
Head	1800 MHz	$\epsilon_r = 40.0 \pm 5\%$	$\sigma = 1.40 \pm 5\% \text{ mho/m}$
Head	1900 MHz	$\epsilon_r = 40.0 \pm 5\%$	$\sigma = 1.40 \pm 5\% \text{ mho/m}$
ConvF X	5.5 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	5.5 $\pm 9.5\%$ (k=2)	Alpha	0.40
ConvF Z	5.5 $\pm 9.5\%$ (k=2)	Depth	2.61

Boundary Effect

Head	900 MHz	Typical SAR gradient: 5 % per mm	
Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	7.8	4.5
SAR _{be} [%]	With Correction Algorithm	0.2	0.4
Head	1800 MHz	Typical SAR gradient: 10 % per mm	
Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	10.6	7.4
SAR _{be} [%]	With Correction Algorithm	0.2	0.2

Sensor Offset


Probe Tip to Sensor Center	2.7	mm
Optical Surface Detection	1.5 $\pm$ 0.2	mm

ET3DV6 SN:1684


November 20, 2002

Receiving Pattern (ϕ), $\theta = 0^\circ$


f = 30 MHz, TEM cell ifi110

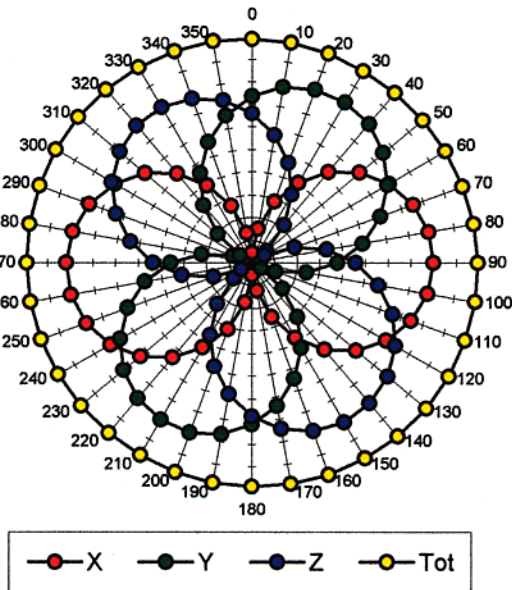


f = 100 MHz, TEM cell ifi110

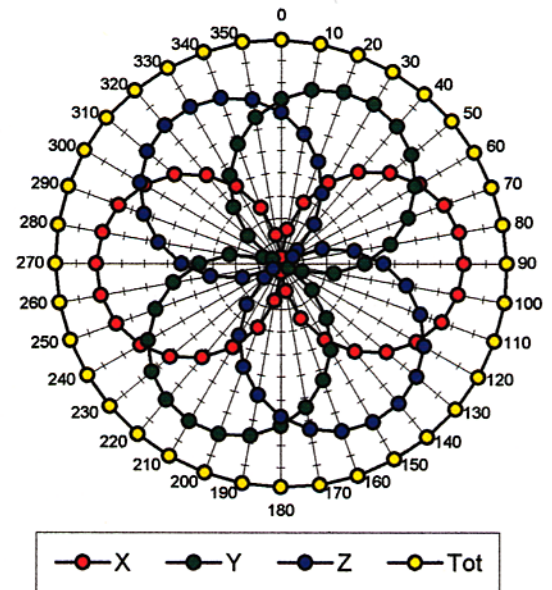

f = 300 MHz, TEM cell ifi110


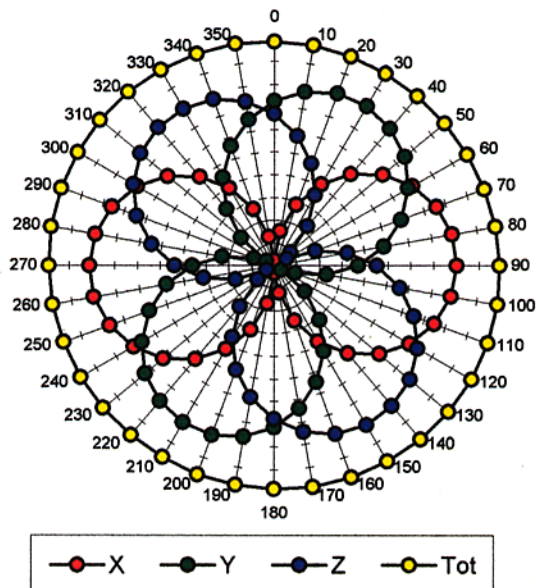


f = 900 MHz, TEM cell ifi110

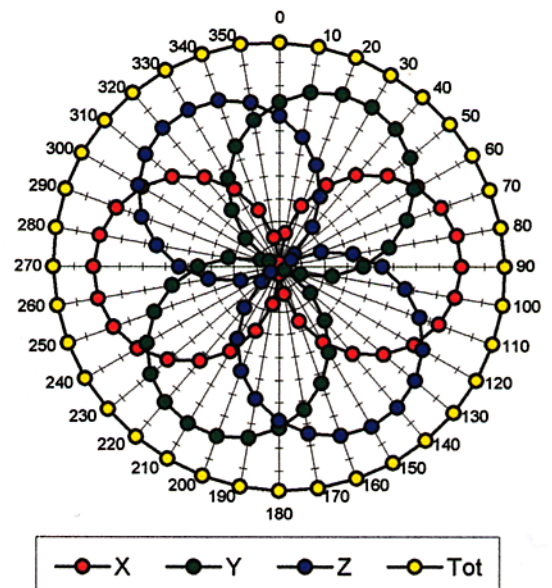


ET3DV6 SN:1684


November 20, 2002

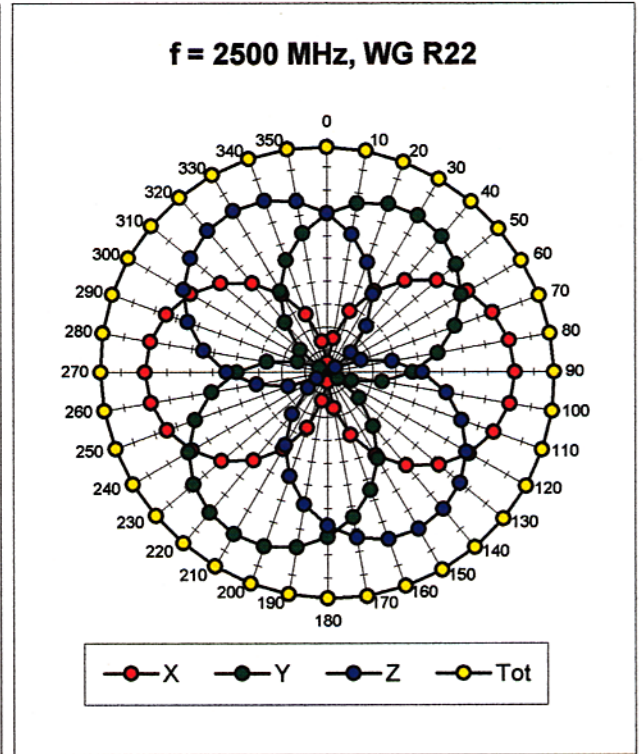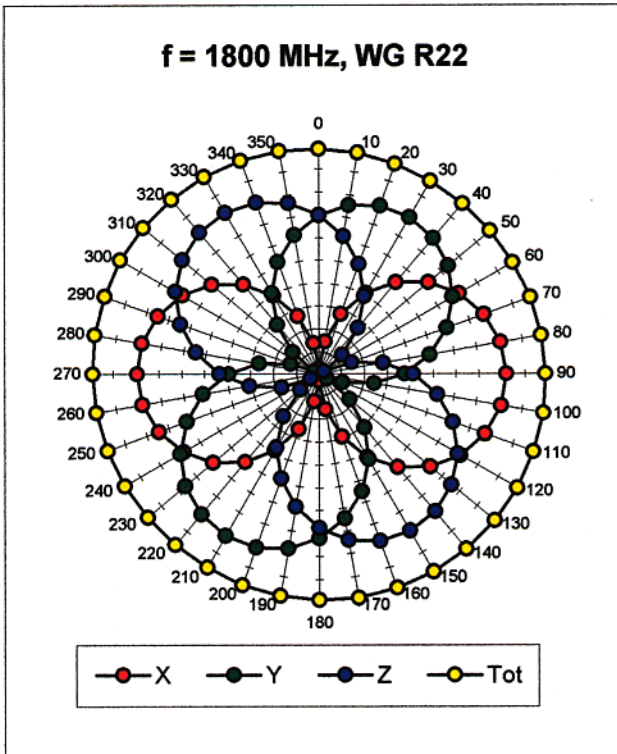

Isotropy Error (ϕ), $\theta = 0^\circ$

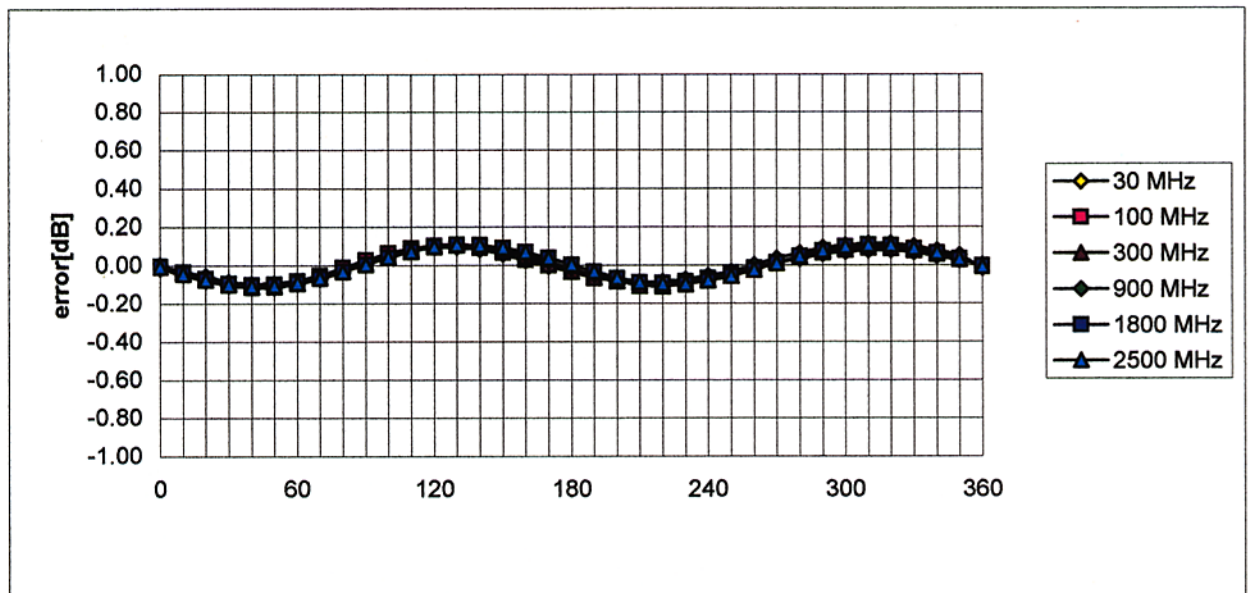

Frequency Response of E-Field
(TEM-Cell:ifi110, Waveguide R22)

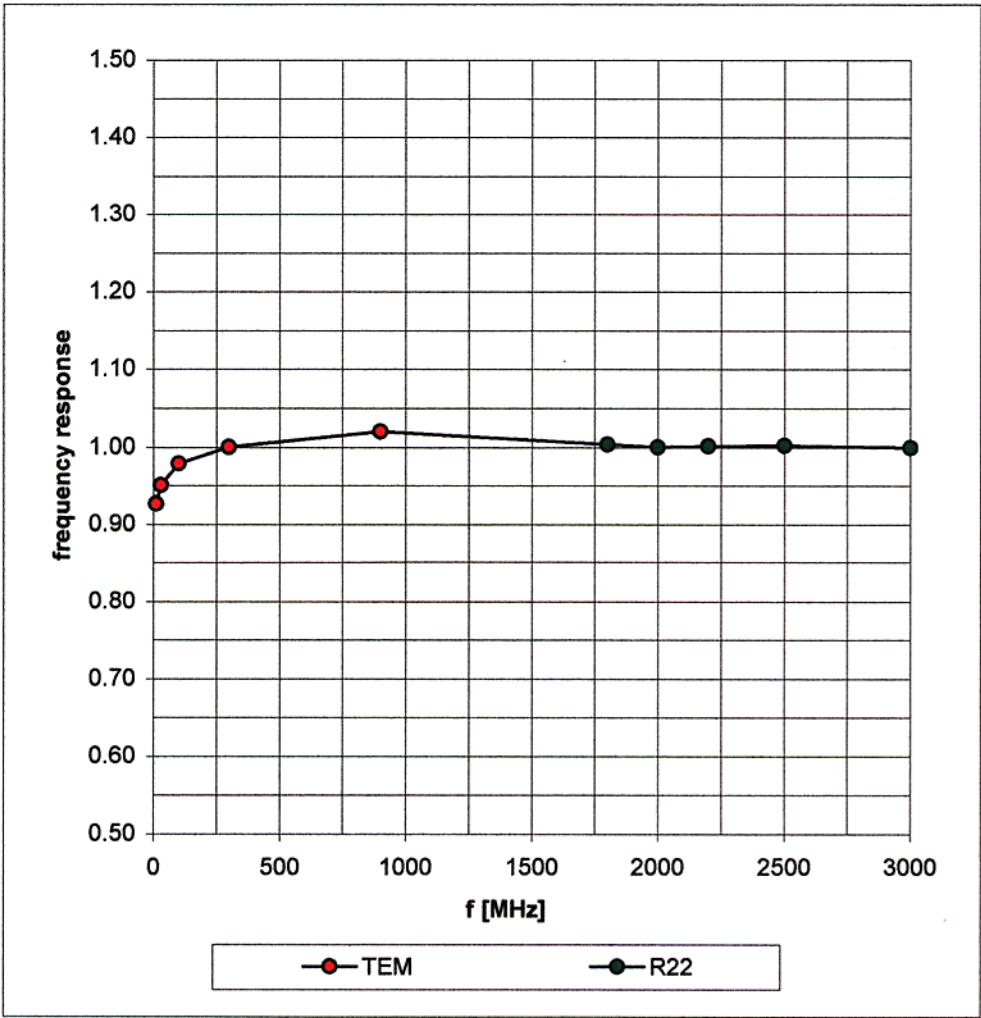

November 20, 2002

Dynamic Range f(SAR_{brain})

(Waveguide R22)

